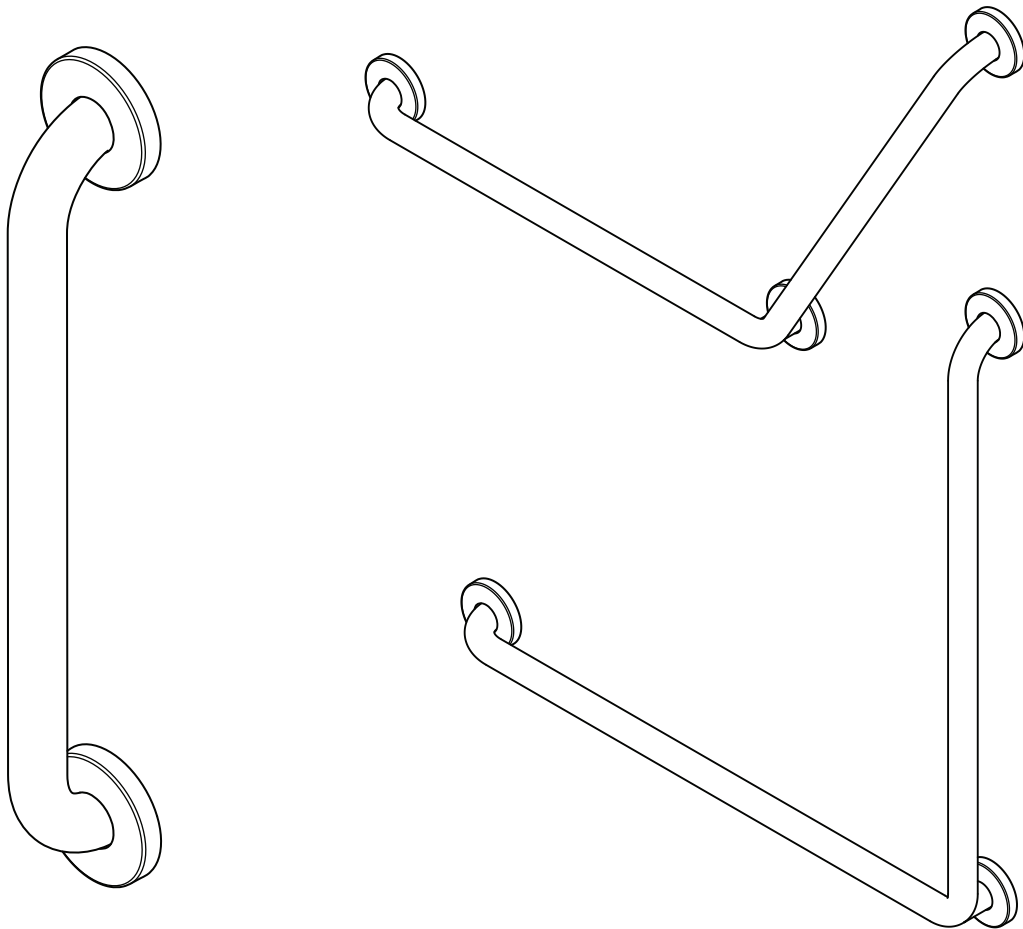
Easy Mount™ Grab Bar

GB-XX
(DC199 - RevE)

healthcraftproducts.com

OVERVIEW

Thank you for choosing the Easy Mount™ Grab Bar. Please read and understand the instructions in this manual; keep manual for future reference. It is your responsibility to see that your Easy Mount™ Grab Bar is properly installed and cared for. Failure to follow instructions in this manual could result in serious injury or death. If you are not equipped to undertake the outlined work, we would recommend that you have your Easy Mount™ Grab Bar installed by a qualified contractor.

The Easy Mount™ Grab Bar is intended for moderate vertical and horizontal load-bearing, providing sitting and standing support for people with disabilities. It should be installed with appropriate fasteners to a structural wall. The maximum user weight for this product is 500lbs / 227kg for bar lengths of 12", 18", 24" and L-shape, and 400lbs / 181kg for lengths of 32", 36", 42" and 120° grab bar.

Note: Product is not intended to support full body weight. The Easy Mount™ Grab Bar is not to be used in any other way than described above.

LIMITED LIFETIME WARRANTY

Products are covered by a limited lifetime warranty against defects in materials and workmanship for the original purchaser. Warranty excludes products that have been damaged through misuse, accidental damage, alteration, normal wear and tear, wood material and stain, or the use of corrosive or abrasive cleaning products.

Buyer hereby indemnifies, agrees to hold harmless and defend HealthCraft Products Inc. from and against any and all liabilities, claims, (founded and unfounded), losses, damages, costs and expenses (including without limitation consequential damages and reasonable professional fees) resulting from buyers specification, application, or improper use of goods described hereon; buyers omission or neglect. HealthCraft Products Inc. does not assume any liability for damage resulting from services performed by others or faulty installation, misuse or misapplication of goods sold by HealthCraft Products Inc. HealthCraft Products Inc. shall not be liable for prospective profits or special, indirect, or consequential damages, or for the cost of any corrective work done without HealthCraft Products Inc. prior approval. HealthCraft Products Inc. total liability hereunder shall in no event exceed the purchase price of the goods specified hereon.

Specifications subject to change without notice.

Manufacturer: HealthCraft Products Inc. 2790 Fenton Road, Ottawa, CA, K1T 3T7. 1-888-619-9992. Patent Pending. © 2026 HealthCraft Products Inc.

DELIVERY CONTENTS

FIGURE A.

1. Grab Bar
2. #10 x 2" wood screw, S/S (x6 for standard GB, x8 for L/120° GB)

INSTALLATION WARNINGS

FIGURE B.

1. Bar must be installed to walls that are structurally sound. Bar may be mounted horizontally across two studs, on an angle across two studs, or vertically on one stud.
- 2a. CAUTION - Do not install if there is no structure behind drywall or if there is no structure for one side.
- 2b. If desired installation location does not have adequate structure, add wood blocking OR
- 2c. To reinforce hollow walls, we recommend Easy Mount™ Anchors (sold separately). Weight capacity 300lbs.
- 3a. CAUTION - Do not install in fiberglass / acrylic shower enclosures with insufficient stiffness and strength.
- 3b. Add wood blocking if required.
- 3c. To reinforce fiberglass, we recommend Easy Mount™ Anchors (sold separately). Weight capacity 300lbs.
4. CAUTION - Do not install into metal studs.
5. CAUTION - Do not install if structure behind walls is unknown.

SUGGESTED PRODUCT LOCATION

FIGURE C.

For horizontal orientation of the grab bar, align the notch on the flange with the center of vertical studs – so that, ideally, 3 screws will contact a stud. To determine the optimal mounting location, we recommend holding the bar in place and walking through the motions of a transfer (sitting, standing, reaching, etc.). ADA recommended bar height of 33-36" / 838-914mm.

INSTALLATION - L SHAPED / 120° GRAB BARS

FIGURE D.

1. Mark location of the corner flange fasteners.
2. Install fasteners, leave enough room for plate.
3. Lower the bar onto the fasteners. Continue installing by following Figure E instructions.

INSTALLATION - GRAB BARS

FIGURE E.

1. Slide flange covers off and out of the way. Hold the bar against the wall, adjust using a level, and mark the screw hole locations.
2. Drill pilot holes for fasteners.
3. Fasten the bar to the wall using appropriate fasteners (min 2 per flange) and sealant.
4. Slide flange covers back onto flanges. Inspect the hardware monthly, and tighten if required.

PRODUCT USAGE

FIGURE F.

CAUTION - Grab bar and bathtub are slippery when wet.

TECHNICAL DATA

FIGURE G.

- *Complies with ADA 609.2 Cross Section - Grab Bars.
- *Complies with ADA 609.3 Spacing - Grab Bars.
- *Complies with ADA 609.8 Structural Strength - Grab Bars.

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