Service Manual D _ EN

Für / for

Badelifter / Seat-Lift SENTA-ST PUR L

Artikel Nr. / Part No. 920101001

Abbildung zeigt / Illustration shows
Aufrichthilfe / Seat lift SENTA-ST PUR L
Art. Nr. / Part No. 920101001

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Dear customer,

we would like to congratulate you on the purchase of a medical product from BEKA-Hospitec.

Please observe the following points:

1. Please read this manual carefully.
2. Please use this unit only for the purpose it is designed for.
3. The operation of this unit is only allowed by trained staff.
4. Maintenance and safety technical control checks must be carried out by authorized companies only.
5. Additional set ups, changes or repairs must be carried out by authorized companies only. We do not take any liability for additional set ups, changes or repairs as well as for the product or product functions, if they were carried out by non-authorized companies.

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The manufacturer of this unit is only responsible for the security and reliability of this unit, if:

- a. the unit is only used for the purpose it is made for,
- b. frequent functional control checks are made.

Caution:

**It is imperative that all
safety instructions be
observed.
(see User / Operating Manual for
SENTA-ST PUR L)!**

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For a safe operating you have to check the SENTA-ST-Seat-Lift PUR L once per day, see the following steps:

- Visual inspection (damages or wear).
- Check that all screws at the SENTA-ST PUR L lift are tight.
- Check that the column work properly.
- Check the Hand control (up / down, width movement).
- Check the emergency lowering (mechanical).
- Check the battery and charge function.
- Check the wheels (roll and turn freely, lock the brakes).

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When the SENTA-ST-Seat-Lift PUR L is not to be used for some time, the emergency stop should be pushed. Please store the SENTA-ST-Seat-Lift PUR L at a temperature (10°C up to 40°C) exceeding freezing point and at a relative humidity (30% up to 75%).

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A "safety check" (SC) must be performed on the SENTA-ST-Seat-Lift PUR L in intervals of **12 months**.

The performance of the safety checks and maintenance must be recorded at the device history book.

Please refer to the manufacturer after the inspection sticker on the back of the unit has expired.

The safety checks are only allowed to be conducted by the manufacturer or by persons authorized by the manufacturer.

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BEKA Hospitec provides you several variants on service agreements for regular safety checks of your BEKA product.

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Please contact at the end of the product life-cycle your BEKA-dealer, who will dispose of the SENTA-ST-Seat-Lift PUR L according to the local regulations and laws for recycling. For an environmentally sound disposal, the company Beka Hospitec GmbH as manufacturer gives more detailed information +49 (0) 641-92220-0.

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7.1 Wings

Picture 1:
Required tools:
1 x 4mm T-key

Picture 2:
Move the locking pin into the "Open" position and insert the two guide pins into the provided holes.

Picture 3:
Turn the locking pin until it falls into the "Locked" position.
Only then safe use can be guaranteed.


Picture 4:
Hang the locking rod into the grate at the back.

Picture 5:
Move the side wings into the desired position.

Picture 6:
Now insert the bolt through the strap holder for fixation and apply screw locking varnish on the bolt (**21 Nm min.**).

Picture 7:
Now screw the strap holder onto the side wing (left).

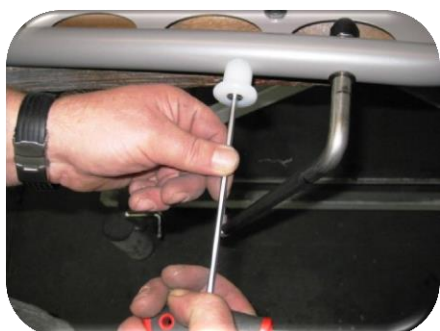Picture 8:
Now screw the strap holder onto the side wing (right).

Picture 9:
The illustration displays the installed side wing incl. the strap holder.

7.2 Handhold

Picture 1:Required tools:

- 1 x T-handle 4mm x 150
- 1 x screw locking varnish (medium strong **21 Nm**)

**Picture 2:**

Apply the screw locking varnish (**21 Nm min**) onto the thread.

**Picture 3:**

Insert the bolt through the handhold from behind and attach the spacer (black) to the screw.

**Picture 4:**

Lightly tighten the top screw onto the handhold.

**Picture 5:**

Now tighten the bottom screw lightly onto the handhold.

**Picture 6:**

Now tighten all screws onto the handhold.

7.3 Column

Picture 1:

Required tools:

- 1 x ratchet
- 1 x socket 13mm
- 1 x box wrench
- 1 x T-handle 6mm
- 1 x Hammer (plastic)
- 1 x Wedge (plastic)
- 1 x Allen key 2mm
- 1 x T-handle 4mm

Picture 2:

Remove the mounting panel with the control box as described at the user manual (14.4).

Picture 3:

Now you can remove the PVC cover.

Picture 4:

Remove all 4 screws (2 x left, 2 x right.) from the holding panel.

Picture 5:

Remove the holding panel.


Picture 6:

Now place the column with the framework onto a worktop.


Picture 7:

Bring the column with the framework into the desired work height.

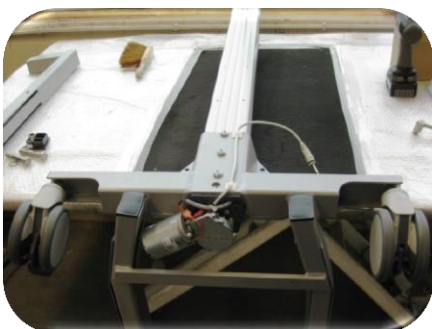Picture 8:

Remove the cable ties with a wire cutter.

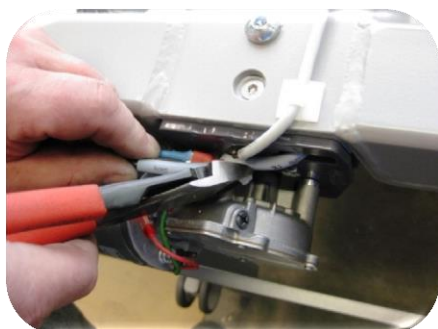Picture 9:

Remove the cable plug connectors from the motor.

Picture 9.1:
Remove the
cable plug
connectors
from the
motor.

Picture 9.2:
Remove the
cable plug
connectors
from the
motor.

Picture 9.3:
Remove the
cable plug
connectors
from the
motor.


Picture 10:
Remove the 3
screws from the
back of the
column.

Picture 10.1:
Remove the 3
screws from the
back of the
column.

Picture 11:
Now remove
the framework
from the
column.

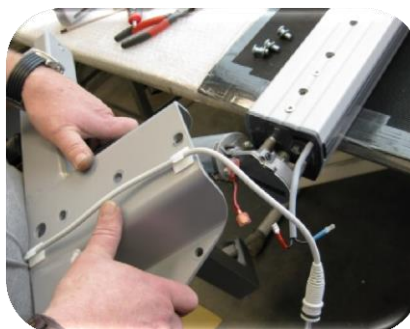Picture 12:
The illustration
shows the
removed column.

7.4 Microswitch change

Picture 1:

Required tools:

- 1 x T-handle 6mm
- 1 x Hammer (plastic)
- 1 x Wedge (plastic)
- 1 x Allen key 2mm
- 1 x T-handle 4mm

Picture 2:

Remove the 2 screws from the cover.

Picture 2.1:

Remove the 2 screws from the cover.

Picture 3:

Carefully remove the cover by lightly hitting it with a plastic wedge and a hammer.

Picture 3.1:

Carefully remove the cover by lightly hitting it with a plastic wedge and a hammer until you can take it off with your hands.


Picture 4:

Carefully remove the cover. Make sure to not tear out the plug connections.

Picture 4.1:

Carefully remove the cover. Make sure to not tear out the plug connections.

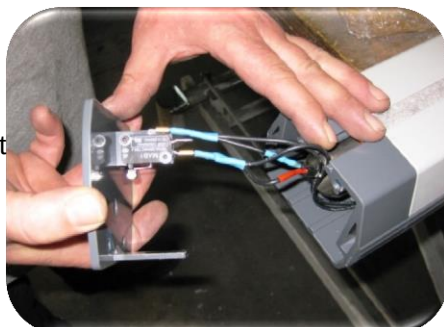
Picture 5:

Remove the 2 screws from the microswitch.

Picture 6:

Remove the 2 plug connections from the microswitch

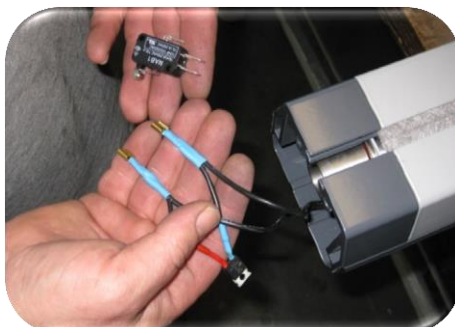
Picture 7:

Take the new microswitch and attach the two plug connections.

Picture 8:
Now carefully insert the cables into the column and attach the cover to where it was before. Make sure to not pinch any cables.

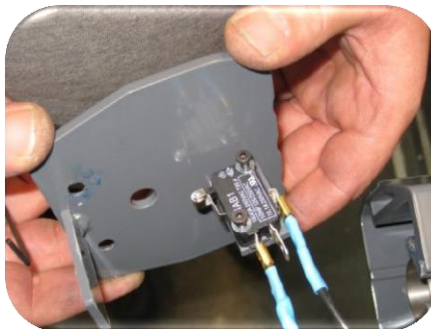7.5 Seat cover

Picture 1:

Required tools:

- 1 x drill, 4.5mm
- 1 x T-handle 2.5mm
- 1 x cordless screwdriver
- 1 x pencil

Picture 2:

Bring the seat into a height that allows easy access.

Picture 3:

Hold the seat cover and mark the 2 holes from the outside with a pencil.

Picture 4:

Here, you can see the 2 hole markings.

Picture 5.1:

Now use the drill to make the 2 holes.

Picture 5.2:

Now use the drill to make the 2 holes.

Picture 6:

Tighten the protective cover with a screw.

Picture 7:

Now you can mark the holes for the other side.

Picture 8:
After marking,
remove the
safety bolt.

Picture 9.1:
Now use the
drill to make
the 2 holes.

Picture 9.2:
Now use the
drill to make
the 2 holes.

Picture 10:
Now you can
attach the seat
cover with all
4 screws.

Picture 11:
The illustration
shows the
installed seat
cover.

7.6 Seat assembly

Picture 1 and 2:

Required tools:

- 1 x Ratchet
- 1 x Socket 8mm
- 1 x Allen key 8mm
- 1 x Screw locking varnish,
medium strong
- 1 x Turbo Grease spray

Picture 3:
Attach the new seat
socket to the lifter
and insert the 2
screws from one
side.

Picture 4:
Now attach the
washer from the
other side.

Picture 5:
Now attach the locking nut.

Picture 6:
Now firmly tighten the screws using the ratchet and the box wrench.


Picture 7:
Now firmly tighten the screws using the ratchet and the box wrench.

Picture 8:
Now bring the thread into shape once again with a knife to remove the old glue from the thread.

Picture 9:
Do this for all 4 threads.

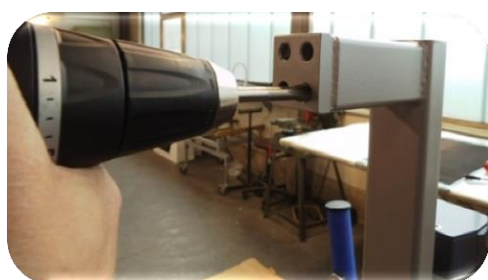Picture 10:
Now apply the new screw locking varnish to all 4 screw threads.

Picture 11:
The screw locking varnish must be at least medium strong (21 Nm).

Picture 12:
The illustration shows the 2 screws needed to attach the seat to the seat socket.


Picture 13:
Now lightly tighten the screw into the seat.


Picture 14:
Now lightly tighten the screw into the seat.

Picture 15:
Place a water level on the seat to adjust it.

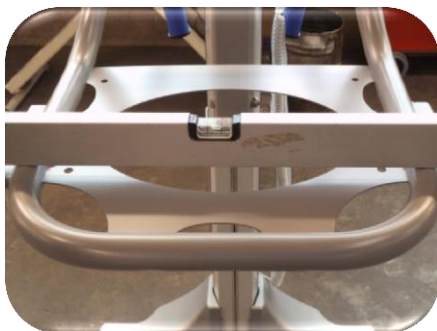Picture 16:
Now adjust the seat as indicated by the water level.

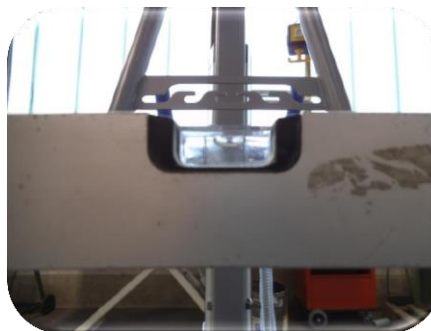
Picture 17:
Firmly tighten the top screw and make sure that the seat stays in position.

Picture 18:
Firmly tighten the bottom screw and make sure that the seat stays in position.

Picture 19:
Check the seat's position once again with the water level.

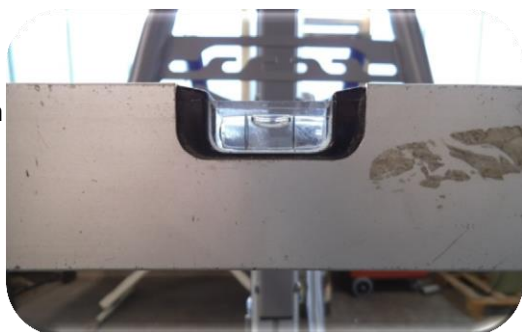Picture 20:
Before attaching the safety arms to the seat, Turbo Grease has to be applied on them.

Picture 21:
Also spray some Turbo Grease on the counterpart for taking in the safety arm.

Picture 22:
Now insert the safety arm into the seat socket. Make sure the two red markings face each other.

Picture 23:
The illustration shows an inserted safety arm.


Picture 24:
Now insert the other safety arm into the seat socket. Make sure the two red markings face each other.

Picture 25:

The illustration shows the two red markings as the safety arm gets inserted.

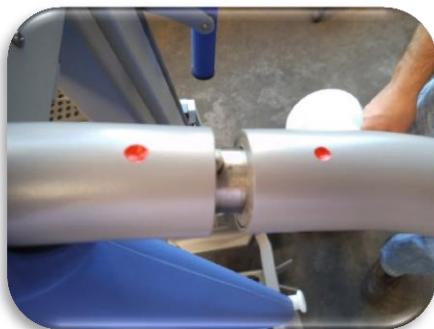Picture 26:

The illustration shows the inserted safety arms.


Picture 27:

Place the water level across the two edges of the arm for adjustment.


Picture 28:

The illustration shows the initial state of the safety arms.

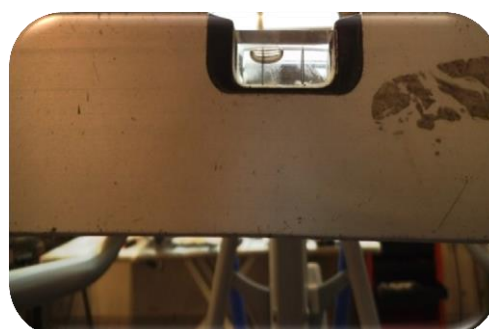Picture 29:

Adjust the safety arm by lightly pressing on the side which is too high until both sides are even.

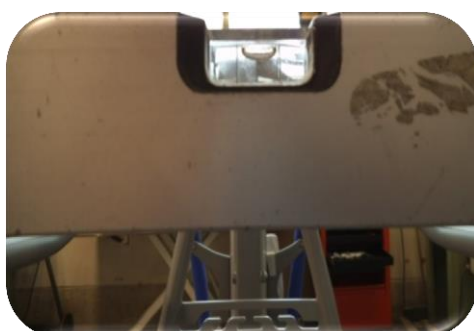7.7 Scale

Picture 1:**Required tools:**

- 1 x Allen key 8mm
- 1 x T-handle 2.5mm
- 1 x Ratchet
- 1 x Socket 8mm
- 1 x Socket 24 mm
- 1 x Socket 30mm
- 1 x Combination spanner 24mm

**Picture 2:**

Remove all 4 screws from the seat cover and take it off.

**Picture 3:**

Remove all 4 screws from the seat's framework and take off the seat.

**Picture 4:**

Remove all 2 screws within the green marked weighing cell.


**Picture 5**

During assembly make sure the arrow on the weighing cell always points DOWN (towards the floor).

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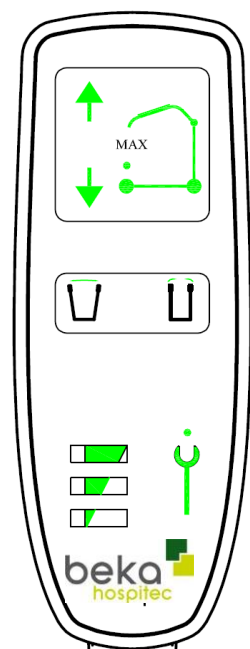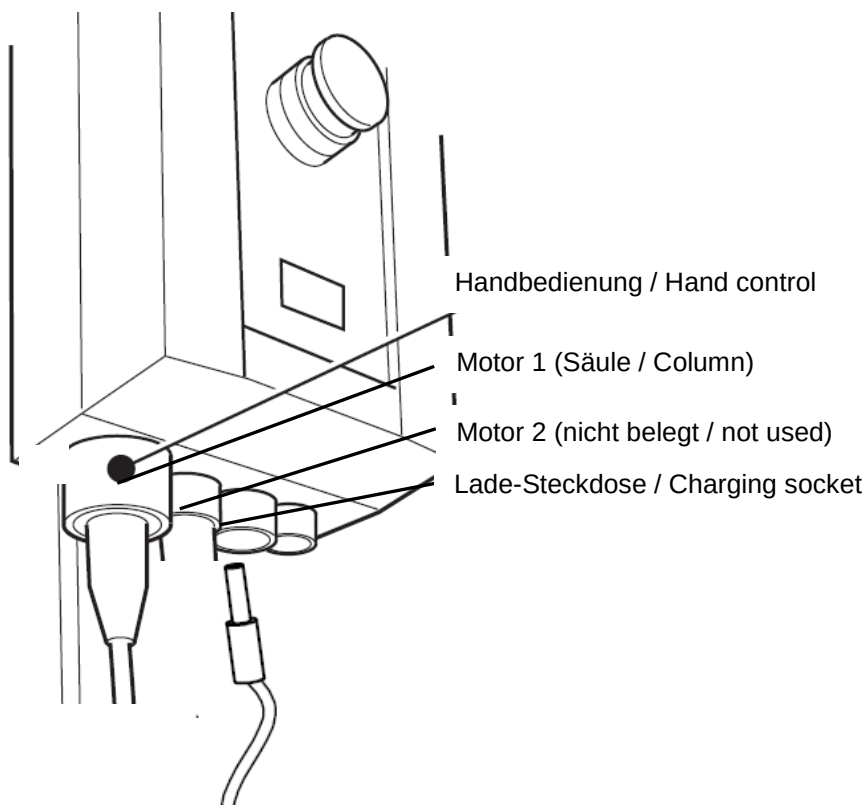
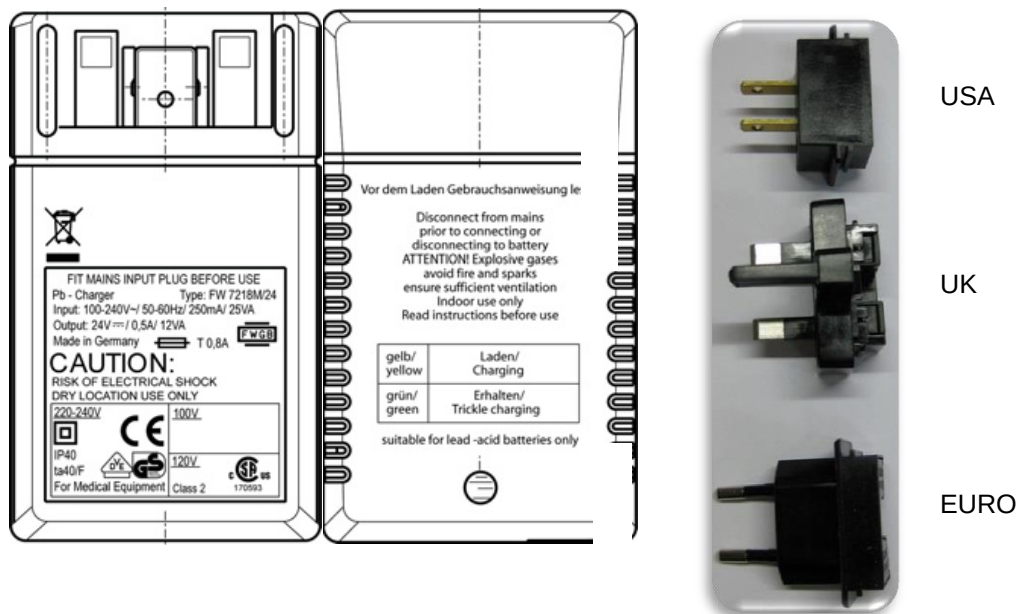
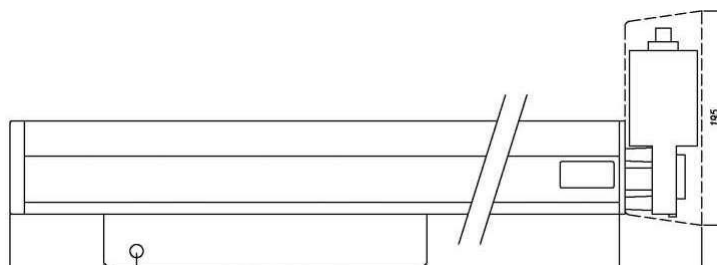
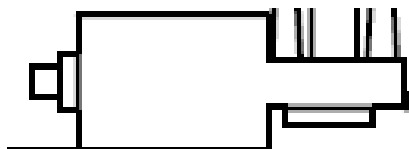
Typ of performance

Test before putting into service (reference value)

Device: Combi-bathing lift SENTA PUR L

SN:

Applied part type: B

Engineer

Paul Feringer

Typ:

Manufacturer: BEKA Hospitec GmbH

Mains connection: DPS

Protection class: I

Accessories:	SN:	YES	NO	Note
1 Scale		☐	☐	
1 Sling		☐	☐	

Visual and Functional Test	Pass	Fail	Note
1. Check mast and Sling suspension (dents, scratches, deformities, unusual surface wear)	☐	☐	
1.1 Check base (scratches, deformities, unusual surface wear, opening and closing)	☐	☐	
1.2 Check castor and castor brakes (turn freely, lock brakes, rolls straight on, cleaned)	☐	☐	
1.3 Screw connections (strength)	☐	☐	
1.4 Inspect the lift arm (lowers and raises freely, screws are secure)	☐	☐	
1.5 Inspect the Flexlink arm to Sling suspension (rotates freely, screws are secure)	☐	☐	
1.6 Inspect the Sling suspension (scratches, deformities, locating pin, rotates freely)	☐	☐	
1.7 Inspect the control box (LED, charging, all cords are in)	☐	☐	
1.8 Check handcontrol (check all functions)	☐	☐	
1.9 Check Emergency Stop (check function)	☐	☐	
1.10 Check Emergency lowering (electrical)	☐	☐	
1.11 Labeling (SN, REF, CE-Sign)	☐	☐	
1.12 Check max. Load (disconnection)	☐	☐	
1.13 Check Emergency lowering (mechanical)	☐	☐	
1.14 Check Moving handle (strength)	☐	☐	
1.15 Check electrical connections and seal up (Channel 2, connector for CH01-Charger)	☐	☐	
1.16 Check Scale (check all functions)	☐	☐	
1.17 Documentation (User Manual, etc.)	☐	☐	

Test Equipment:	SN:	Calibration Date:

Functional test:	set value	measured value	Complies	
			YES	NO
Run the lift with maximum load (automatic disconnection of supply)	185 kg ± 4 kg 408 lbs ± 8,8 lbs	kg / lbs	☐	☐

Overall assessment:	Next Service Date
☐ No safety or functional deficiencies were detected ☐ No direct risk, deficiencies detected may be corrected on short term ☐ Equipment does not comply – Taking out of service – is recommended	MM / YY /

Evaluation:	DATE	Signature Test Engineer

FB 19/01_EN Rev.03

PIE: Permanent installed equipment

NPS: Non-Detachable Power Supply Cord

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14.2 Troubleshooting

Problem with lifter	Relief
Lift does not work up and down.	a) Check that the emergency stop button is not pushed. b) Check that all cables are connected correctly. c) Check that the battery is charged. d) Check that the battery contacts are not defective.
Lift stops in the top position.	a) Check that the emergency stop button is not pushed. b) Check that the battery is charged. c) Use the emergency lowering (mechanical) to tower the patient.
Actuator "stalls" during lift.	a) Battery is low. Charge battery. b) The lifting capacity is exceeded (max. patient weight).
The electrical positioning adjustment does not respond.	a) Check if the hand control is connected. b) Check the plug-in connection on the control box.
Audible "beep" is heard from the control box.	Battery is low. Charge battery.
Hand control does not respond.	a) Check the connector on hand control cable. b) Check the battery condition (replace with a fully charged battery).
Raise and lower buttons on hand control do not respond.	a) Check that the emergency stop button is not pushed or pressed. b) Check the battery condition (replace with a fully charged battery). c) Check that the cables are properly connected at the control box.
If you hear unusual sounds.	Contact BEKA Service.
If the lift is damaged.	Contact BEKA Service.
Problem with charger	Relief
Battery charging does not work.	a) Check that the battery contacts are ok. b) Check the power plug.
"Power on" light on charger is not lit.	a) Check if the charger is connected to the wall receptacle. b) Check the fuse socket.
Charger is plugged in, but "Power on" light is not lit.	a) Check that there is power to the wall outlet. b) Remove the battery and check for damage.
Problem with battery	Relief
Battery is properly seated but no lights are visible.	Contact BEKA Service.
Indicator light does not go off after several hours of charging time.	Battery needs replacing. Contact BEKA Service.
Battery indicates it is fully charged when in the charger, but when placed in the lift, will only do a few lifts.	Battery needs replacing. Contact BEKA Service.

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Please contact the customer service of your local dealer or the manufacturer, if the SENTA-ST-Seat-Lift PUR L does not function properly and you could not eliminate the error with the relief actions according to section 12.2.

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Electrical data

Control unit

- Protection class:	class II
- Application part:	Type B
- Operation mode:	Non-stop operation
- Power supply:	24 V 50-60 Hz
- Power input:	max. 3,5 Amp.

Electric motor data

- Power supply:	24 V
- Frequency:	50-60 Hz
- Power consumption:	10 Amp.
- Power-on time max.:	2 Min ON and 18 Min OFF
- Protection type:	IPX4

Environmental conditions

Operation

- Temperature range:	10 °C up to 40 °C
- Relative air humidity:	35 % up to 75 %, non-condensing
- Atmosphärischer Druck:	800 – 1060 hPa

Storage and transport

- Temperature range (Lift):	-40 °C up to 70 °C
- Temperature range (Battery):	-15 °C up to 40 °C
- Relative air humidity:	10 % up to 80 %, non-condensing
- Atmospheric pressure range:	500 – 1100 hPa

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Kein Hausmüll
No domestic waste
Do not put the treatment system, and packaging materials
into the domestic waste.

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The declaration of conformity of the SENTA-ST-Seat-Lift PUR L will be provided by the manufacturer BEKA Hospitec GmbH upon your request.

15.4 Zubehör Lifter-System SENTA-ST PUR L / Accessories SENTA-ST-Seat-Lift PUR L

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Spare parts:

Part. No.	Description
K9200740	Manual switch
Z9200001	Middle element for lifter for lying patients, PUR
Z9200011	Seat element for hygiene opening
Z9200021	Back support, PUR
Z9200031	Cover for seat, PUR
Z9200041	Lying element for lifter for lying patients, PUR
K9300270	Steering roller double 125 mm back, lockable in position
K9200350	Steering roller double 100 mm front
K9200541	Control box CBJ
K9200390	Battery pack

Accessories:

Part. No.	Description
920103001	SENTA-ST PUR L side attachment wings
920104000	SENTA-ST PUR L digital scales SL
920103100	Side grid
920103120	Handle for lifter for lying patients
920103200	Foot support SENTA-ST PUR L PUR