Grand Tourer

by Travel Buggy

Portable Power Wheelchair

TABLE OF CONTENTS

NOTICES	3	TRANSPORTING	18
SAFETY INFORMATION	4	BATTERY MAINTENANCE	20
BATTERY SAFETY INFORMATION	5	CHARGING	22
SPECIFICATIONS	7	MAINTENANCE,CARE,TROUBLESHOOTING	23
UNBOXING	9	LIMITED WARRANTY	24
UNFOLDING	10	CONTACT & SUPPORT	25
FOLDING	11		
SET UP	12		
CONTROLLER	14		
HOW-TO FEATURES			
BACKREST HEIGHT ADJUSTMENT	16		
ARMREST HEIGHT ADJUSTMENT	16		
JOYSTICK FORWARD AND BACKWARD ADJUSTMENT	16		
FOOTREST HEIGHT ADJUSTMENT	16		
FREE WHEEL	17		

⚠ CAUTION:Read all precautions and instructions in this manual before using this equipment.
Improper use or maintenance can void the warranty.Keep this manual for future reference.

NOTICES

Thank you for choosing Grand Tourer by Travel Buggy, a convenient and portable powerchair.

It is important to read through this manual before operating your Grand Tourer to become familiar with its features as well as its limitations and safety information.

You will find the operation and care of your Grand Tourer very simple and easy, though should you have any questions, please contact Travel Buggy.

OUR CUSTOMER HOTLINE

Canada: 1.888.56.BUGGY(28449)

USA: 1.855.56.BUGGY(28449)

EMAIL

contact@travelbuggy.com

DISCLAIMER:

The notes, warnings, and cautions contained within the manual and marked by the triangular Caution Symbol ⚠ should be given special care. Users should also pay special attention to information marked in this manual beginning with **NOTICE**.

Travel Buggy, its affiliates and partners cannot be held responsible for personal injury or property damage resulting from the unsafe or improper operation or maintenance of the Grand Tourer Power wheelchair.

As part of our ongoing product improvement initiative, Travel Buggy reserves the right to change specifications and design without notice. As a result, there may be minor differences between your Grand Tourer and the photos, illustrations and instructions in this manual.

SAFETY INFORMATION

- DO make certain that the seat is in a locked position when you are operating your GRAND TOURER
- DO reduce your Grand Tourer speed when you are driving around a corner
- DO proceed with caution as you approach the downgrade of a ramp or other incline and when driving near a raised surface or unprotected ledges and drop offs (curbs, porches, stairs etc.)
- DO avoid exposure to rain snow, ice, salt, or standing water, whenever possible
- DO maintain and store in a clean and dry condition
- DO keep your hands away from the tires when driving. Be aware that loose fitting clothing can become caught in drive tires
- DO keep tools and other metal objects away from battery terminals. Contact with tools can cause electrical shock
- DO follow the recommended weight capacity (300lbs /136 kg)
- DO follow the recommended incline slope rating (10 degrees)
- DO check the battery charge to ensure there is enough charge to get to your destination before using.
- DO ensure your Grand Tourer power is OFF when entering and exiting the power wheelchair
- DO NOT drive your Grand Tourer across the side of a hill or diagonally up or down a hill, and do not stop (if possible) while driving up or down an incline
- DO NOT operate your new GRAND TOURER for the first time without completely reading and understanding the owner's manual
- DO NOT carry additional passengers on your GRAND TOURER
- DO NOT use your GRAND TOURER as a seat in a moving vehicle
- DO NOT allow unsupervised children to play near the power wheelchair
- DO NOT operate your Grand Tourer while under the influence of alcohol, as this may impair your ability to drive safely
- DO NOT sit in the Grand Tourer while in a moving vehicle, as an accident or sudden stop may cause you to be thrown from the power wheelchair and compromise the seat lock latch
- ALL servicing must be performed by qualified personnel only. Unauthorized modifications to power wheelchairs, its power system or its control system may compromise the performance and safety of the system. UNDER NO CIRCUMSTANCES should you attempt any repair or service that is not covered in this manual.
- ALWAYS ensure all safety warning labels are in place and legible.

⚠ WARNING: The operator should read all the warning signs and instructions here and on the power wheelchair before use.

⚠ WARNING: You should not operate your power wheelchair on public streets and roadways. Be aware that it may be difficult for traffic to see you when you are seated on your power wheelchair. Obey all local pedestrian traffic rules. Wait until your path is clear of traffic and then proceed with extreme caution.

BATTERY SAFETY INFORMATION

WHEN USING THE BATTERY

1. Lithium-Ion battery packs may get hot, explode or ignite and cause serious injury if exposed to abuse conditions. Be sure to follow the safety warnings listed below:
 - DO NOT place the battery near fire or heat the battery
 - DO NOT install the battery backwards so the polarity is reversed
 - DO NOT connect the positive and negative terminals of the battery to each other with any metal object(such as a wire)
 - DO NOT carry or store the battery together with necklaces, hairpins, or other metal objects
 - DO NOT pierce the battery with nails, strike the battery with a hammer, step on the battery or otherwise subject it to strong impacts or shocks
 - DO NOT solder directly on to the battery
 - DO NOT expose battery to water or salt water, or allow the battery to get wet
2. DO NOT disassemble or modify the battery. The battery contains safety and protection devices, which if damaged, may cause the battery to generate heat, explode or ignite
3. Do NOT place the battery in or near fire, on stoves or other high temperature locations. Do not place the battery in direct sunlight, or use or store batteries inside cars in hot weather. Doing so may cause the battery to generate heat, explode or ignite. Using the battery in this manner may also result in a loss of performance and a shortened life expectancy.
4. If the device is to be used by children the caregiver should explain the contents of this document to the children and provide adequate supervision to ensure the device is being used appropriately
5. When the battery is worn out, insulate the terminals with adhesive tape or similar materials before disposal
6. Immediately discontinue use of the battery if, while using, charging or storing the battery, the battery emits an unusual smell, feels hot, changes color, or shape, or appears abnormal in any way. Contact Travel Buggy if any of these problems are observed
7. DO NOT place battery in microwave oven, high pressure containers or on induction cookware

BATTERY SAFETY INFORMATION



CHARGING PRECAUTIONS

1. BE SURE to follow the warnings listed below while charging the battery. Failure to do so may cause the battery to become hot, explode or ignite and cause serious injury.
 - When charging the battery, ONLY USE THE SPECIFIED CHARGER PROVIDED WITH THE PRODUCT
 - DO NOT attach the battery to a power supply plug or directly to a car's cigarette lighter
 - DO NOT place the battery in or near fire, or in direct sunlight. Heating the battery can damage the safety circuitry, which can cause additional heating, rupture or ignition of the battery
2. DO NOT continue charging the battery if it does not recharge in the timelines specified in the operating manual (4-6 hours or until the charging brick turns green)
3. THE TEMPERATURE range over which a battery can be charged is 0°C to 45°C (32°F to 113°F). Charging the battery at temperatures outside of that range may cause severe damage to the battery or reduce battery life expectancy

WHEN DISCHARGING THE BATTERY

1. DO NOT discharge the battery using any other device except for the specified device. When the battery is used in devices other than the specified device, it may damage the battery or reduce its life expectancy. The device may cause abnormal current flow, which can cause the battery to become hot, explode or ignite and cause serious injury.
2. THE TEMPERATURE range over which the battery can be discharged is -20°C to 60°C (-4°F to 140°F). Use of the battery outside this temperature range may damage the performance of the battery or may reduce its life expectancy.

SPECIFICATIONS

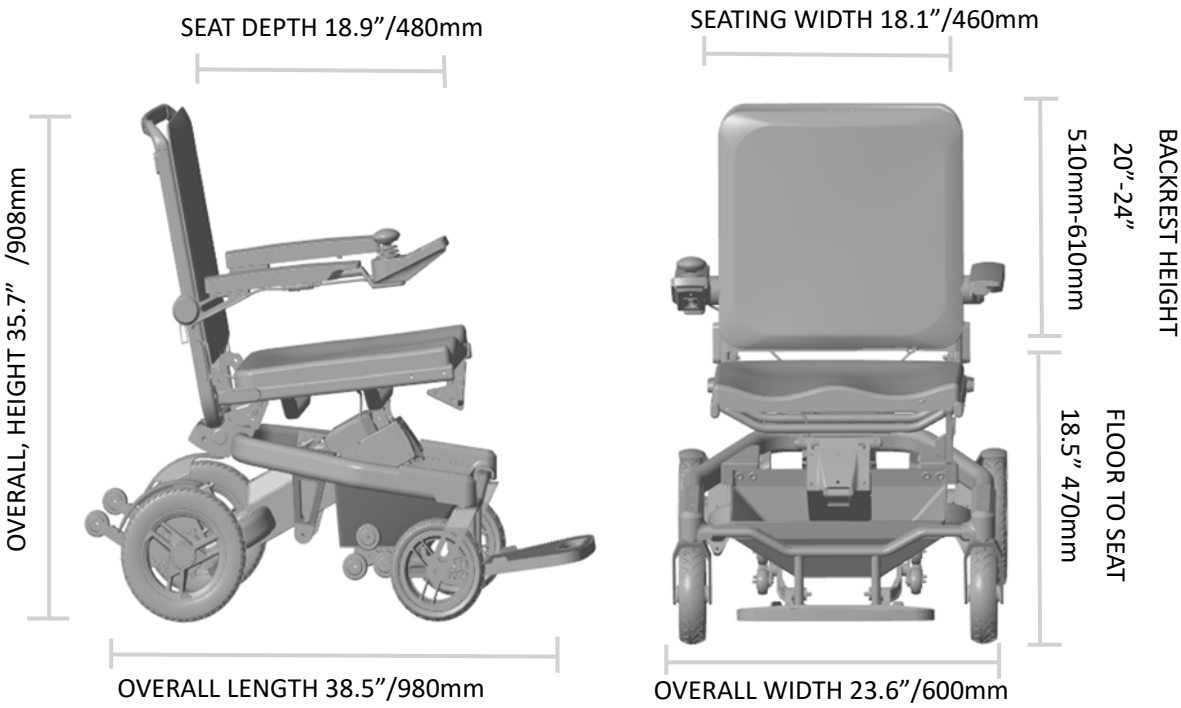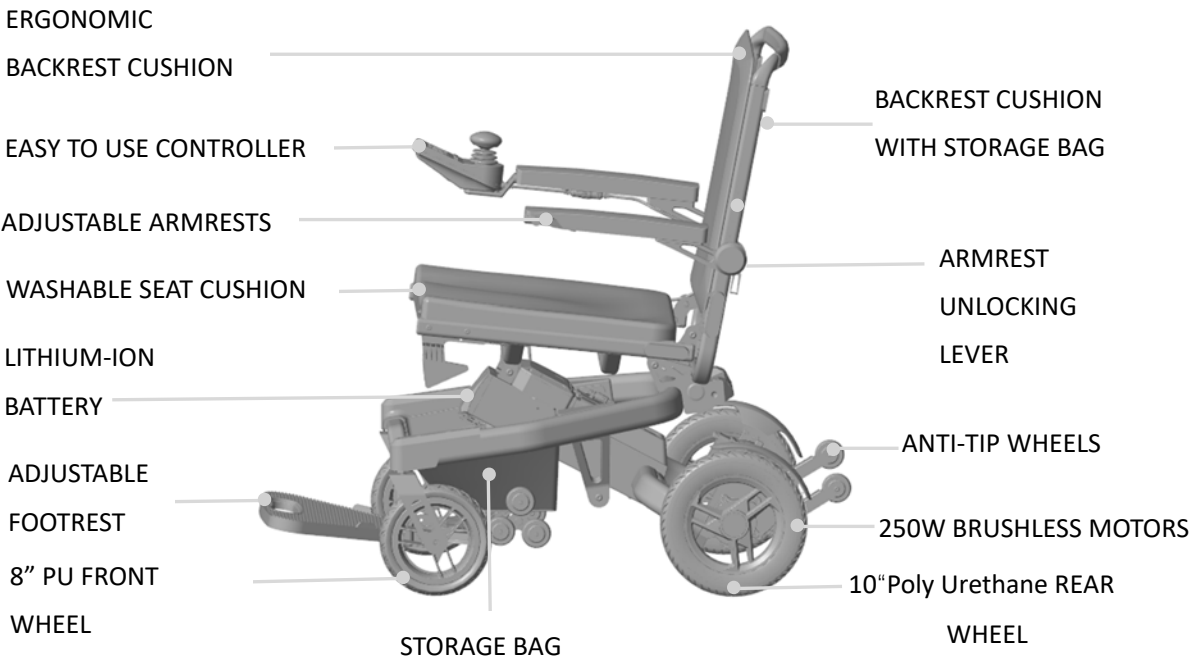SPECIFICATIONS

PERFORMANCE SPECIFICATIONS

DRIVING RANGE: 15.53MI/25KM

MAXIMUM SPEED: 3.73MPH/6KPH

DRIVE SYSTEM: TWO MOTOR DRIVE SYSTEM

MOTOR TYPE:250 WATT BRUSHLESS

INCLINE CAP ABILITY:10°

TURNING RADIUS:31.5"/800mm

STOPPING DISTANCE:3.28FT/1m
(AT FULL SPEED)

OBSTACLE CLEARANCE:1.57"/40mm

BRAKING SYSTEM: ELECTROMAGNETIC

MAX WEIGHT CAPACITY:300LBS/136KG

DIMENSIONAL SPECIFICATIONS

UNFOLDED DIMENSIONS: 38.5"*23.6"*35.7"
980mm*600mm*908mm

CHAIR GROSS WEIGHT: 55LBS / 25KG

CHAIR NET WEIGHT: 50.7LBS /23KG
(NO BATTERY)

SEAT CUSHION SIZE:
18" (W)*18(D)"/450mm(W)*450mm(D)

FOLDED DIMENSIONS: 29.5" *23.6"*13.8"
750mm*600mm*350mm

SEATING WIDTH:18.11"/460mm

FRONT WHEEL DIAMETER: 8"

REAR WHEEL DIAMETER: 10"

FLOOR TO SEAT HEIGHT: 18.5"/470mm

BACK SUPPORT HEIGHT:
19.7"/21.65"/23.62" 500mm/550mm/600mm

OTHER SPECIFICATIONS

BATTERY TYPE: Lithium-Ion (24VX12AH)

BATTERIES PROVIDED:ONE (1)

BATTERY WATT HOUR:288WH

BATTERY CHARGE CYCLES: ~500

TIRES: PU

CHARGING TIME:6 TO 7 HOURS


CHARGER INPUT:100-240V,50/60HZ
(INTERNATIONAL)

CHARGER OUTPUT:24V 2.0A

LAP BELT: PROVIDED

HORN: BUILT-IN

UNBOXING

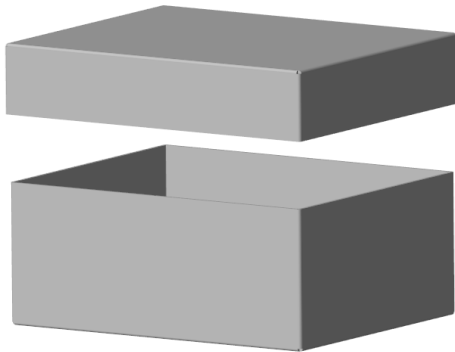
- To remove the powerchair, orientate the box, so it is laid down on the ground
- Grab the top and begin to pull your Grand Tourer out horizontally parallel with the ground.

- Locate the components provided with your power wheelchair

SEAT CUSHION

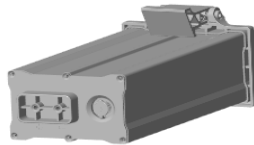LITHIUM BATTERY

JOYSTICK

BATTERY CHARGER


TOOL KIT

TOOL

UNFOLDING

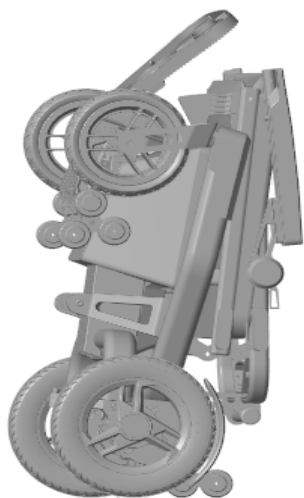1-Lift the footrest and pull the safety lock upward to unlock the chair.

2-Unfold the wheelchair outward until you hear a clicking sound indicating it is locked in place.

3-Lower the anti-tip wheel bracket, complete the unlocking of the wheelchair.

⚠ NOTICE: Failure to fully extend the power chair into its locked position may result in injury while driving.


FOLDING

1-Lift both the armrests and the anti-tip wheel brackets upward.


2-Pull the unlocking handle and push the wheelchair forward.

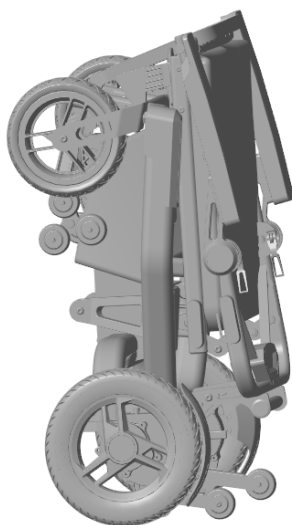
3-Until you hear a clicking sound indicating it is locked in place.

SET UP

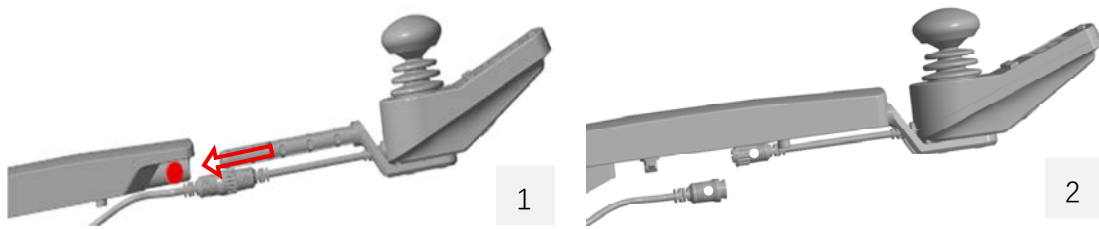- Slide the joystick controller into the armrest tube until a comfortable arm length is achieved for the user (image 1).
- Align the grey points of the joystick connection cables and press them together to plug the joystick controller in (image 2).
- Rotate the coupling over the connection point to secure the connection point
- Make sure the quick-release pin locks into the tube to complete the fixing action.



SET UP

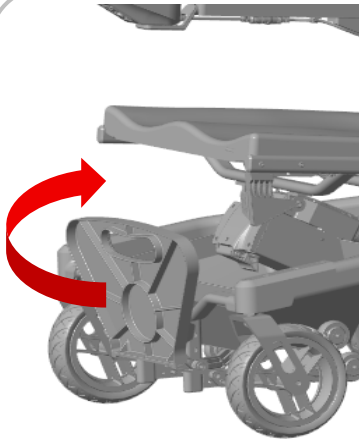- Before sitting in the electric wheelchair, be sure the power button is OFF and the rear brake levers are in the locked position
- GRAND TOURER offers the user 2 optional entry alternatives, direct front entry or side entry
- For a front entry, we recommend flipping up the foot plate to create an easier sitting experience
- DO NOT step on the footplate when entering or exiting the GRAND TOURER as it may cause the power chair to move in an unsafe manner


- To enter on the side, lift the arm rest vertical to provide the required clearance
- To enter on the side, lift the arm rest vertical to provide the required clearance

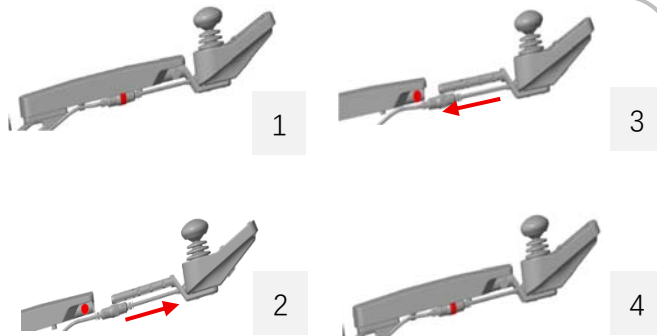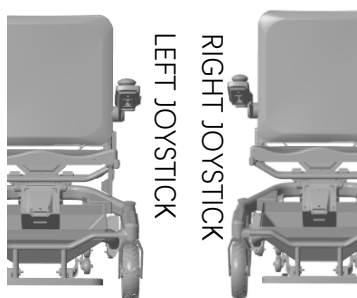- Joysticks can be positioned on left or right arm rest
- Rotate the joystick connector to loosen it (Image 1).
- Pull up the quick-release pin to remove the joystick (Image 2).
- Insert the joystick fixing plate into the armrest tube on the other side. Adjust to the desired position and lock the quick-release pin (Image 3).
- Insert the joystick connector (Image 4) and rotate the locking ring to secure it. Ensure the quick-release pin is fully locked into the groove to complete the installation.
- Move the cable to the new side and secure it with ties. Make sure it doesn't block folding or movement, then reconnect the joystick.

CONTROLLER

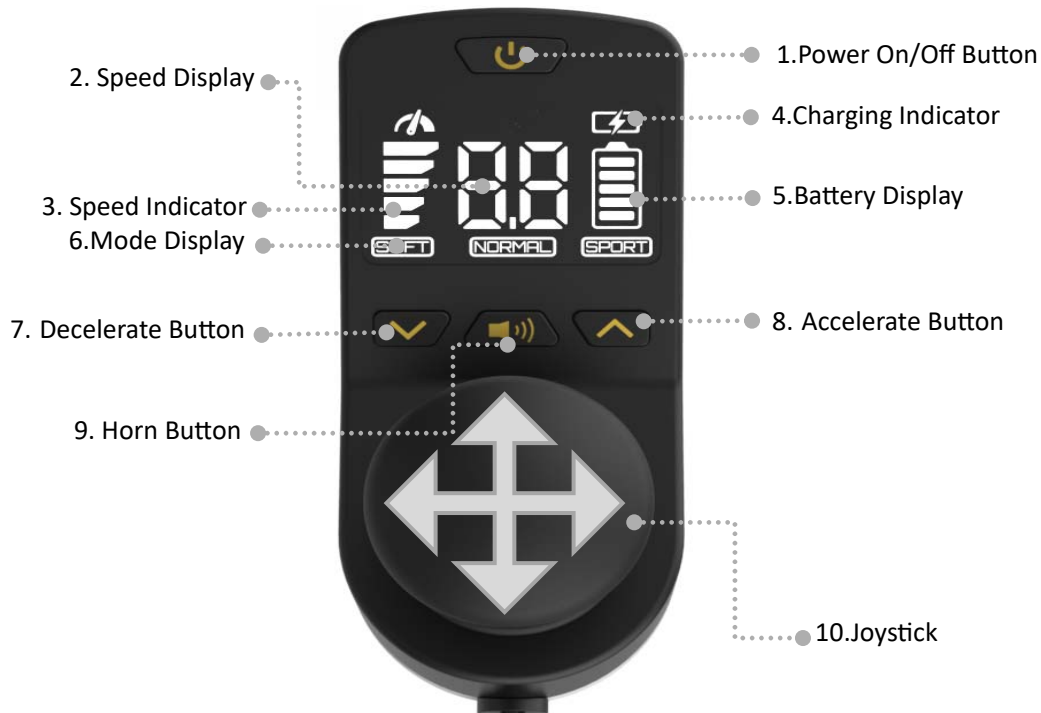1. **Power On/Off Button** — Press this button briefly to power on the electric wheelchair and hold it down to power off.

2. **Speed Display** — Displays the current driving speed of the wheelchair in real time. As you accelerate or decelerate, the numbers on the screen will change accordingly, helping you easily monitor your current speed.

3. **Speed Indicator** — One bar represents the slowest speed, and five bars represent the fastest speed.

4. **Charging Indicator** — When the device is connected to the charger, the light will blink during charging; The indicator light will remain steadily on after charging is complete.

5. **Battery Display** — Displays the real-time battery level of the electric wheelchair.

6. **Mode Display** — There are a total of three modes: SOFT, NORMAL, SPORT.

7. **Decelerate Button** — Press this button to decrease speed. Hold the brake button down, and the driving mode of the wheelchair will cycle from SPORT to NORMAL, and finally to SOFT. You can adjust based on your preference.

8. **Accelerate Button** — Press this button to increase speed, with a maximum speed of 6 kilometers per hour. Hold the accelerate button down, and the driving mode will change from SOFT to NORMAL, and finally to SPORT.

9. **Horn Button** — Press this button to activate a beeping sound from the electric wheelchair to alert the surrounding environment.

10. **Joystick** — Controls the direction and speed of the electric wheelchair; push the joystick in the direction you wish to go.

CONTROLLER

11. **Type-C Power Supply Port** — Use a power cable with a 5V/2A output specification to connect mobile devices for charging.

12. **Wheelchair Charging Port**— You can charge the electric wheelchair here without removing the battery.

HOW-TO FEATURES

1

2

Backrest Height Adjustment

By loosening the four set screws on the top of the backrest (Image 1), the height of the backrest can be freely adjusted. The adjustable range is 100mm (Image 2).


3

Armrest Height Adjustment

By loosening the quick-release handle on one side of the armrest, you can freely adjust the height of the armrest easily, the adjustable range is 100mm, based on your preference (Image 3).


4

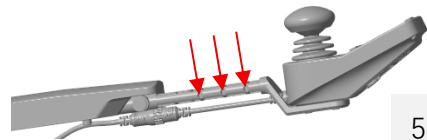5

Joystick Forward and Backward Adjustment

Using the quick-release button located at the base of the armrest (Image 4), you can easily pull it to quickly adjust the armrest to three different positions (Image 5). The total adjustable range is 50mm.

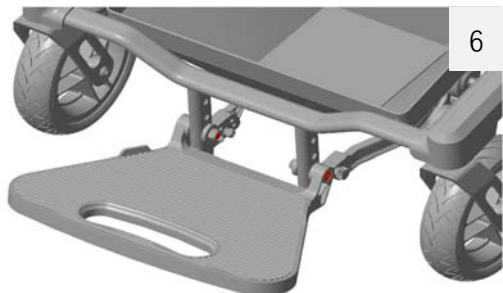6


Footrest Height Adjustment

By loosening the two screws located at the footrest (Image 6), you can adjust the height of the footrest to three different levels. The adjustable range is 50mm.

⚠ WARNING: Always make certain the wheels are properly inserted back into position and that both the teeth of the locking plate and the locking pin are in place before using your Travel Buggy.


HOW-TO FEATURES

PUSH MODE

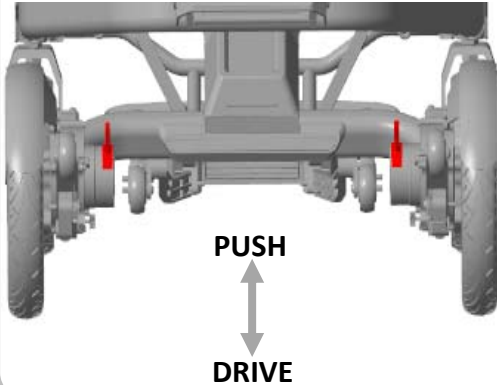
- In order to have a care giver push the Grand Tourer, it must be put in PUSH mode
- This mode disengages the drive mechanism to the wheels, making them free to push manually

LOWER ANTI-TIP WHEELS: Please ensure the anti-tip wheel brackets are in the lower position during use to help reduce the risk of tipping and improve overall safety.

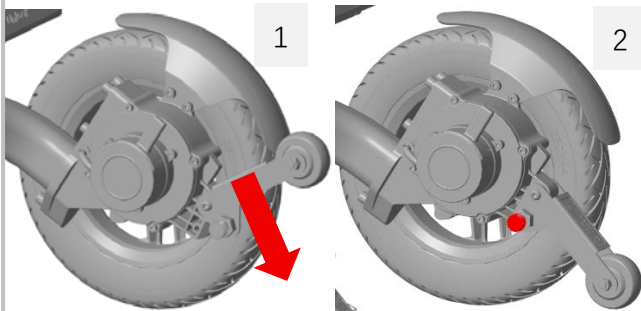- Press the anti-tip wheel bracket downward; the spring pin will automatically slide into the built-in round recessed hole on the motor (Image 1).
- securely locking the anti-tip wheel into position (Image 2).

FOLD ANTI-TIP WHEELS

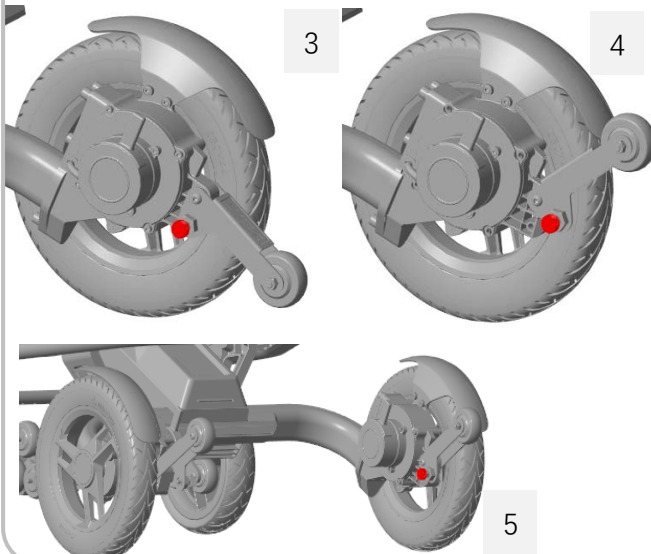- Pull the anti-tip wheel unlock button (Image 3). The anti-tip wheel will automatically fold upward via the internal tension spring.
- It will then automatically lock into place through the spring mechanism within the unlock button (Image 4&5)

TRANSPORTING

Towing Function Guide

After folding the wheelchair, you can easily tow your electric wheelchair by following these steps:

1. Flip up the footrest (Image 1).
2. Grab the hollow part of the footrest (Image 2), and using the rear wheels, you can easily move your electric wheelchair.


TRANSPORTING


Trunk Loading Guide

After folding the wheelchair, you can easily load it into the trunk of your car by following these steps:



1. Move the electric wheelchair to the rear of your car's trunk (Image 1), ensuring that the crawling wheels are positioned close to the trunk.
2. For user convenience, the climbing wheels are designed with two adjustable positions to accommodate both sedans and full-size SUVs. You can select the appropriate angle based on the height of your vehicle's trunk.

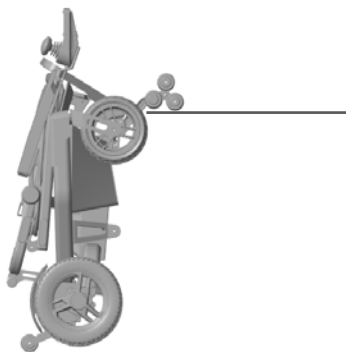3. Pull the unlock pin to open the crawling wheels (Image 2, Image 3). The unlock pin has two different angles to accommodate different trunk heights. After adjusting to the first angle, you can continue adjusting to a larger angle.

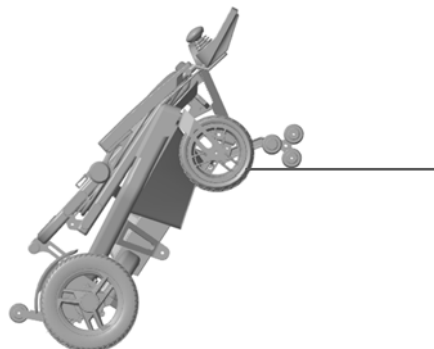3. Grasp both sides of the bottom of the backrest tube (Figure 8.2.4) so that the crawling wheels contact the flat surface of your car. Then, hold both sides of the bottom of the backrest frame (Figure 8.2.5) and lift your electric wheelchair upwards. Push it towards the interior of your car's trunk to successfully load it into the trunk.

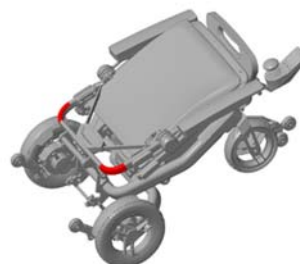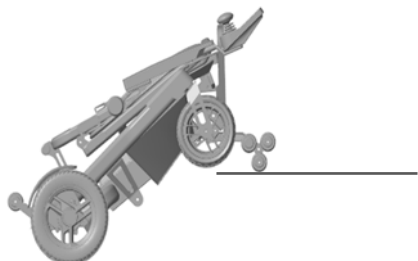BATTERY MAINTENANCE

Lithium-Ion rechargeable batteries require routine maintenance and care in their use and handling. Read and follow the guidelines in this document to safely use Lithium-Ion batteries and achieve the maximum battery life span.

OVERVIEW

Do not leave batteries unused for extended periods of time, either in the product or in storage. When a battery has been unused for 6 months, check the charge status and charge or dispose of the battery as appropriate.

The typical estimated life of a Lithium-Ion battery is about 3 years or ~500 charge cycles, whichever occurs first. One charge cycle is a period of use from fully charged, to fully discharged, and fully recharged again. Use a three-year life expectancy for batteries that do not run through complete charge cycles.

Rechargeable Lithium-Ion batteries have a limited life and will gradually lose their capacity to hold a charge. This loss of capacity (aging) is irreversible. As the battery loses capacity, the length of time it will power the product (run time) decreases.

Lithium-Ion batteries continue to slowly discharge(self-discharge) when not in use or while in storage. Routinely check the battery's charge status.

Use only the supplied battery in your product.

Observe and note the time that a new fully charged battery provides for powering your product. Use this new battery run time as a basis to compare run times for older batteries. The running time of your battery will vary depending on the product's configuration and the applications that you run.

Routinely check the battery's charge status.

Carefully monitor batteries that are approaching the end of their estimated life.

Consider replacing the battery with a new one if you note either of the following conditions:

- The battery run time drops below about 80% of the original run time.
- The battery charge time increases significantly.

If a battery is stored or otherwise unused for an extended period, be sure to follow the storage instructions in this document. Not follow the instructions, and the battery has no charge remaining when you check it, consider it to be damaged. Do not attempt to recharge it or to use it. Replace it with a new battery.

CHARGING

Always follow the charging instructions provided with your product. Refer to your product's user manual and/or online help for detailed information about charging its battery.

Use only the supplied charger to charge your battery.

BATTERY MAINTENANCE

STORAGE

- Charge or discharge the battery to approximately 50% of capacity before storage.
- Charge the battery to approximately 50% capacity at least once every six months.
- Remove the battery and store it separately from the product.
- Store the battery at temperatures between 5 °C and 20 °C (41 °F and 68 °F).

HANDLING PRECAUTIONS

- Do not disassemble, crush, or puncture a battery.
- Do not short the external contacts on a battery.
- Do not expose a battery to fire or water.
- Do not expose a battery to temperatures above 60°C (140 °F).
- Keep the battery away from children.
- Avoid exposing the battery to excessive shock or vibration.
- Do not use a damaged battery.
- If a battery pack has leaking fluids, do not touch any fluids.
- Dispose of a leaking battery pack (see Disposal and Recycling in this document)
- In case of eye contact with fluid, do not rub eyes. Immediately flush eyes thoroughly with water for at least 15 minutes, lifting upper and lower lids, until no evidence of the fluid remains. Seek medical attention.
- Immediately discontinue use of the battery if, while using, charging or storing the battery, the battery emits

an unusual smell, feels hot, changes color, or shape, or appears abnormal in any way. Contact Travel Buggy if any of these problems are observed.

TRANSPORTATION

- Always check all applicable local, national and international regulations before transporting a Lithium-Ion battery.
- Transporting an end-of-life, damaged, or recalled battery may, in certain cases, be specifically limited or prohibited.

DISPOSAL AND RECYCLING

Lithium-Ion batteries are subject to disposal and recycling regulations that vary by country and region. Always check and follow your applicable regulations before disposing of any battery. Contact Rechargeable Battery Recycling Corporation (www.rbrc.org) for U.S.A. and Canada, or your local battery recycling organization.



Many countries prohibit the disposal of waste electronic equipment in standard waste receptacles.

Place only discharged batteries in a battery collection container. Use electrical tape or other approved covering over the battery connection points to prevent short circuits.

CHARGING

Battery Charging Methods

Grand Tourer is equipped with dual charging ports, located beneath the controller (Image 1) and at the bottom of the battery (Image 2). You can choose either charging method according to your needs.

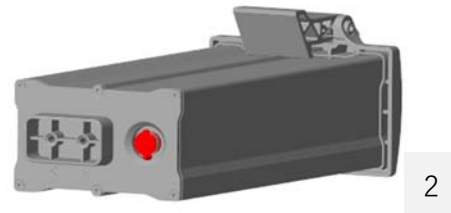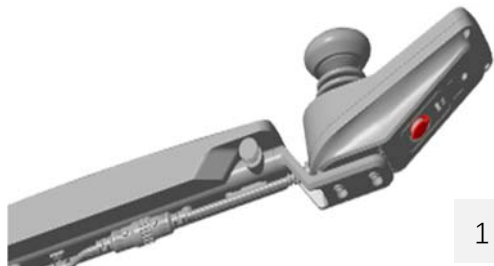- Fully Charged
- Charging

- When charging, the indicator light on the power charger is red.
- After 6 to 7 hours of charging (depending on the remaining battery level), the indicator light will change from red to green, indicating that the battery is fully charged.
- To check your battery level, turn ON the joystick controller power.
- There are a total of five WHITE battery indicator lights.
- From top to bottom, the lights will gradually turn off as the battery drains, indicating the remaining charge.
- Recharge the battery before the last WHITE light turns off, as this indicates the battery is fully depleted.

MAINTENANCE, CARE, TROUBLESHOOTING

It is important to perform regular maintenance of your Grand Tourer for longevity and warranty eligibility. Even a little cleaning now and again will go a long way to keeping your device looking and functioning like new.

Certified repair shops with trained professionals who specialize in electric wheelchairs will be qualified to perform the service and maintenance on your Grand Tourer electric wheelchair should it be required. Please be sure to keep a log of the service and maintenance performed. Your shop or technician may require parts, which Travel Buggy has available.

We recommend you inspect your Grand Tourer two to three times per year for in the following areas:

Frame: Tighten loose screws as necessary using the tool kit included with your Grand Tourer

Wheels: Make sure they rotate easily

Significant resistance may require new bearings

Armrests: Look for excessive wear or if they are loose

Replacement parts and more major maintenance issues should be done by a professional. We do not recommend that you do intense types of maintenance work for the Grand Tourer yourself.

LUBRICATION

- According to use condition, please clean or lubricate the front axle, front fork and rotating parts after 180 days
- Molybdenum disulfide grease is recommended
- Your GRAND TOURER Power Wheelchair has electric hub motors that do not require lubrication, please do not lubricate the motors

Question	Possible Causes	Answer
Why is there NO POWER to the Grand Tourer?	<ul style="list-style-type: none">● Controller system power is not connected.● Controller system CPU and joystick are not connected.● Battery power is too low.	<ul style="list-style-type: none">● Connect the Battery.● Ensure all connections between Joystick, CPU and Battery are securely tightened.● Charge the Battery.
Why does my power chair beep and not move?	<ul style="list-style-type: none">● Freewheel mode engaged.● Capacity is reached.● Other error code.	<ul style="list-style-type: none">● Ensure the rear motor switches are in POWER mode, rear wheels should be locked.● Contact Travel Buggy if other error code present.

LIMITED WARRANTY

The purchaser of this product is entitled to a limited warranty as offered by Travel Buggy and its affiliates on the following components and timeframes from the date of purchase:

5 YEAR WARRANTY - Main Frame

2 YEAR WARRANTY - Motor

2 YEAR WARRANTY -Control System

1 YEAR WARRANTY -Battery

90 DAY LIMITED WARRANTY - On wear parts, such as the tires, armrest padding, joystick button covers, joystick knobs, seat cushion, and backrest padding.

This warranty is valid for the replacement of dysfunctional parts only. Any parts under warranty will be replaced and shipped to your door. Any service and labor fees, if applicable, to replace parts under warranty must be paid by the user.

Due to its straightforward design, most parts can be easily exchanged by the end user without a professional service tech required. However, it is always recommended you seek professional help with maintenance and service, to make sure the work is done properly.

The warranty only covers NON-DELIBERATE Damage.

The warranty does not cover:

- Products damaged by user negligence.
- Products were damaged accidentally.
- Products damaged intentionally.
- Products that have been subjected to negligence.
- Products that have been subjected to abuse.
- Products that have been improperly stored.
- Products that have been improperly handled.
- Products that have been improperly operated.
- Products that have experienced general misuse.
- Products that have been modified in an unauthorized, unapproved way.

The warranty is non-transferable and only valid for the original wheelchair purchaser.

Travel Buggy will, at its sole discretion, replace any warrantable item that fails within the warranty period.

CONTACT & SUPPORT

You will find the operation and care of your Grand Tourer very simple and easy, but should you have any questions please contact Travel Buggy.

OUR CUSTOMER HOTLINE

Canada: 1.888.56.BUGGY (28449)

USA: 1.855.56.BUGGY (28449)

EMAIL

contact@travelbuggy.com

Additional images and videos are available on our website.

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