
User manual

Rabbit Up

R82® by Etac

M1062 - 2026-03-30 - Rev. 10.0

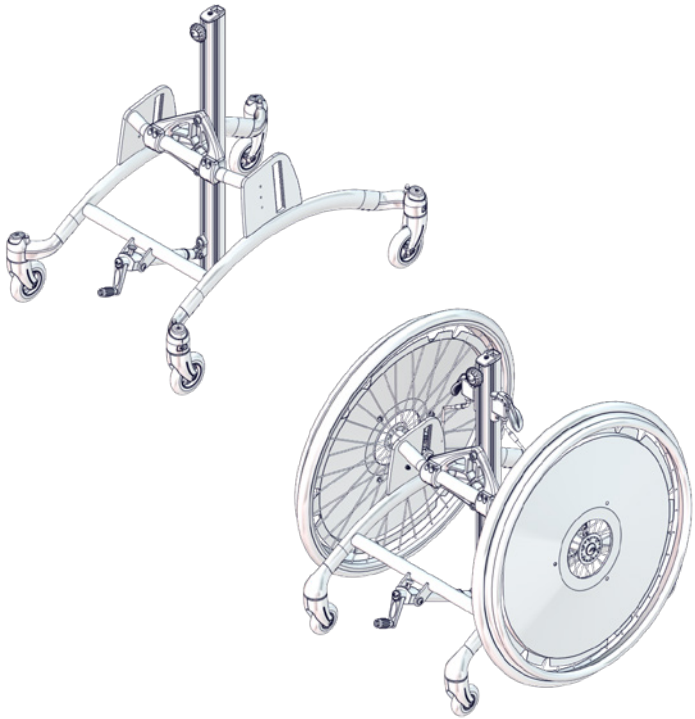www.etac.com

MD CE UKCA

	4
	5
	6
	7
	8
	9
	10
	12
	13
	14
	15
 mm (inch) kg (lb)	20
 mm (inch) kg (lb)	21
	22

ENGLISH	24
DEUTSCH	28
NEDERLANDS	34
DANSK	40
NORSK	44
SVENSKA	48
SUOMI	52
ESPAÑOL	56
FRANÇAIS	62
ITALIANO	68
PORTUGUÊS	74
PORTUGUÊS DO BRASIL	80
РОССИЯ	86
POLSKI	92
ČESKY	98
中國	103
日本	108
ΕΛΛΗΝΙΚΗ	114


1.


RABBIT UP 1	RABBIT UP 2
<40 kg <88.1 lb	<50 kg <110.2 lb

RABBIT UP 3	RABBIT UP 4
<60 kg <132.2 lb	<70 kg <154.3 lb


A.

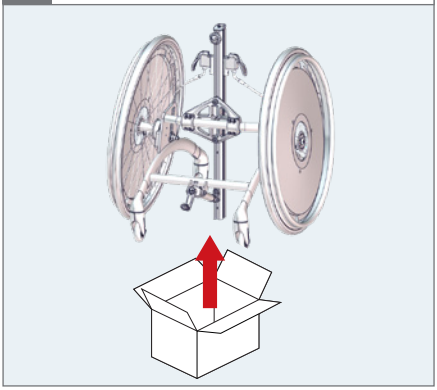
B.

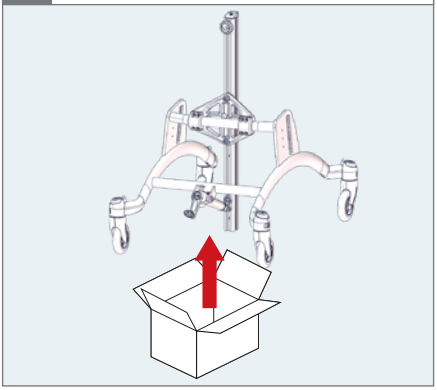
1A.

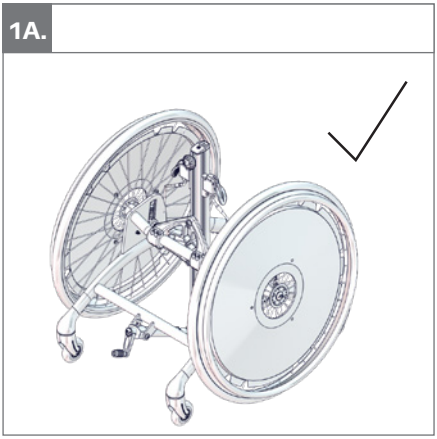
1B.

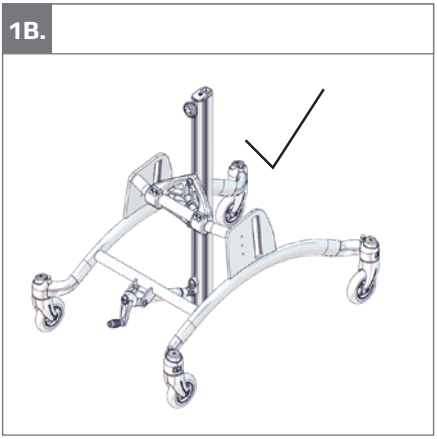
1C.

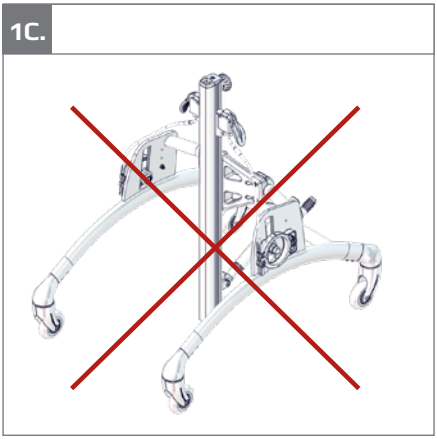
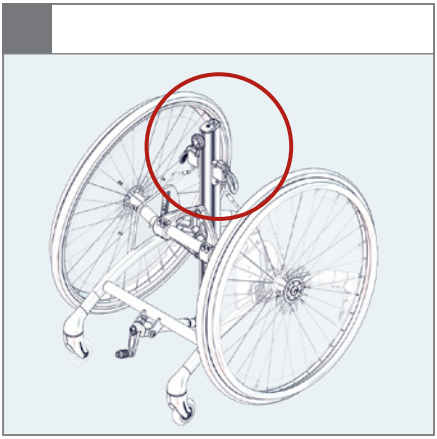
1.

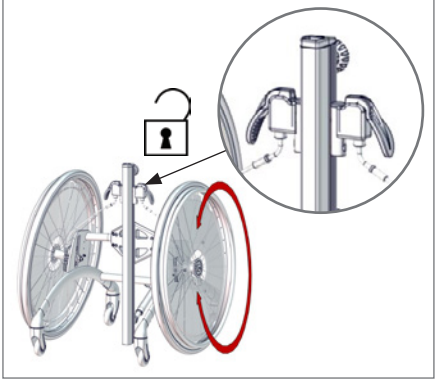
2.

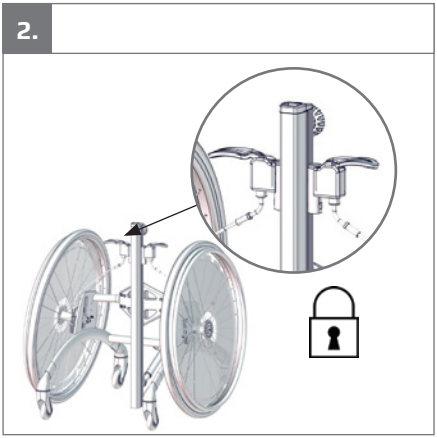
3.

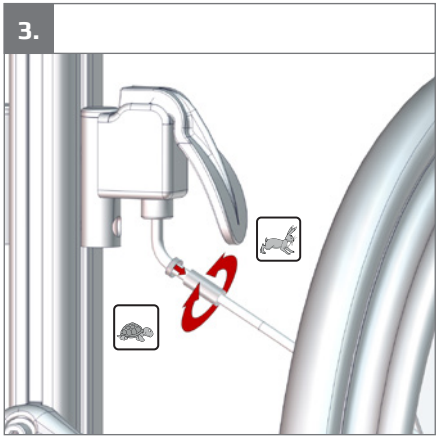
4.

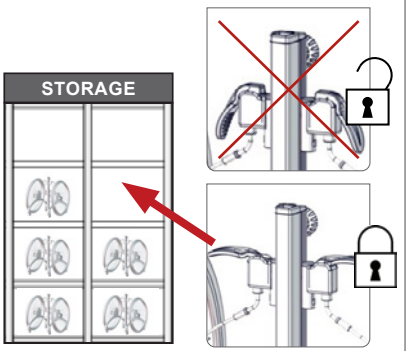
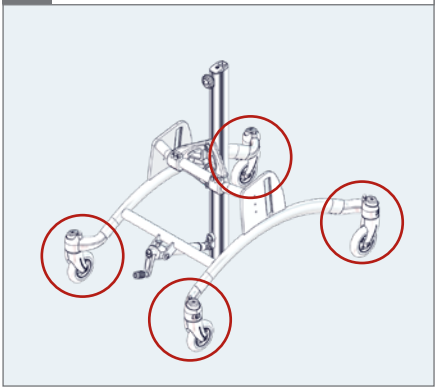
1.

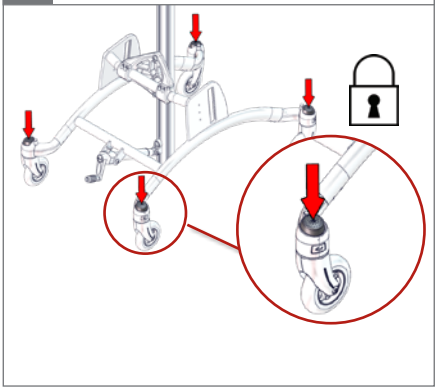
2.

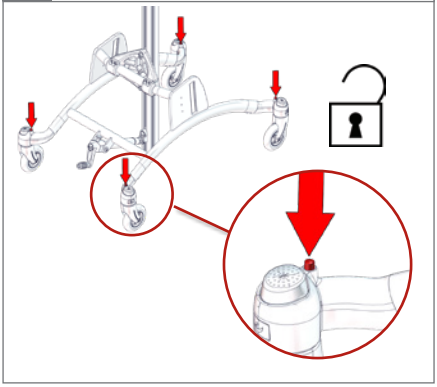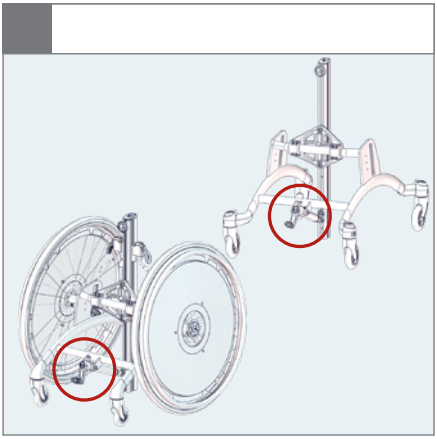
1.

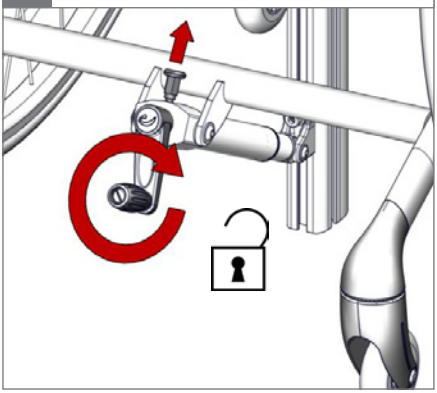
2.


3.

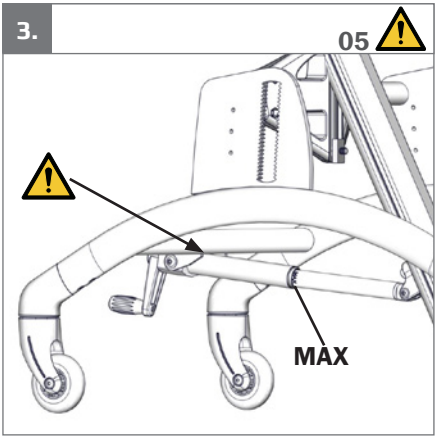
4.

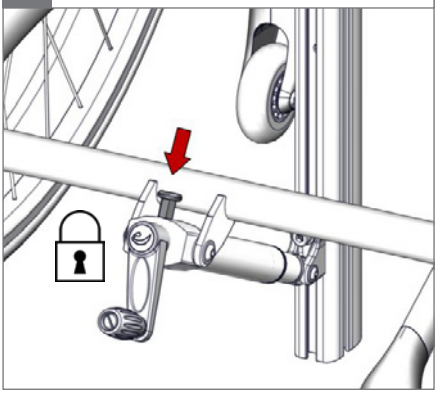
5.

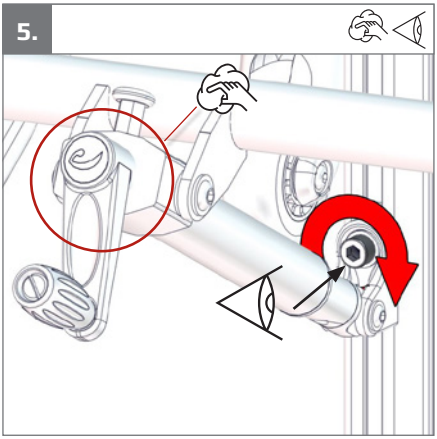
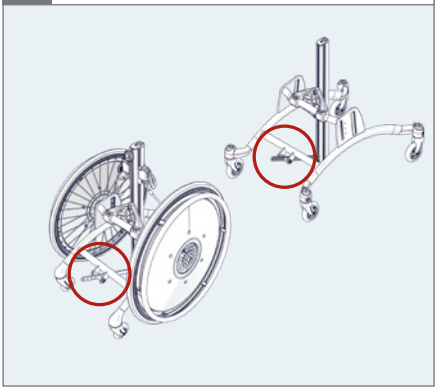
1. ⬡ 6 mm 13 mm

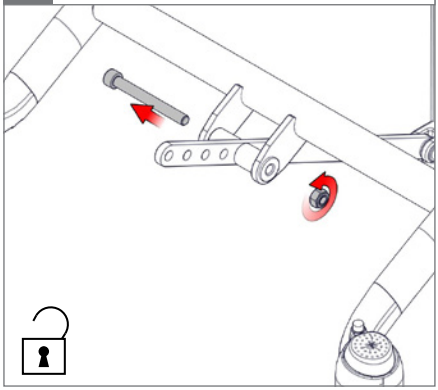
2. 0-20°

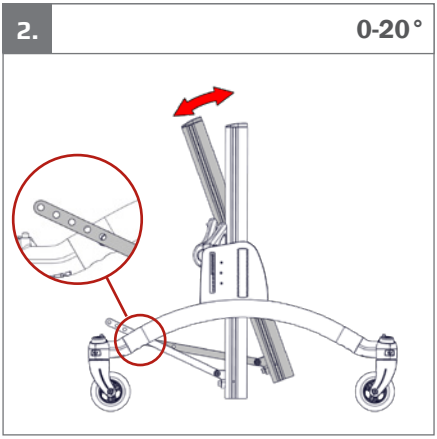
3.

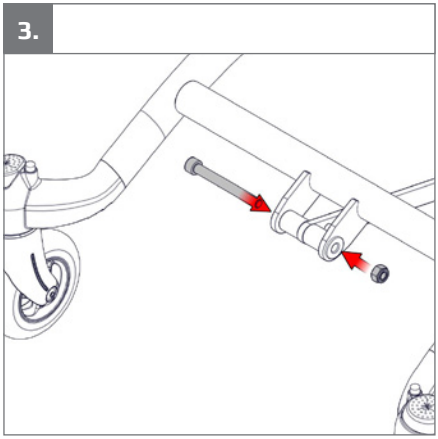
4. ⬡ 6 mm & 13 mm

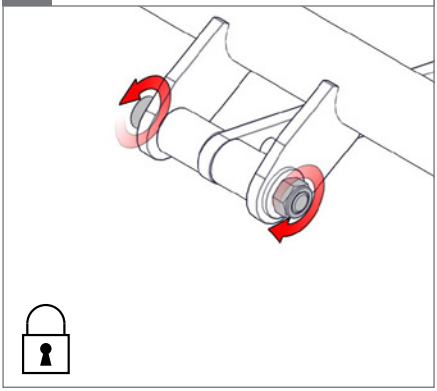
5.

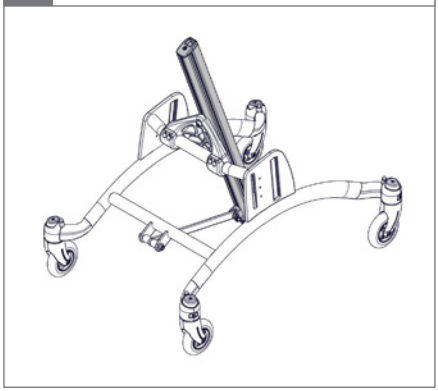
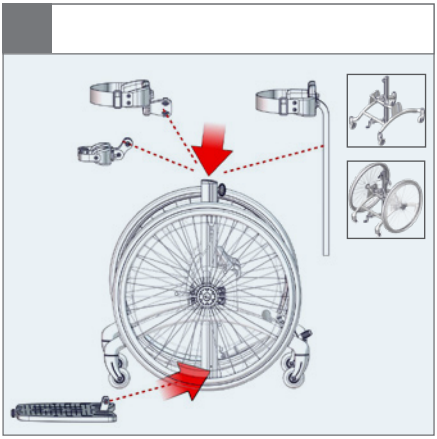
1. ⬡ 5 mm

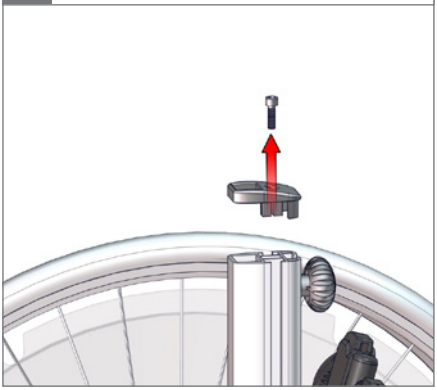
2. ⬡ 5 mm

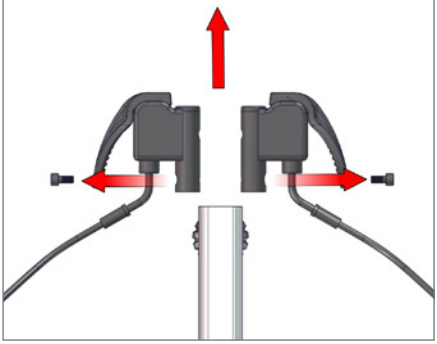
3. ⬡ 6 mm


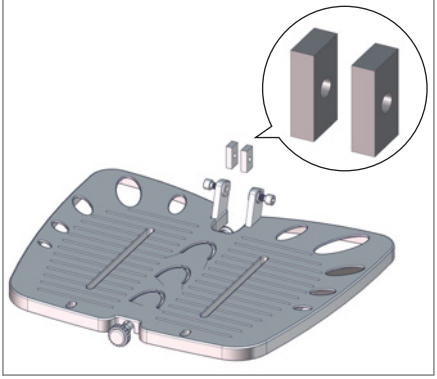
4.


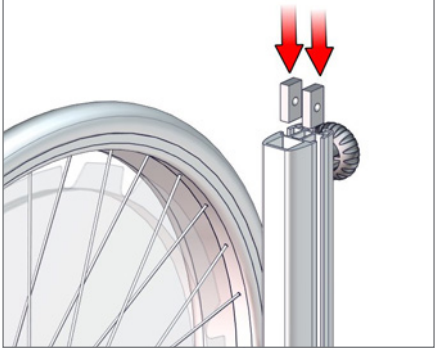
5.

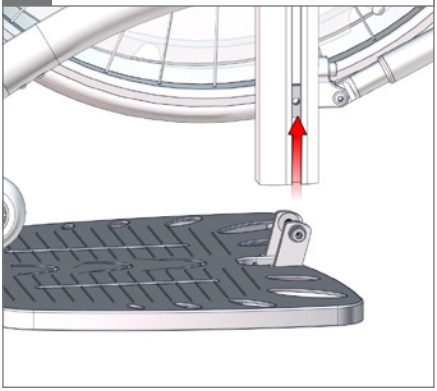
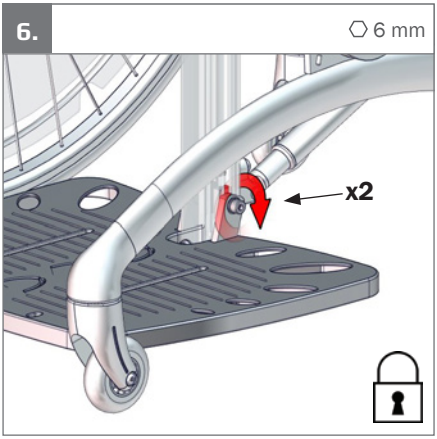
6.
6 mm
x2

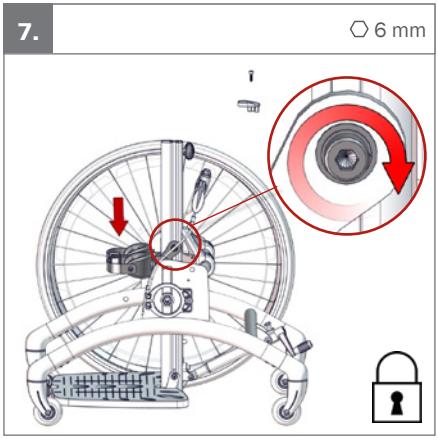
7.
6 mm

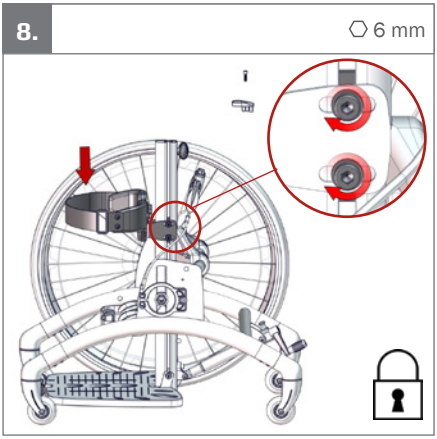
8.
6 mm

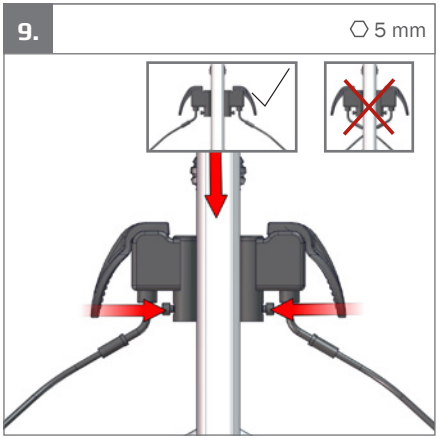
9.
5 mm

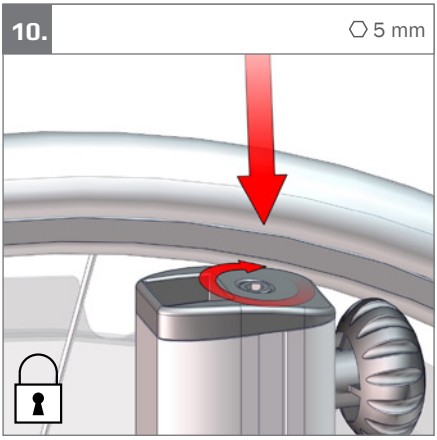
10.
5 mm

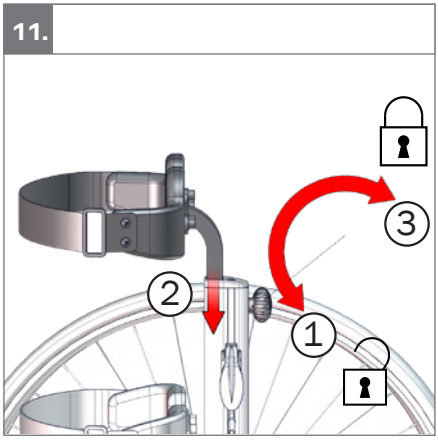
11.
3
2
1

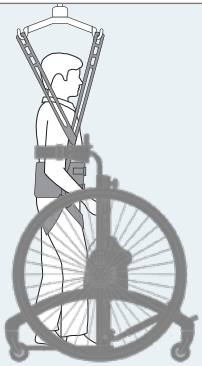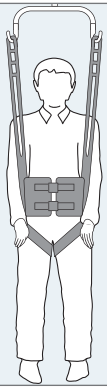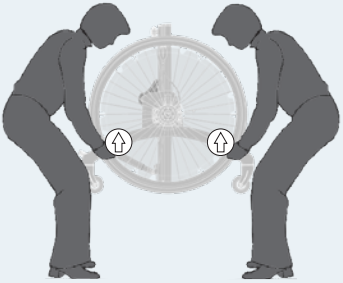
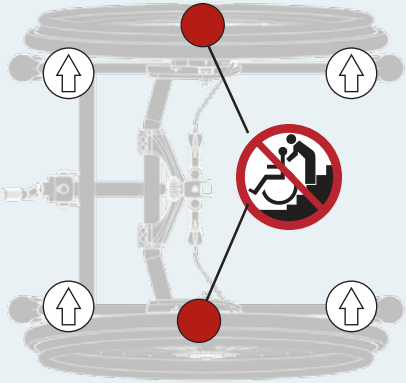
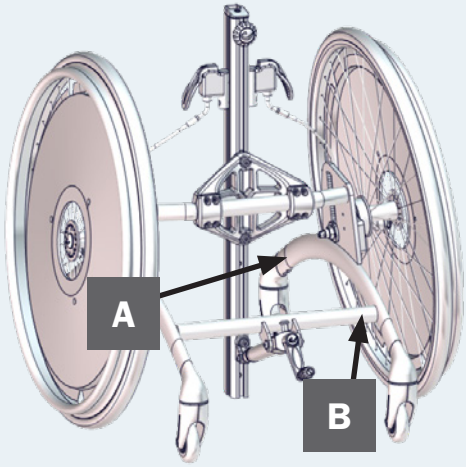
A
B

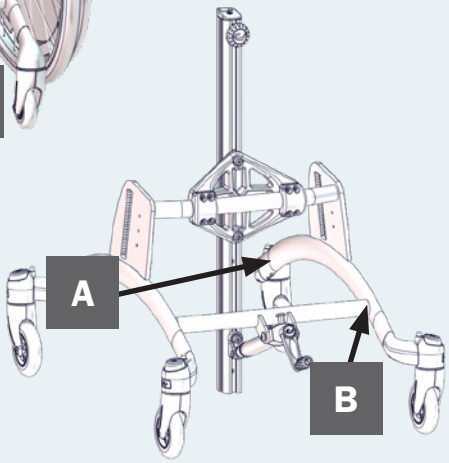
A
B

A

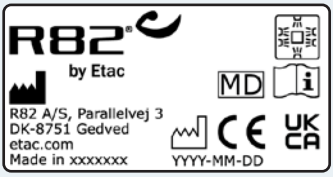R82
by Etac
MD
R82 A/S, Parallelvej 3
DK-8751 Gedved
etac.com
Made in xxxxxxx
YYYY-MM-DD
UK
CA

B

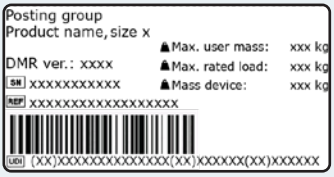Posting group
Product name, size x
Max. user mass: xxx kg
DMR ver.: xxxx
Max. rated load: xxx kg
SN xxxxxxxxxxx
Mass device: xxx kg
REF xxxxxxxxxxxxxxxxx
UDI (XX)XXXXXXXXXXXXX(XX)XXXXXX(XX)XXXXXX

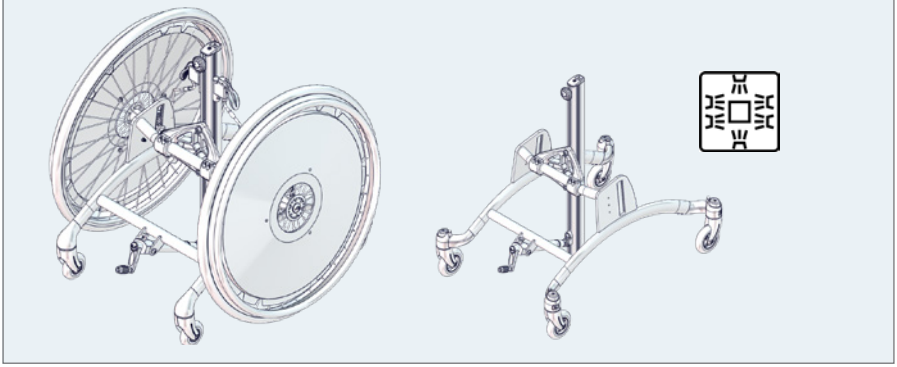EN: Warning: This symbol appears in the User Guide along with a number referring to the instructions below. It draws attention to situations where the product or the user or carer's safety may be in danger
01 = Entrapment Hazards
02 = Tilting Hazards
03 = Transportation in motor vehicles
05 = Squeezing Hazard

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





mm (inch) kg (lb)

	1	2	3	4
 Overall width, with driving wheels	635 mm (25")	680 mm (26¾")	730 mm (28¾")	720 mm (28¼")
 Inside width	325 mm (12¾")	375 mm (14¾")	425 mm (16¾")	425 mm (16¾")
 Length, with driving wheels	620 mm (24½")	750 mm (29½")	800 mm (31½")	900 mm (35½")
 Overall height	670–1000 (26½–39¼")	850–1180 (33½–46½")	990–1310 (39–51½")	1200–1530 (47¼–60¼")
 User height, floor to armpit	500–740 mm (19¾–29¼")	750–1080 mm (29½–42½")	830–1190 mm (33¾–46¾")	960–1420 mm (37¾–56")
 Length, center bar	500 mm (19¾")	700 mm (27½")	800 mm (31½")	1000 mm (39¼")
 Angle, center bar	0°–30°	0°–30°	0°–30°	0°–30°
 Foot plate angle	(15°) – (-10°)	(15°) – (-10°)	(15°) – (-10°)	(15°) – (-10°)
 Floor to foot plate	60 mm (2¼")	60 mm (2¼")	60 mm (2¼")	60 mm (2¼")
 Camber	8°	6°	6°	4°
 Total mass with driving wheels	13.5 kg (29.75 lb)	15.5 kg (34.1 lb)	16 kg (35.25 lb)	18.5 kg (40.75 lb)
 Maximum user mass	40 kg (88.1 lb)	50 kg (110.2 lb)	60 kg (132.2 lb)	70 kg (154.3 lb)
 Max rated load	55.5 kg (122.3 lb)	67.5 kg (148.8 lb)	78 kg (171.9 lb)	90.5 kg (199.5 lb)
 Static stability sideways, uphill & downhill	10°	10°	10°	10°

mm (inch) kg (lb)

	1	2	3	4
 Overall width	515 mm (20¼")	550 mm (21¾")	600 mm (23½")	600 mm (23½")
 Inside width	325 mm (12¾")	375 mm (14¾")	425 mm (16¾")	425 mm (16¾")
 Length	700 mm (27½")	796 mm (31¼")	796 mm (31¼")	926 mm (36½")
 Overall height	670–1000 mm (26½–39¼")	850–1180 mm (33½–46½")	990–1310 mm (39–51½")	1200–1530 mm (47¼–60¼")
 User height, floor to armpit	500–740 mm (19¾–29¼")	750–1080 mm (29½–42½")	830–1190 mm (33¾–46¾")	960–1420 mm (37¾–56")
 Length, center bar	500 mm (19¾")	700 mm (27½")	800 mm (31½")	1000 mm (39¼")
 Angle, center bar	0°–30°	0°–30°	0°–30°	0°–30°
 Foot plate angle	(15°) – (-10°)	(15°) – (-10°)	(15°) – (-10°)	(15°) – (-10°)
 Floor to foot plate	60 mm (2¼")	60 mm (2¼")	60 mm (2¼")	60 mm (2¼")
 Minimum Turning diameter	850 mm (33½")	950 mm (37½")	980 mm (38½")	1085 mm (42¾")
 Mass device	9 kg (19.8 lb)	10.5 kg (23.1 lb)	11 kg (24.2 lb)	13.5 kg (29.3 lb)
 Maximum user mass	40 kg (88.1 lb)	50 kg (110.2 lb)	60 kg (132.2 lb)	70 kg (154.3 lb)
 Max rated load, no driving wheels	51 kg (112.4 lb)	62.5 kg (137.7 lb)	73 kg (160.9 lb)	85.5 kg (188.4 lb)
 Static stability downhill & uphill	10°	10°	10°	10°
 Static stability sidewards	10°	10°	10°	10°

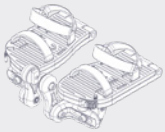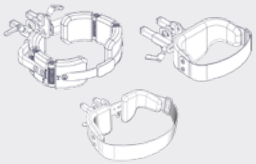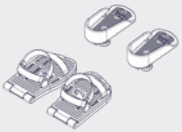	9996097635
	9996097636
	9996097637
	9996097638
	9996097639
	9996097642
	9996097643
	9996097645

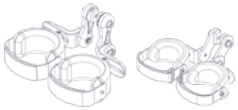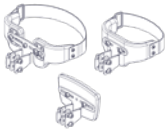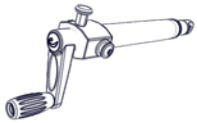	9996097646
	9996097647
	9996097248
	9996097649
	9996097651
	9996097652
	9996097656
	9996097655



9996097810

9996098152

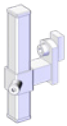9996098153

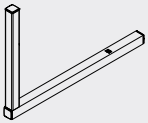9996098211

ENGLISH

General

We appreciate your choice of a new device from R82 - a worldwide supplier of technical aids and appliances for children and teenagers with special needs. To fully benefit from the options offered by this device please read this user manual before use and save it for future reference.

Intended purpose

The family of R82 Standing devices consist of assistive devices intended for the alleviation of or compensation for a functional impairment due to an injury or disability. The devices are designed for an individual lacking the ability to maintain the standing position without the use of an assistive device.

Intended user group

Selection of an appropriate device size is dependent on body weight and physical dimensions of the user. Secondary users of the device are caregivers providing assistance to the primary user as well as clinicians/ technicians who set up the device.

The device is intended for children and adults with a need for postural support to enable the standing position due to impaired mobility and functional ability and not on a specific health condition or age.

Intended environment

The device is intended for indoor use, in home environments or institutions and outdoor on even level and hard ground. It is not intended for wet areas.

Intended application

The user is the person standing in the device. The caregiver is the person operating the device.

The cumulative sum of time of use (exposure) will be either prolonged (more than 24 hours but less than 30 days) or long-term (exceeding 30 days). The use of the device will be discontinued multiple times during the day when the user will be lying down, sitting in another assistive device or a conventional product.

The device is intended for use with intact skin only. Depending on the configuration and the accessories used, the contact will often be with the skin of the head and neck region, shoulder region, upper extremities and lower extremities of the user and caregiver.

The device is intended for refurbishment and reuse.

Indications

The assistive devices in the R82 Standing Family are indicated for children and adults who need postural support to enable supported standing due to a persistent, significant, or substantial impairment of mobility, i.e., disability. The indications for using the devices in the R82 Standing Family are based on functional ability and not a specific health condition or age.

Contraindications

There are no known contraindications.

Precautions

Precautions for use are based on known and possible risks as identified in the R82 Standing devices User Manuals and detailed in the R82 Standing device Risk Management as follows.

Cautions:

Special attention should be given to the considerations about the configuration of the device and the use of accessories for users with limitation in their:

- Sensory functions and the sensation of pain, consider extra checks of the pressure caused by the support provided by the device and its accessories.
- Movement function (e.g. experience involuntary movements).
- Ability to maintain a standing position, use the device with accessories (e.g. chest or back support).
- Ability to maintain the position of their head, consider using accessories or to use a device that provide support for positioning of the head.

Enclosed tools:

- 5 mm Allen key
- 6 mm Allen key
- 13 mm spanner

CE Declaration of Conformity

This device conforms to the requirements of the Medical Devices Regulation (2017/745).

The CE mark must be removed when:

- rebuilding the device
- using the device outside the intended purpose and declaration of conformity
- using other than original R82 spare parts and fittings.

Combinations

When combining the device with other devices, the combination must not change the intended purpose of the devices or modify the devices in such a way that compliance with the applicable requirements may be affected. The person / company responsible for the combination of the devices must ensure that the applicable requirements are fulfilled.

Applicable requirements examples:

- Material safety (biocompatibility, flammability), mechanical safety (stability, squeezing hazards, inter connection), electrical safety.

Disposal

When a device reaches the end of its lifetime, it should be separated by material type so that the parts can be recycled or disposed of properly. If required, please contact your local dealer for a precise description of each material. Please contact your local authorities to be informed of the options available for recycling in an environmentally friendly way.

Accessories and spare parts

The devices from R82 can be supplied with a variety of accessories which comply with the needs of the individual user. Spare parts can be ordered on request. Find the specific accessories on our website or contact your local dealer for further information.

R82 Warranty

www.etac.com/support

Presales information and instructions

Presales information and instructions are available at www.etac.com, under your device, in the documents section.

Safety

Signs, symbols and instructions placed on this product must never be covered or removed and must remain present and clearly legible throughout the entire lifetime of the product. Immediately replace or repair illegible or damaged signs, symbols and instructions. Contact your local dealer for instructions.

Any adverse event in relation to the device must be reported to your local dealer and the national competent authority in a timely manner. The local dealer will forward information to manufacturer.

Caregiver

- Read the user manual carefully before use and keep it for future reference. Incorrect use of this product may cause serious injury to the user.
- This product contains small parts which could present a choking hazard if removed from their intended placing.
- Settings and adjustments to the product or accessories are only to be made by persons who have read the user manual and who has been given the necessary instructions by a product advisor.
- Users who have difficulties reading the user manual can have someone read the user manual aloud to understand how to use the product.
- Always use the correct lifting techniques and aids for this purpose.
- Never leave the user unattended in this product. Ensure permanent supervision by a caregiver.
- Repair/replacements must only be made using new original R82 spare parts and fittings and carried out in accordance with the guidelines and service intervals prescribed by the supplier.
- If there is any doubt as to the continued safe use of your R82 product or if any parts should fail, stop using the product immediately and contact your local dealer as soon as possible.
- Ensure that moving parts does not squeeze or trap body parts when adjusting them.
- Ensure that body parts are not trapped in openings for mounting of supports/ accessories.

Environment

- Observe the surface temperature on the product, before placing an user in it. This goes especially for users with insensitive skin as they cannot feel the heat. If the surface temperature is above 41 °C then allow the product to cool before use.
- When using the device on inclines, please refer to the measurements table for allowed inclines.

User

- If the user is close to the maximum user mass and/or has a lot of involuntary movements, e.g. rocking, a larger size with a higher maximum user mass or a different R82 device should be considered.

Product

- Carry out all positional adjustments on the product and accessories and ensure that all knobs, screws and buckles are securely fastened before use. Keep all tools out of the reach of children.
- As a minimum, use the product with a support placed at the top of the pelvis or higher.
- Activate the brakes before placing the user into the product.
- Ensure the stability of the product is maintained before placing the user in it.
- Ensure the wheels are fully-functional and secure before use.
- When changing wheels or camber parts always refer to the mounting guide 9996097649 for correct placement and configuration. The mounting guide is always available at www.etac.com.
- Inspect the product and all of its accessories and replace any worn out parts before use.
- This device is tested for resistance to ignition at level non-ignition in accordance with applicable standards.
- The product might loose its flame resistant characteristics if using other than original upholstery and/or accessories.

Transportation

- Transporting the user across doorsteps or other small obstacles with the user safely strapped into the device is allowed.
- Transporting the user in the device is only allowed with the device in its transport position.
- The transport position is with the device in an upright position.
- Transport direction: back and forth. Do not transport the product sideways.
- This product is not to be occupied in a motor vehicle. The user must be positioned in a vehicle seat and the product must be stored in a luggage compartment.
- The user must not occupy the product while it is being transported into or out of a vehicle by a lift.

Service & maintenance information

Service and maintenance of a medical device is wholly the responsibility of the owner of that device.

Failure to service and maintain a device in accordance with the instructions may invalidate the device's warranty. Furthermore, failure to service and maintain a device may compromise the clinical condition or safety of users and/or their caregivers. Do not perform service and maintenance while the user is in the product. Contact the local dealer for assistance, if needed, in setting up, using or how to perform service and maintenance of the device.

Lifetime

The lifetime of this product in normal use is 8 years if all maintenance and servicing is carried out in accordance with the manufacturer's instructions and demonstrably recorded.

Reuse

This product and its accessories/components are suitable for reuse.

Refurbishment for reuse

Before reuse or changing the user, follow the service checklist, maintenance and washing instructions and carry out regular inspections.

Service interval

The frequency of inspection may be altered depending upon use and wear. It is recommended that the product is inspected yearly, every time the product is re-issued for use and after a long-term storage (more than 4 months). The inspection should be carried out by a person who understands the use of the product.

Service checklist:

Check and adjust the following as per their intended use:

- Wheels are moving freely and that all wheel locks/brakes are easily activated. Adjust if required
- Spoke guards is without cracks and secured to the wheel with the clips mounted on each guard
- Castors are moving freely and all wheel locks/brakes are easily activated
- The column angle is adjustable and it is posible to lock it in position
- Nuts and bolts on the product are tightened
- Oil the swing-away parts
- Check all fittings and buckles for signs of worn out parts
- Product for cracks or signs for worn-out parts
- All labels on the product are intact

Repair or replace any damaged or worn out parts.

Maintenance

Before every use it is recommended to:

Wipe off any residue and dirt from the product, using a cloth with warm water and a mild detergent/soap without chlorine and let it dry before use.

Carry out a visual inspection, for any damaged or worn-out parts.

Washing

This device can be washed by hand or in a machine depending on configuration.

Machine wash

Dismount all accessories and wash them separately.

This device can be washed at 60°C using a mild detergent for 10 minutes in a washing machine designed for medical devices. Use the machine's drying function to dry the device.

Tilt/rotate the device after washing to facilitate water drainage from within the device.

Manual wash

Dismount all accessories and wash them separately.

This device and its accessories can be washed manually. Use a mild detergent/soap and let all parts dry before use.

Fabric and covers

The cover can be washed at max. 60°C using a mild detergent in a washing machine.

Disinfection

The product can be disinfected with a 70% disinfectant IPA solution. It is recommended to wipe off any residue and dirt from the product, using a cloth with warm water and a mild detergent/soap without chlorine and let it dry before disinfection.

Materials

- Aluminium
- Steel
- Plastic
- Rubber
- Stainless steel
- Polyurethan
- Wood

Surface treatment

The following surface treatment processes were used to protect against corrosion:

- Lacquered surfaces = Polyester powder coating or ED-coating
- Non-lacquered aluminium parts = Anodized coating
- Non-lacquered steel surfaces = Galvanized

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Etac Ltd.
Gelbart AG