




Adjustable bath safety step  
assembly instruction


**7100**

048503710002

Marche réglable pour baignoire  
instructions d'assemblage



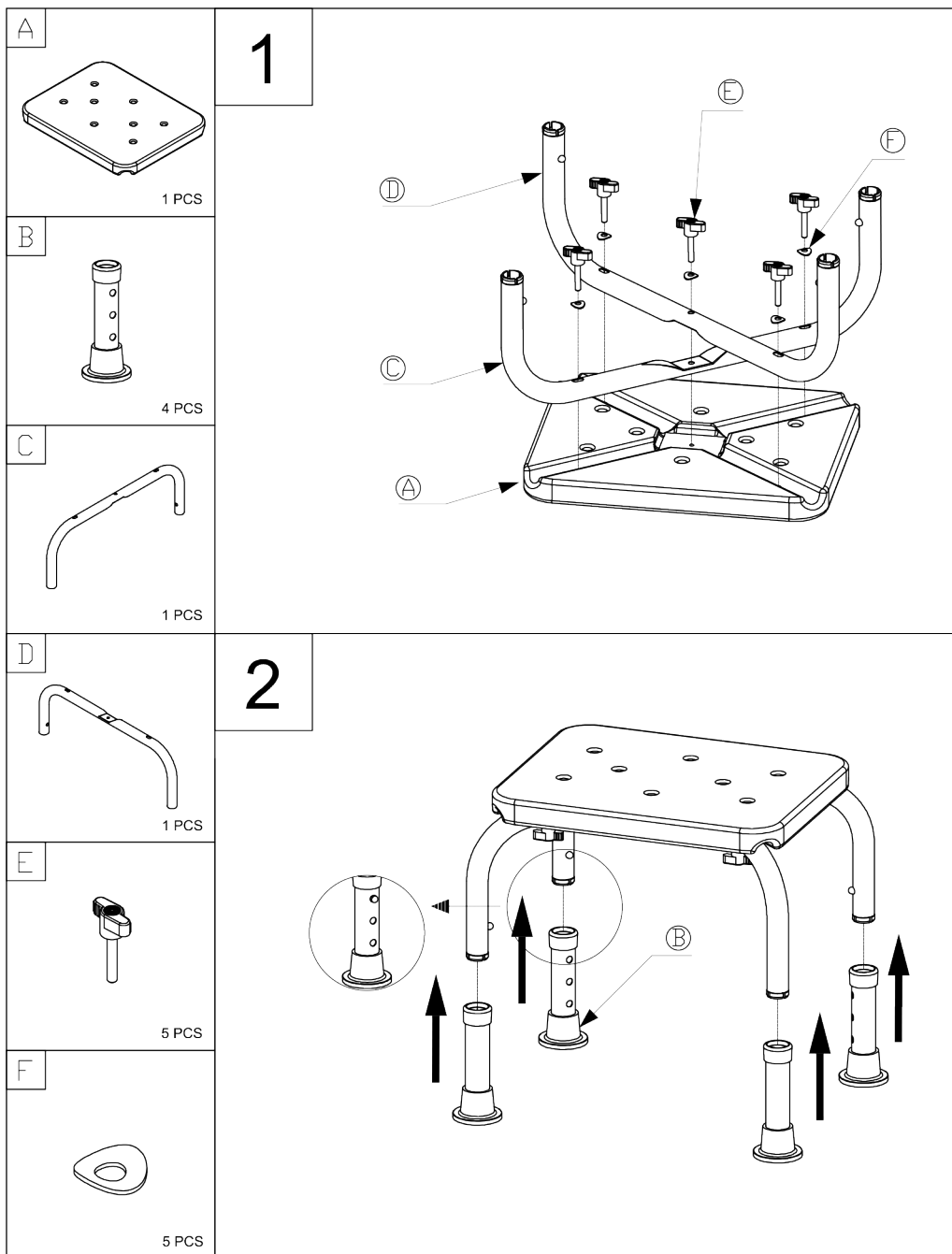

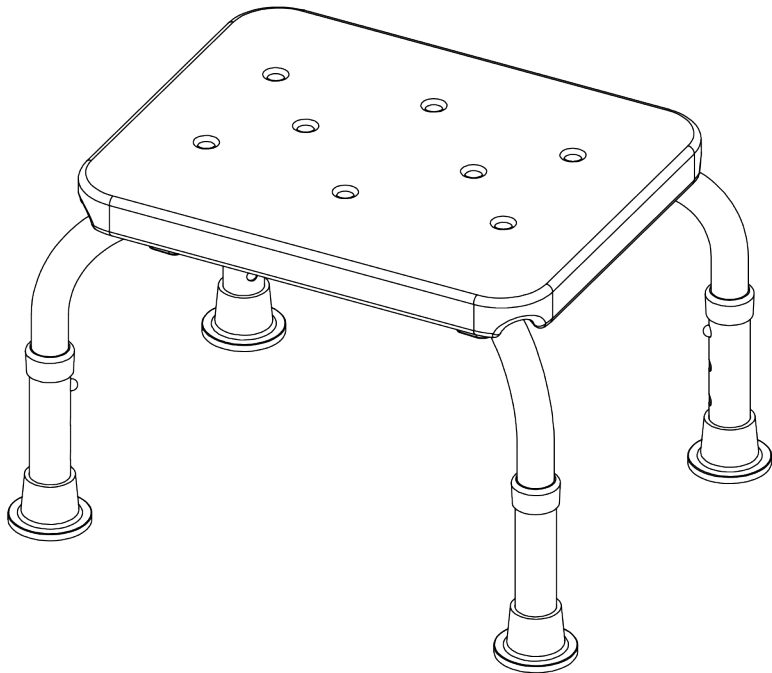

- 1) Place padded step piece face down on a flat surface
- 2) Position top and bottom frame pieces over step
- 3) Fasten frame pieces to seat by attaching washer and screw through center screw hole
- 4) Then fasten four remaining washers and screws
- 4) Fasten all loosely at first, then tighten once all screw holes are engaged
- 5) Insert leg extensions into ends of frame assembly

- 1) Placez la marche matelassée orientée vers le bas sur une surface plane
- 2) Positionnez les parties supérieure et inférieure du cadre sur la marche
- 3) Fixez les parties du cadre au siège en passant une rondelle et une vis par le trou de vis central
- 4) Fixez alors les quatre rondelles et vis restantes
- 5) Fixez-les toutes sans serrer au début, puis serrez davantage une fois que tous les trous de vis sont en prise
- 6) Insérez les extensions des pattes dans les extrémités du cadre

# PCP<sup>®</sup>

## Personal Care Products



PCP, division of SAI,  
 Ripley, OH 45167