

AVERO Premium plus



beka DIFFERENT
BY DESIGN




Watch product
video!



Functionality meets wellness.

The AVERO Premium plus lifting bathtub was developed in close cooperation with nurses and occupational therapists. The ergonomic tub body with the lowered foot area and high-quality materials enable a relaxation bath with a wellness character. The tub body and the foot construction allow the use of all modern care lifts on the market.

The AVERO Premium plus stands freely in the room without anchoring to the floor and is therefore freely accessible from 3 sides. A control panel with symbols prevents incorrect operation and simplifies the work of the nursing staff. Available in different bath lengths as model 175 or model 190.

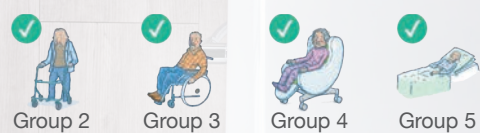

Group 2

Group 3

Group 4

Group 5




Operating panel with temperature display



Accessible from 3 sides with all modern lifters



Great arm and shoulder room

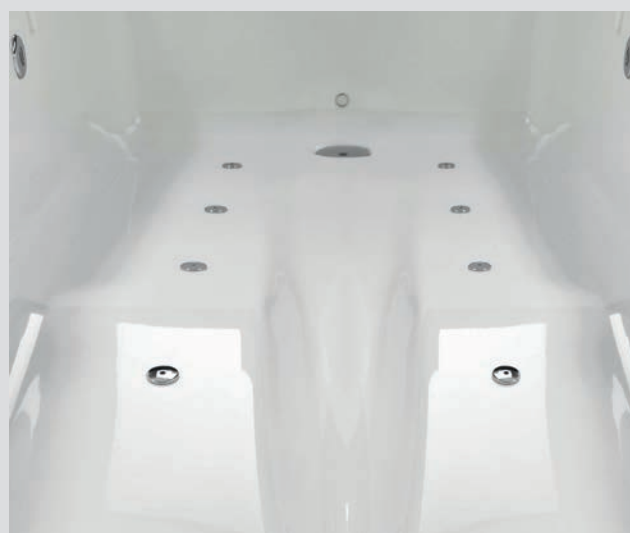

The lowered foot area supports ergonomic lying and prevents slipping



Large armrests with high quality handles offer safety and comfort to bathers



Optional: integrated disinfection system



Lowered foot area & special back profile

The lowered leg area offers bathers reliable support. Together with the distinctive armrests and the molded back profile, the components form the basis for the ergonomic tub body. This ensures a relaxed and safe bathing process with sufficient hold, especially when used with seat lifts.



Intuitive operating and control unit

All functions of the bath can be easily controlled using the control panel integrated in the console. The water temperature of the bath spout, bath water and hand shower is always displayed.



Wide tub body offers plenty of space

Residents of all sizes enjoy sufficient freedom of movement inside the body of the AVERO Premium plus.



Technical Data

Total dimensions* (L x W)	2.335 mm x 850 mm
Inner tub body dimensions* (L x W x H)	1.780 mm x 575/720 mm x 430/490 mm
Outer tub body dimensions* (L x W x H)	1.934 mm x 800 mm x 520 mm
Height adjustment	690 mm - 1.040 mm
Capacity*	ca. 220 litres
SWL	300 kg
Power connection	230 V, 50 Hz (120 V, 60 Hz)

* Specifications for AVERO Premium plus, model 190

Product Information

Electrical Height Adjustment
Easy operation via control panel with intuitive symbols
Safety modules prevent incorrect operation of the tub
Temperature display for bath spout, hand shower and bathing water
Ergonomic inner tub body
Lowered foot area
Armrests with handles for more comfort and safety
Seat with hygiene opening
Scalding protection for hand shower and bath spout

Available Options

-  Aroma system
-  Colored-light system
-  AirSpa system
-  Music-sound system
-  Bluetooth-music-sound system
-  RAL colors available
-  Disinfection system
-  Automatic tub filling
-  Emergency function in case of power failure
-  Bathtub shortener



BEKA Hospitec GmbH | Am Rübenmorgen 3 | 35582 Wetzlar-Dutenhofen
Tel +49 641 922 20-0 | E-Mail info@beka-hospitec.com
www.beka-hospitec.com

