
HANDICARE 1100 | SERVICE MANUAL

MANUAL DE SERVICIO
KUNDENDIENSTANLEITUNG
SERVICEHANDLEIDING
MANUAL DE SERVICE
MANUALE PER SERVIZIO

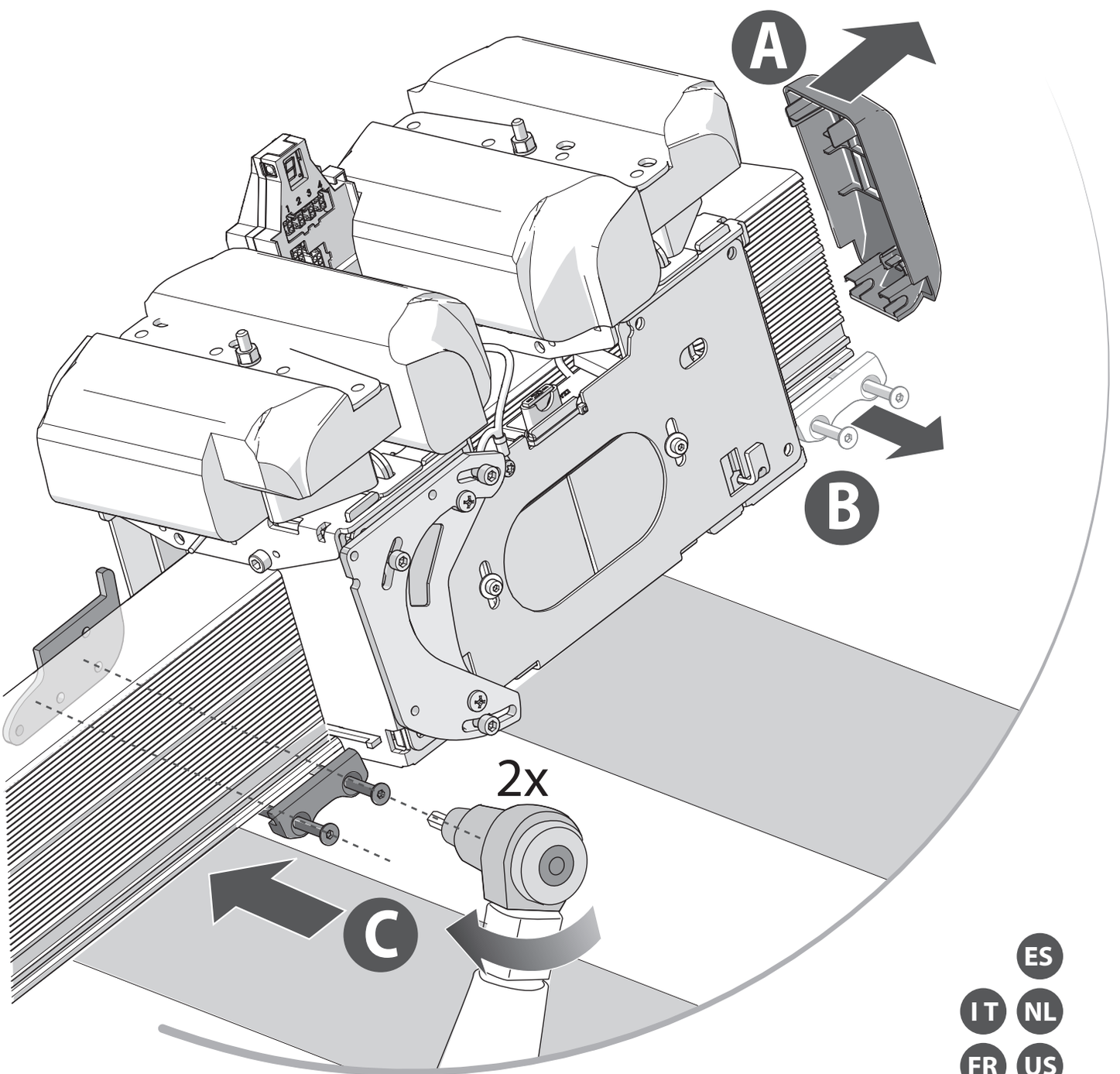ES
IT NL
FR US
DE EN

handicare

AA17833, Revision 21

EN *All rights reserved. No part of this publication may be reproduced and/or published by means of print, photocopy, microfilm or in any other way or form, without prior written permission from the manufacturer. This also applies for the accompanying drawings and diagrams. The manufacturer reserves the right to change parts at any time, without prior or direct notification to the customer. The contents of this manual may also be changed without prior warning. This manual is intended for the standard models of the lift. The manufacturer therefore cannot be held liable for any damage resulting from specifications that differ from those of the standard model of the lift supplied to you.*
This manual has been compiled with the best intentions and every possible care but the manufacturer cannot accept responsibility for any errors in this manual or for the consequences thereof.

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EN
The Dutch manual is the original language

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
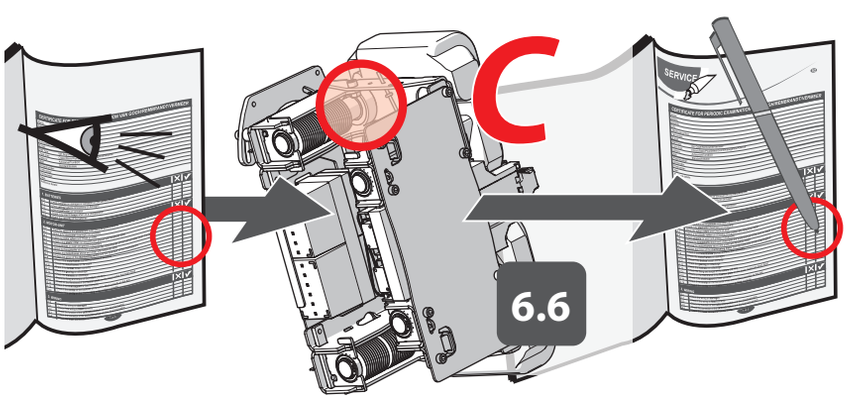
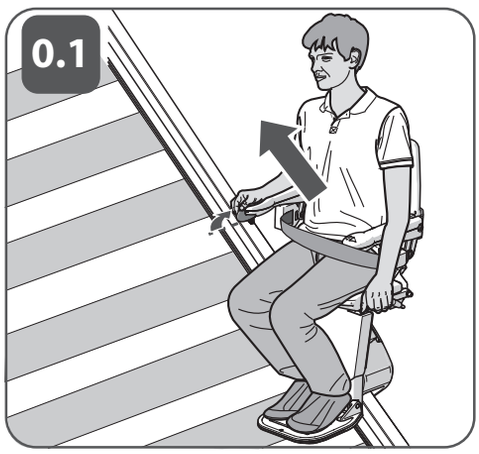
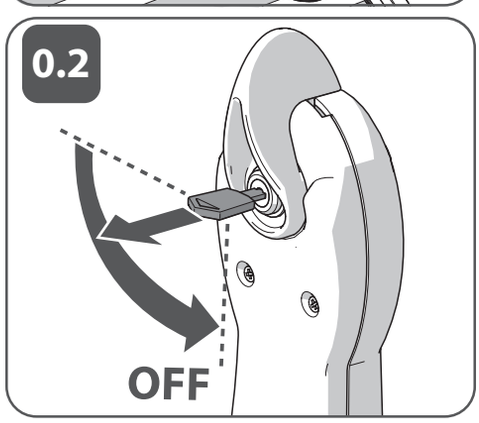
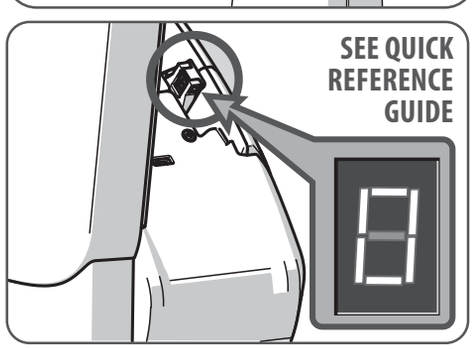
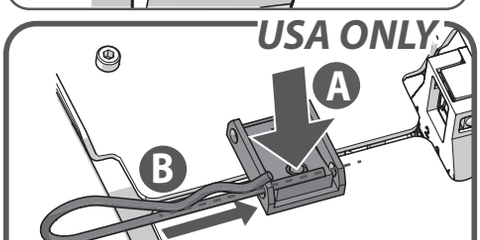
Turn off and lock the switch, before the service work is started

When installing a new part on a existing lift always replace the part by a part with the same article number.

TOOLS

REQUIRED TOOLS & DOCUMENTS

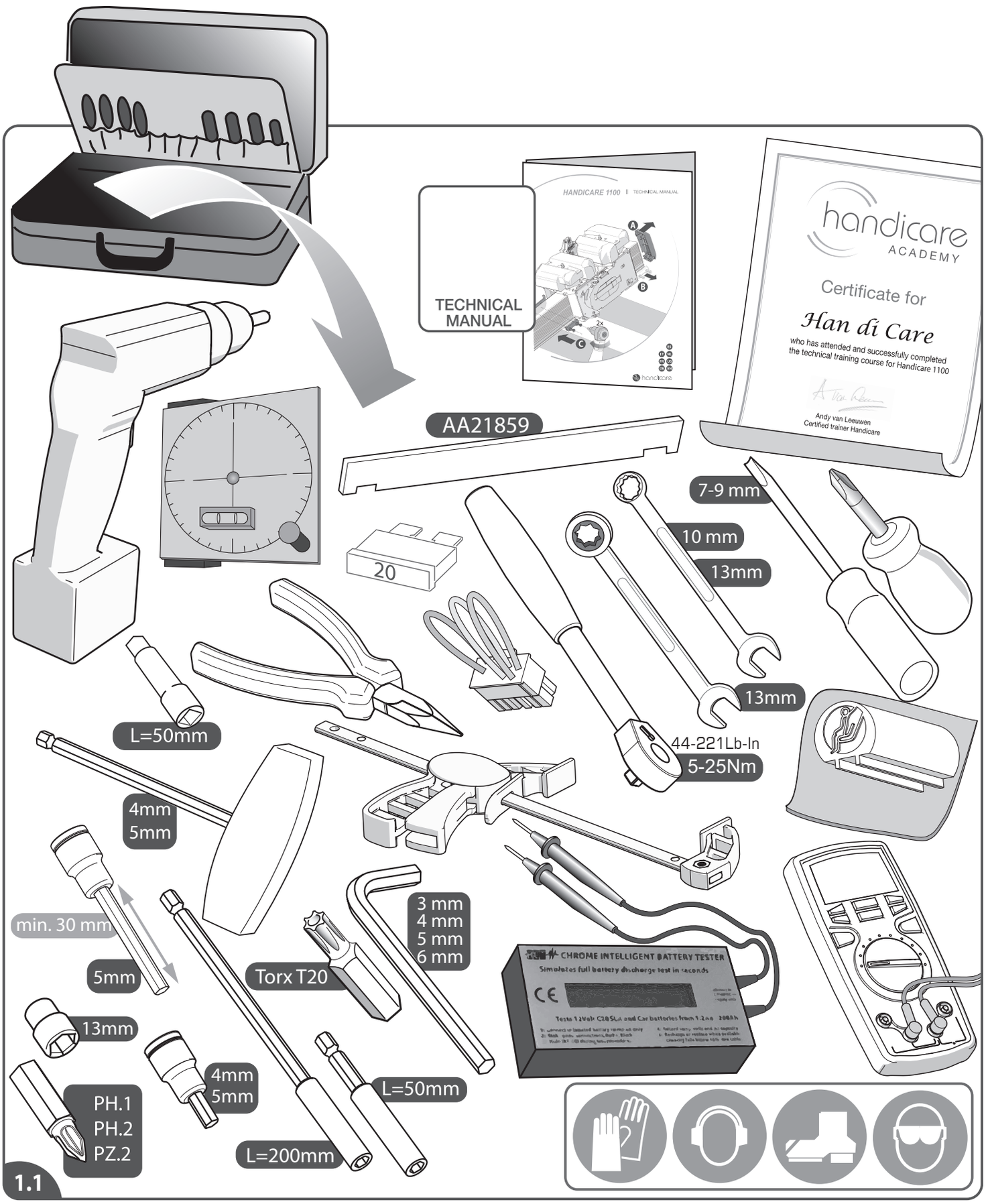EN

Only a technician certified for this lift should install this lift (for more detailled information about the required competence see attachment G of the installation manual**). The latest instructions for the lift and instructions for options are available for download at www.handicarepartners.com. Routine, Periodic, and Acceptance inspections and test in the USA shall be done according to ASME A18.1.**

Ensure that tools are sufficient and in good condition. The minimum tools required in addition to the Installation kit are shown below:

Allen keys {3mm, 4mm, 5mm}	Shorting connectors	Torque spanner
Sockets {10mm,13mm}	Pliers	T-bar Allen keys {4mm, 5mm}
Drill	Spanners {10mm, 13mm}	Battery Tester
Screwdriver with fine flat head	Torx bit {T20}	Installation kit
Phillips bit {PH2, Pz2}	Bit extenders {200mm, 50mm}	Spirit level

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




1 SERVICE LIFT
RAIL

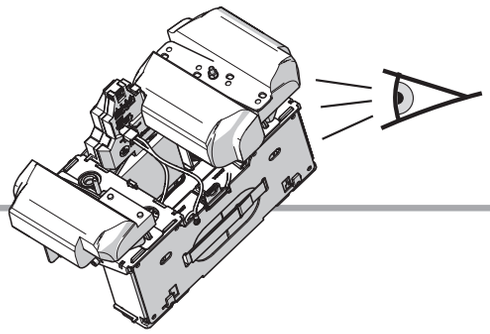
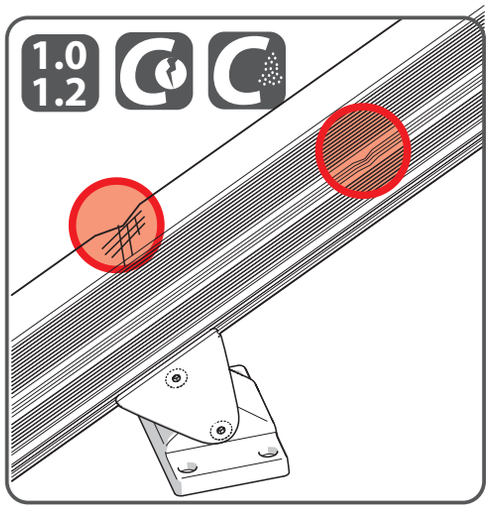
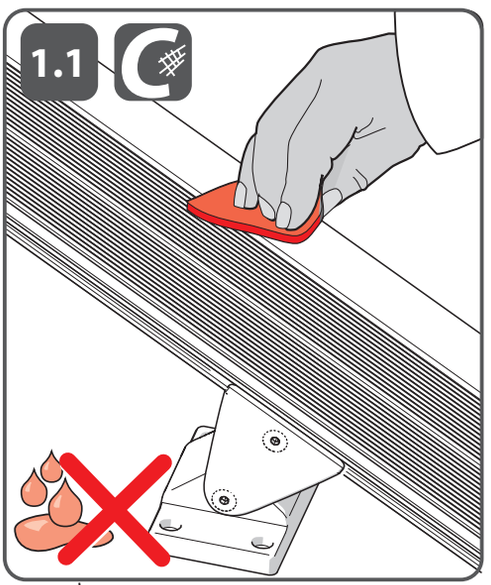
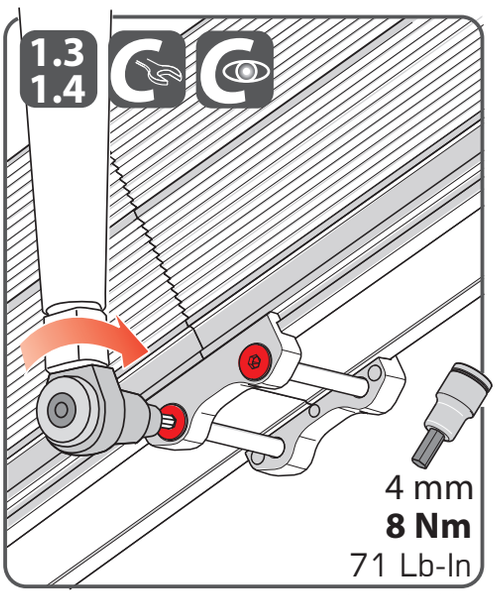
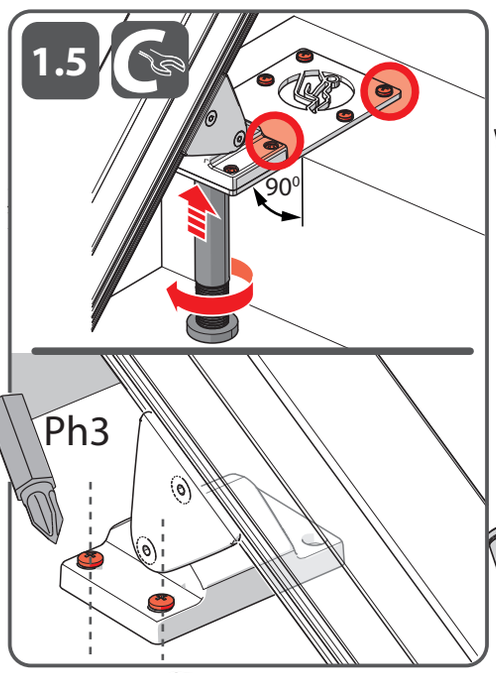
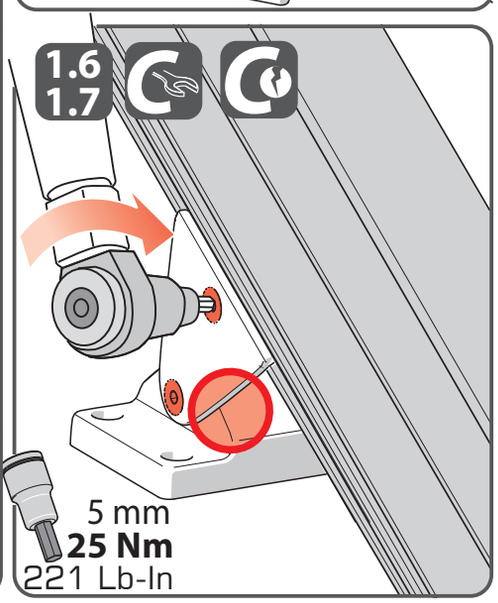
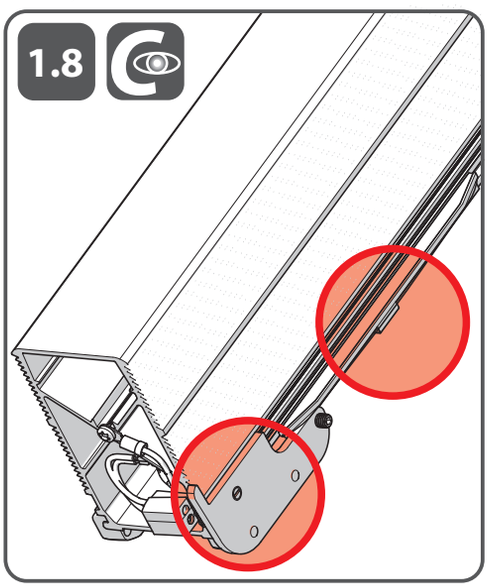
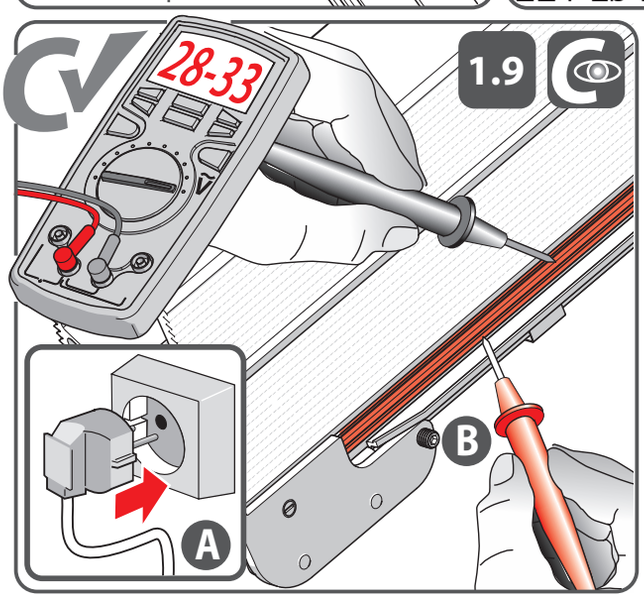
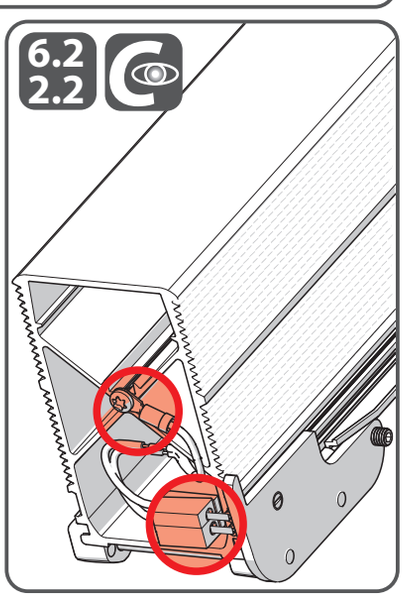
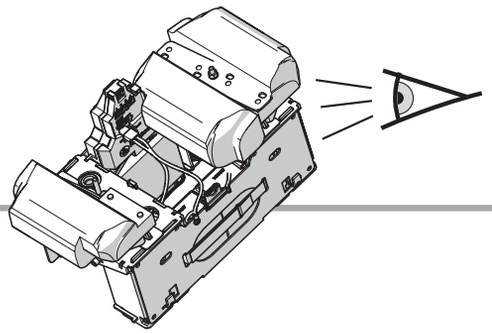
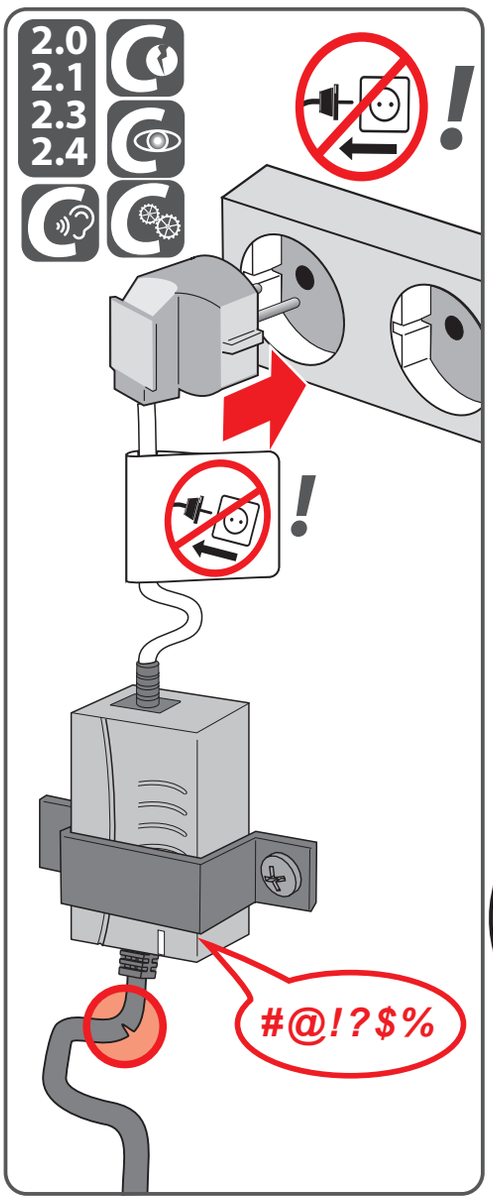
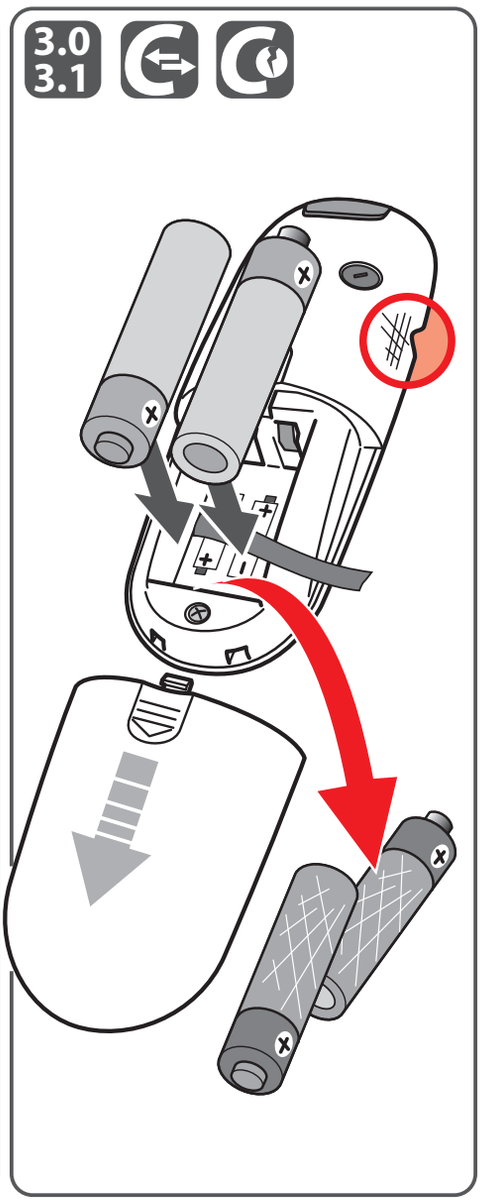
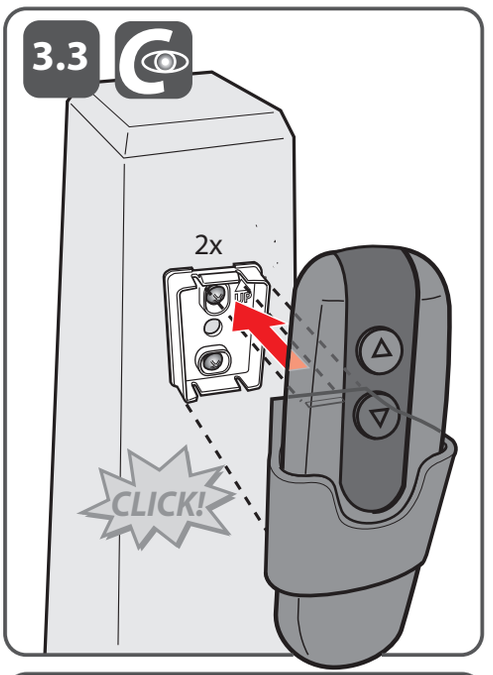
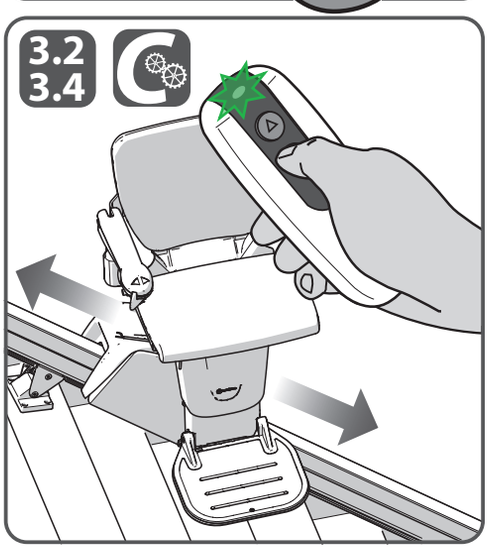
FOR EXAMINATION OF THE NEXT SECTION THE CHAIR AND THE COVERS OF THE UNIT SHOULD BE REMOVED AND THE UNIT SHOULD BE TAKEN OF THE RAIL (SEE ATTACHMENT A)

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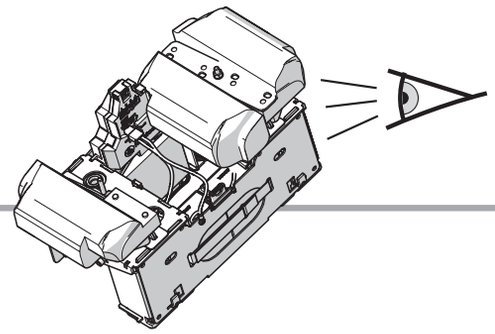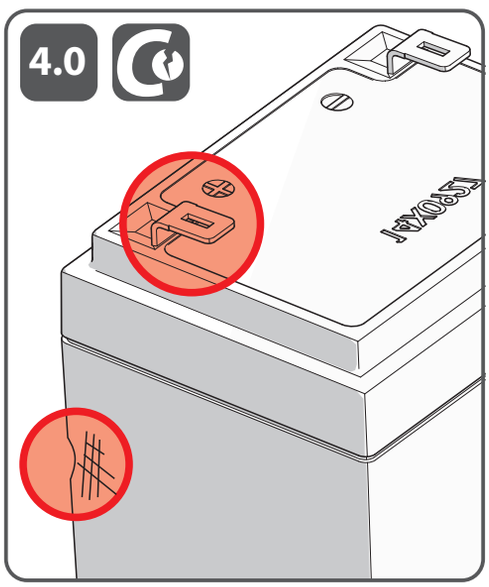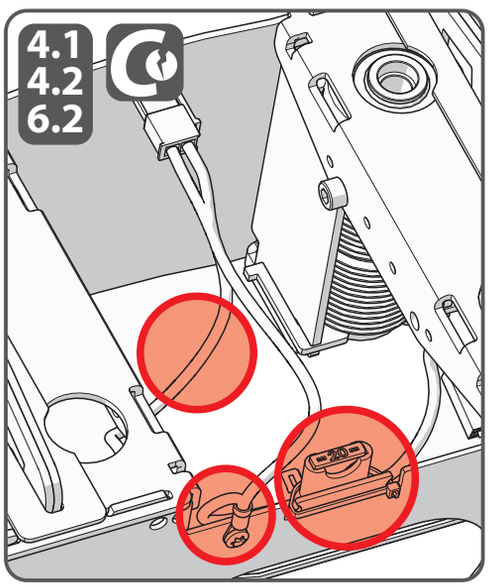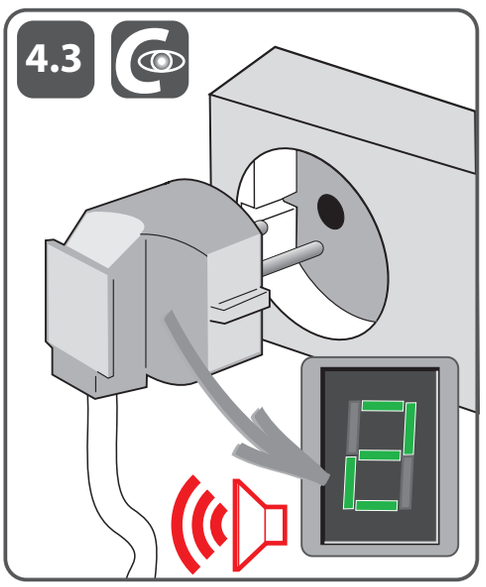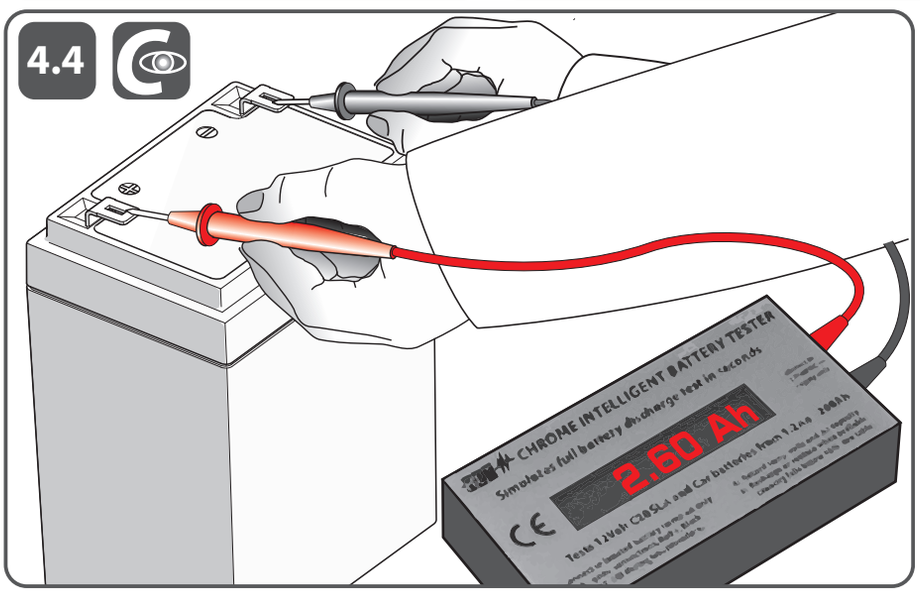
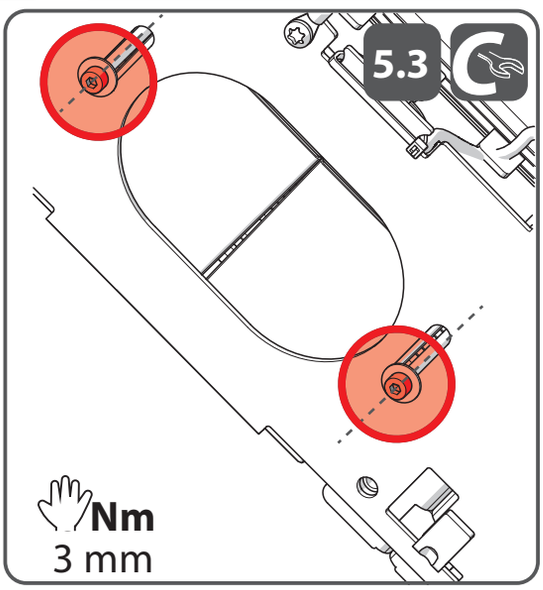
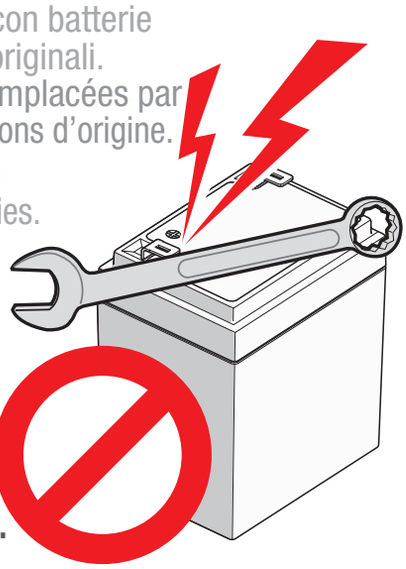
***Liquimolly* Battery clamp grease**
***RS* Pro contact treatment grease 503-335**


In a humid environment, the spray below can be applied to prevent oxidation of the battery contacts.
En un ambiente húmedo, se puede aplicar el siguiente spray para evitar la oxidación de los contactos de batería.
In vochtige omgeving kan onderstaande spray worden toegepast om oxidatie van de accucontacten te voorkomen.
In einer feuchten Umgebung kann das untenstehende Spray aufgetragen werden, um eine Oxidation der Akkukontakte zu verhindern.
Dans un environnement humide, le spray ci-dessous peut être appliqué pour éviter l'oxydation des contacts de batterie.
In un ambiente umido, lo spray sottostante può essere applicato per prevenire l'ossidazione dei contatti di batterie.

5 SERVICE LIFT


POWER PACK

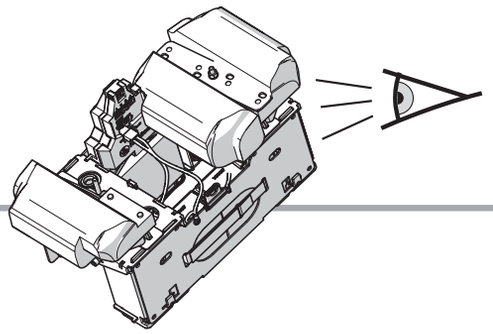5.1
5.2
5.3

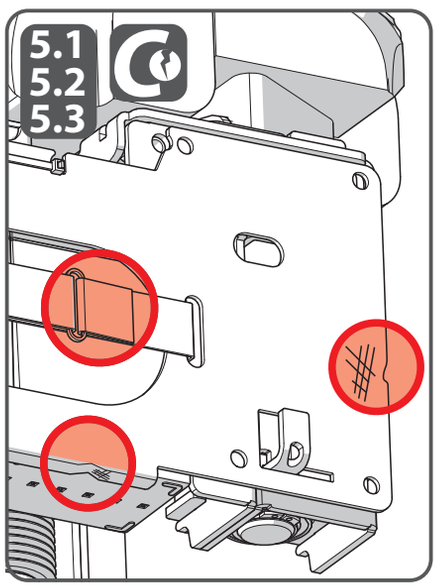
5.4

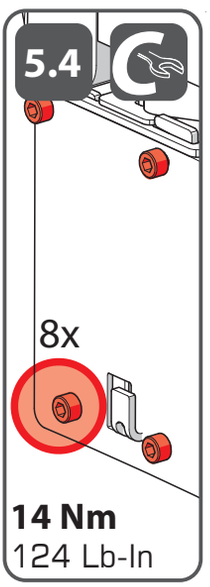
5.12


14 Nm
124 Lb-In

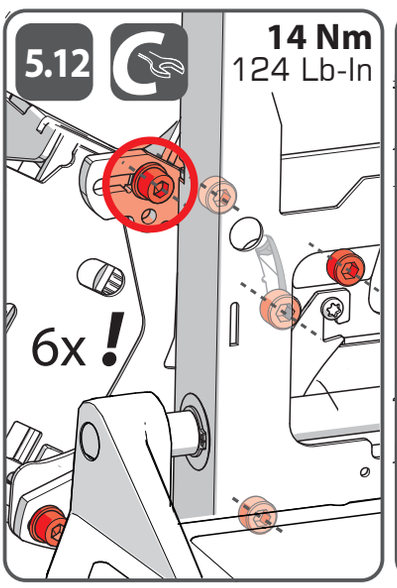
5.5

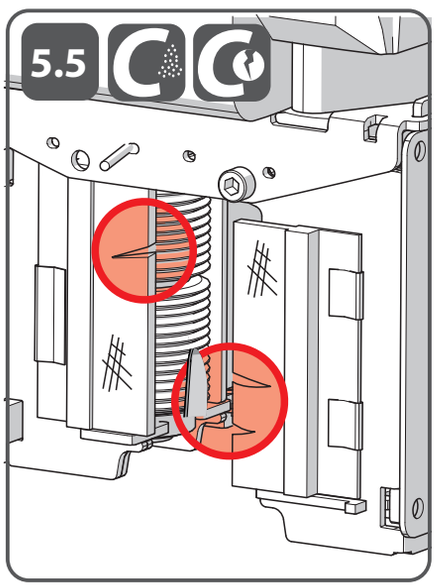5.6

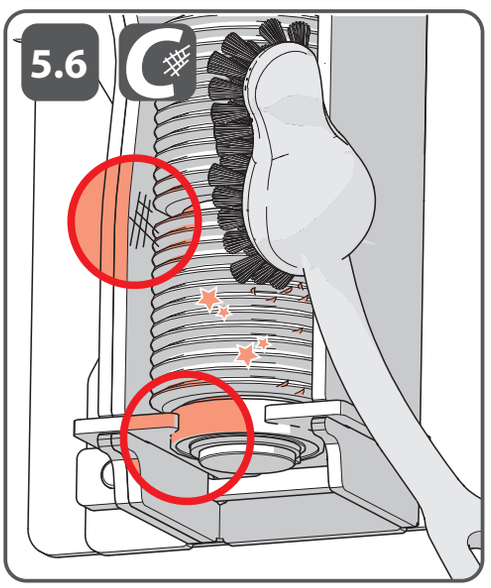
5.7

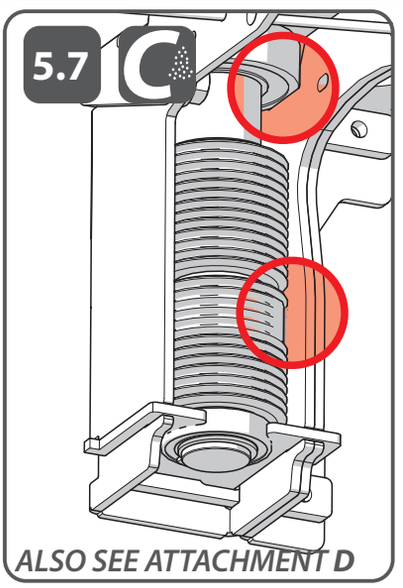
ALSO SEE ATTACHMENT D


5.8

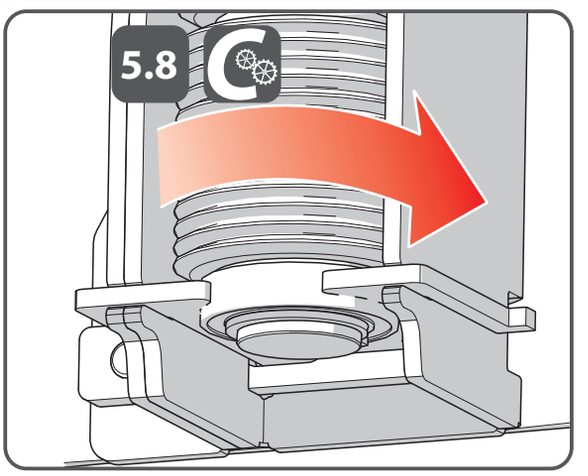
5.9

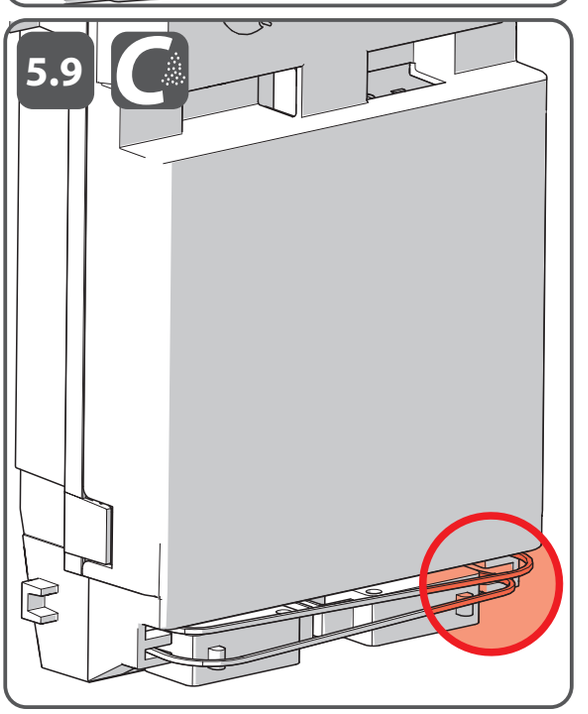
5.10
5.11


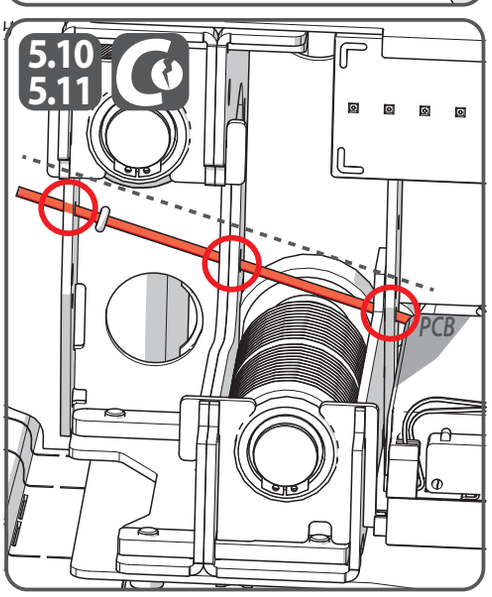
5.11
5.13
5.14


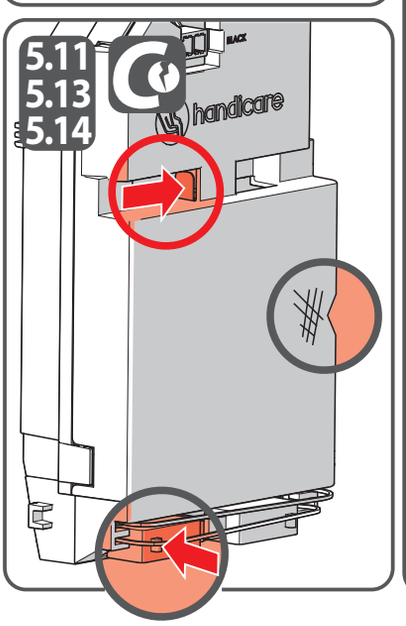
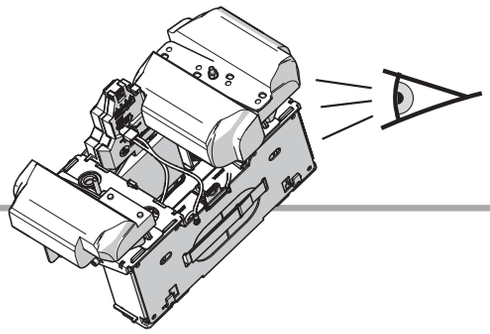
FOR EXAMINATION OF THIS SECTION REINSTALL ALL REMOVED PARTS OF THE LIFT (SEE ATTACHMENT B)

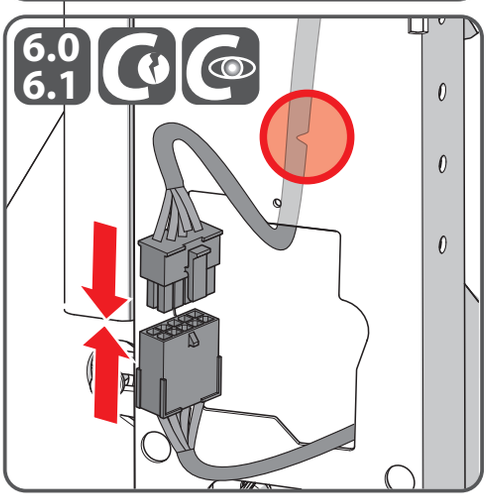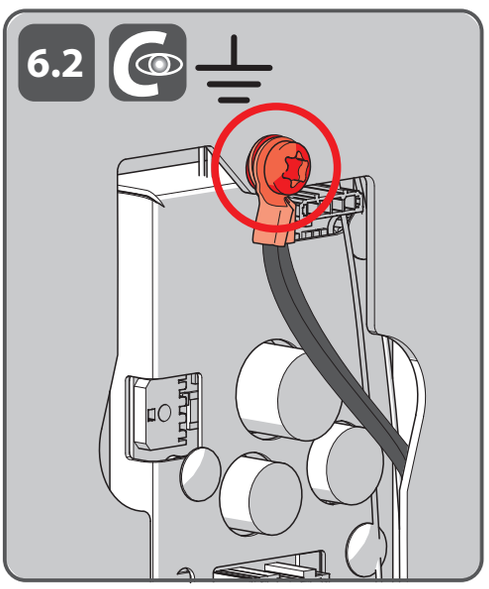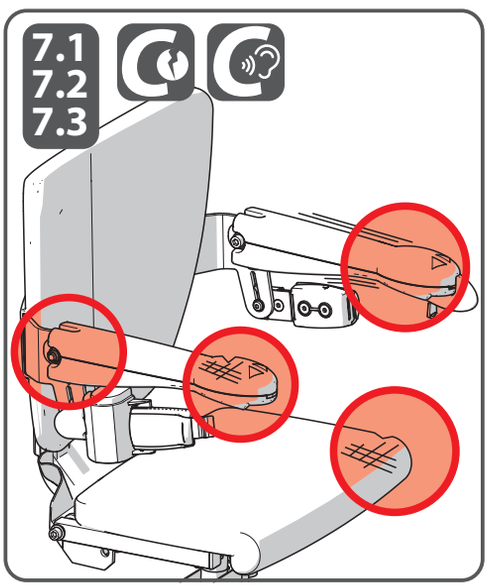
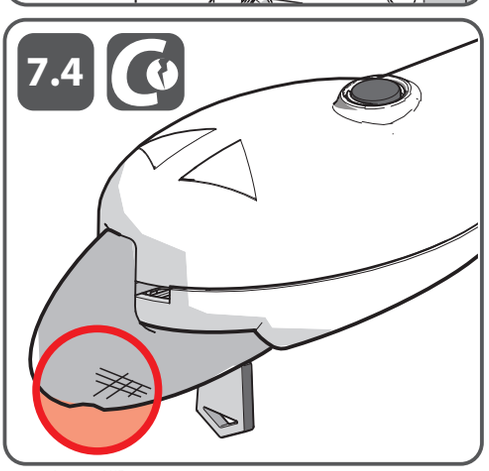
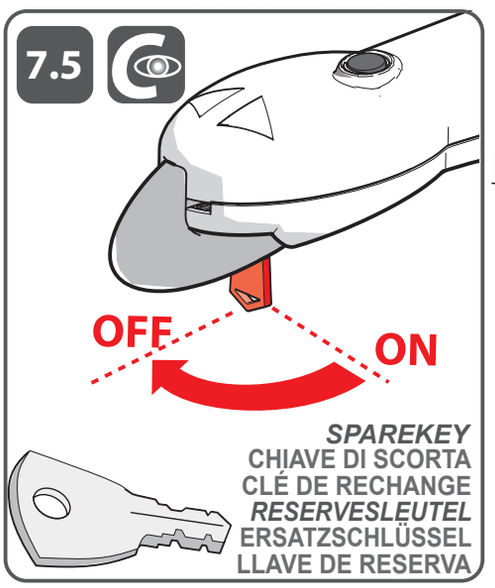
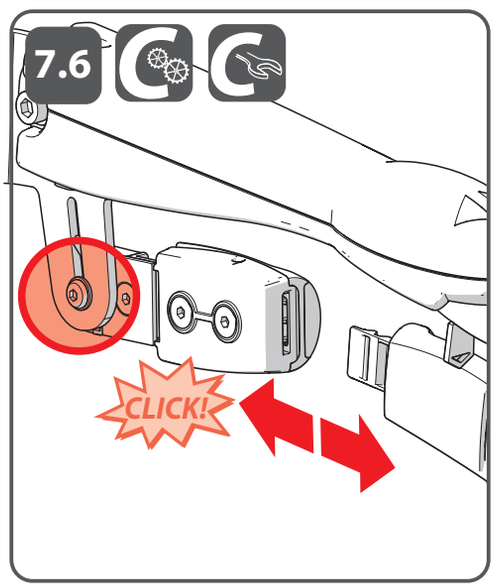
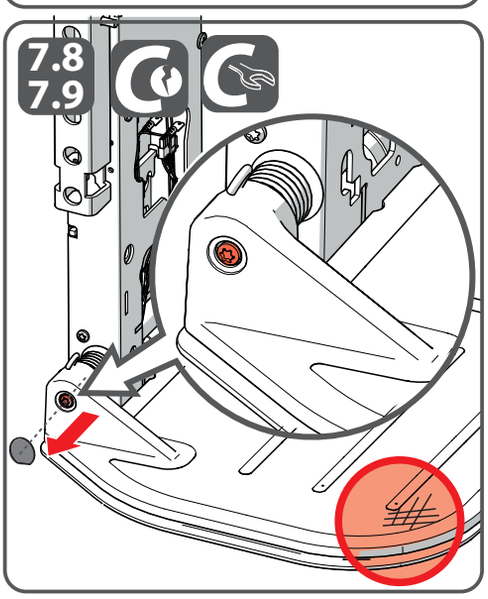
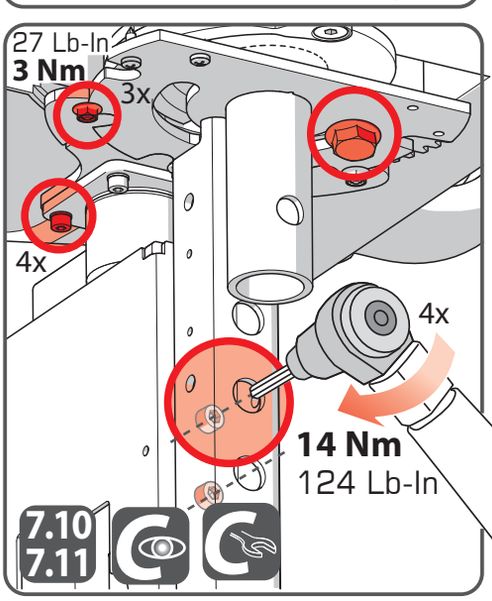
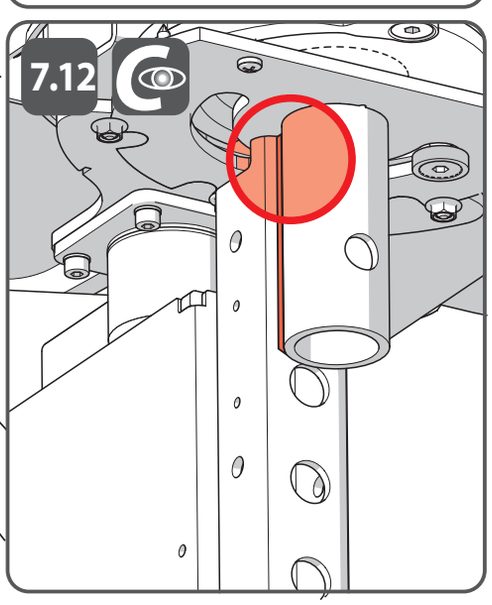
7 TEST SAFETY/FUNCTIONS


KEYSWITCH / JOYSTICK

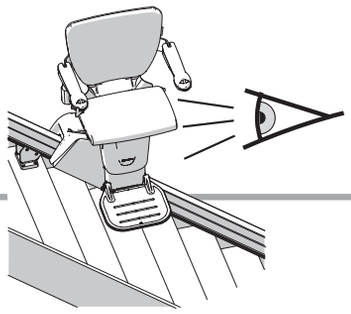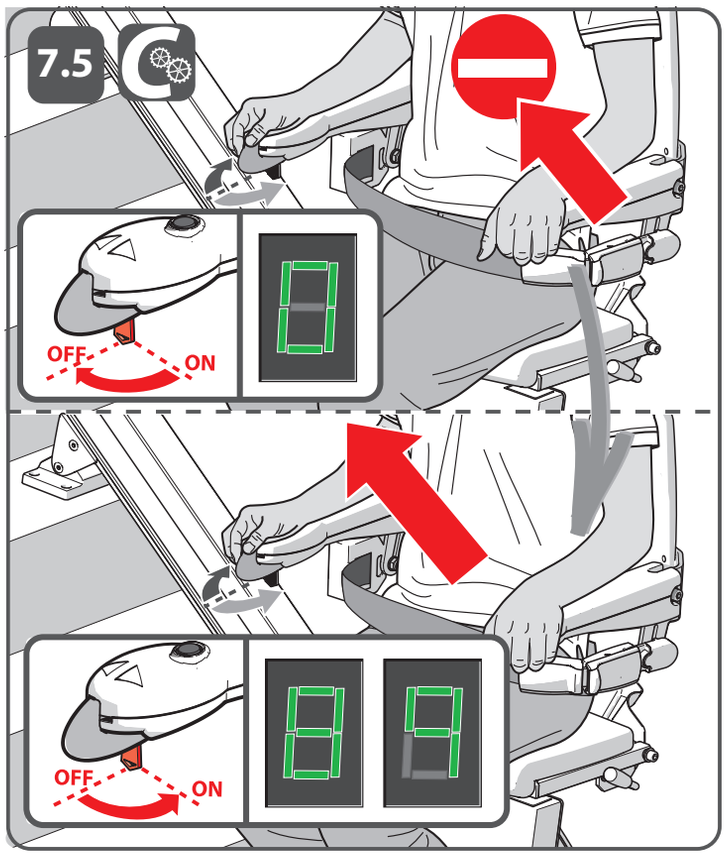
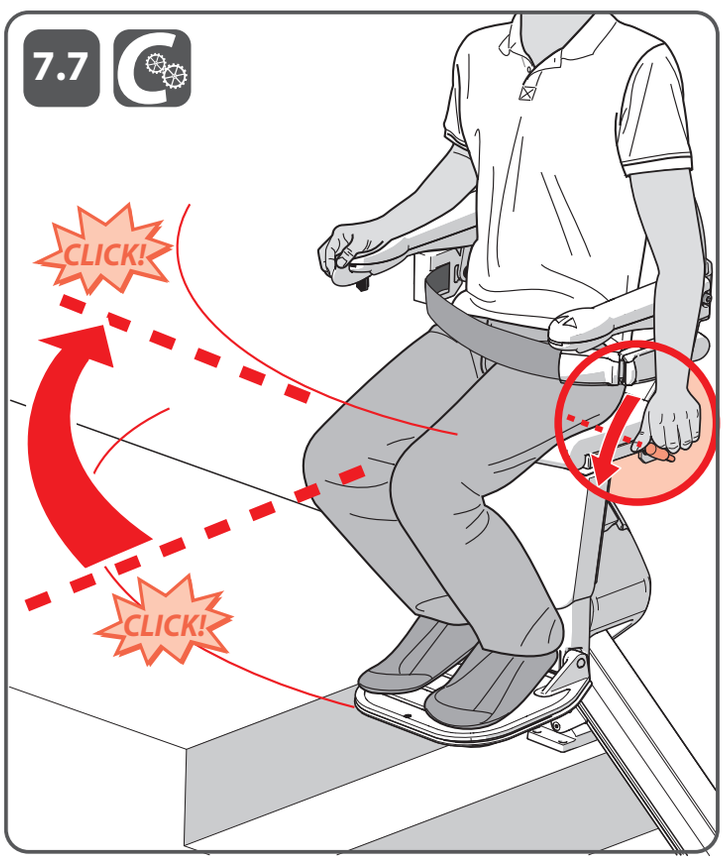
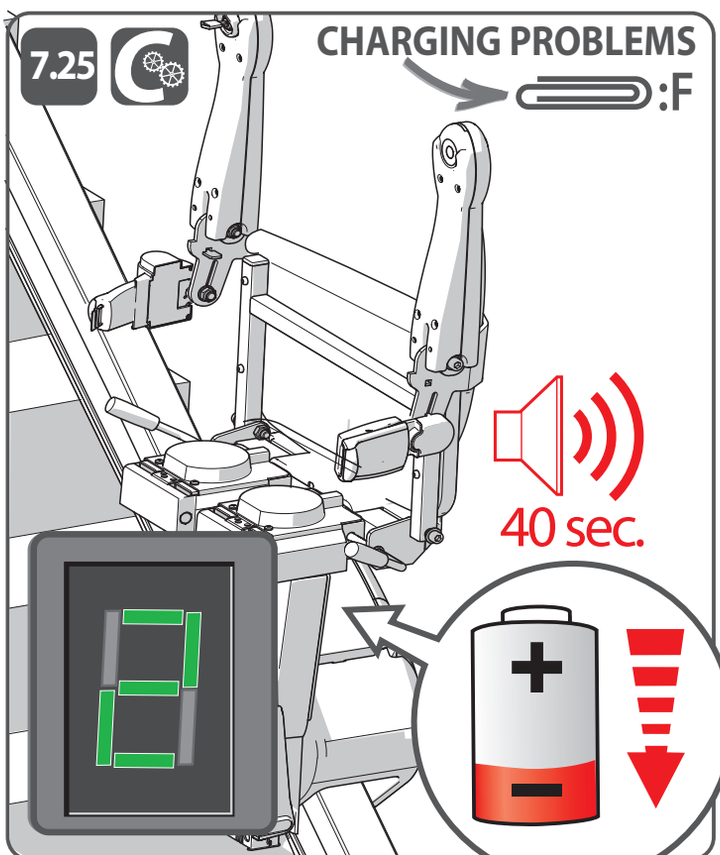
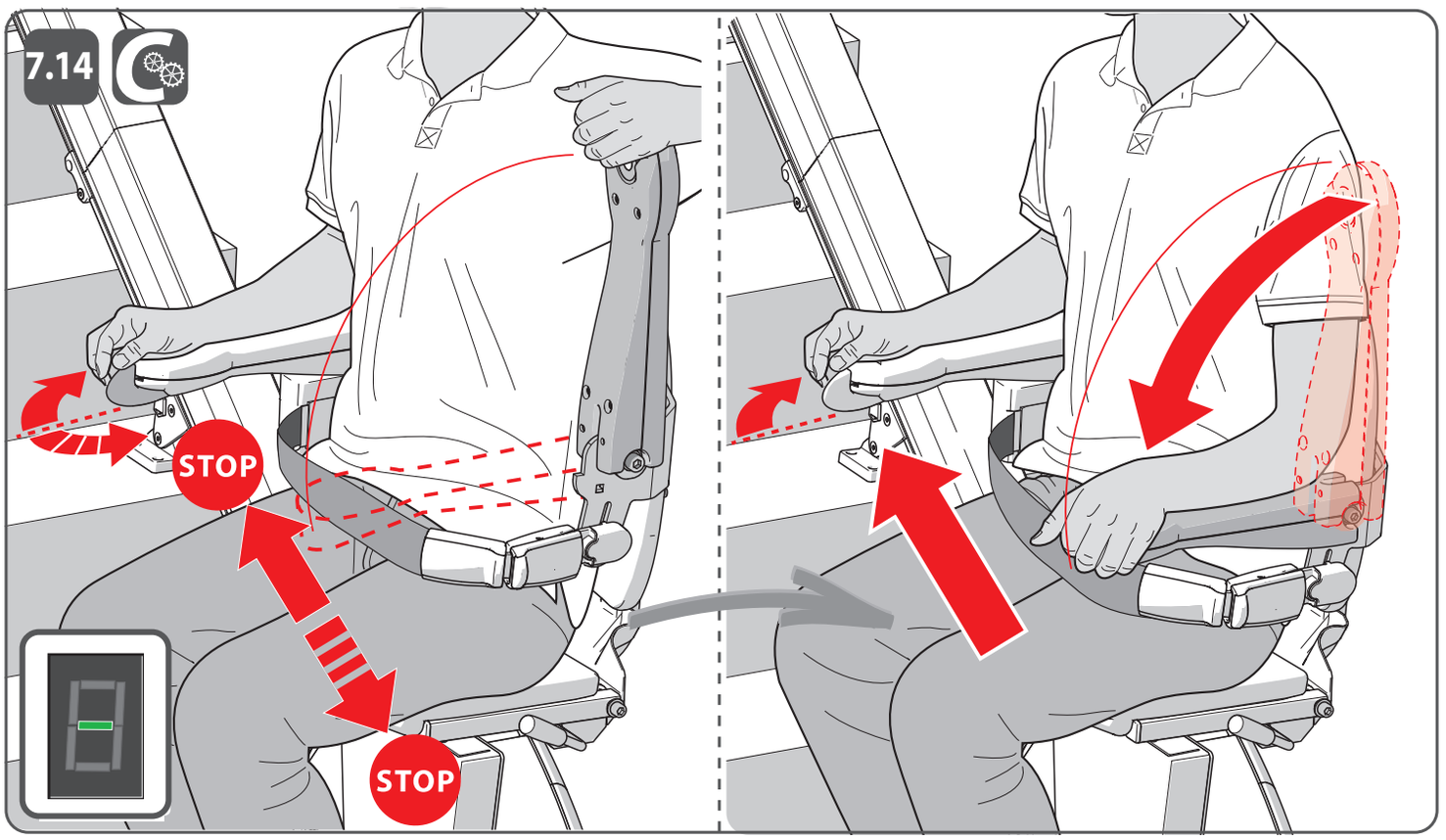
7.14
STOP
STOP

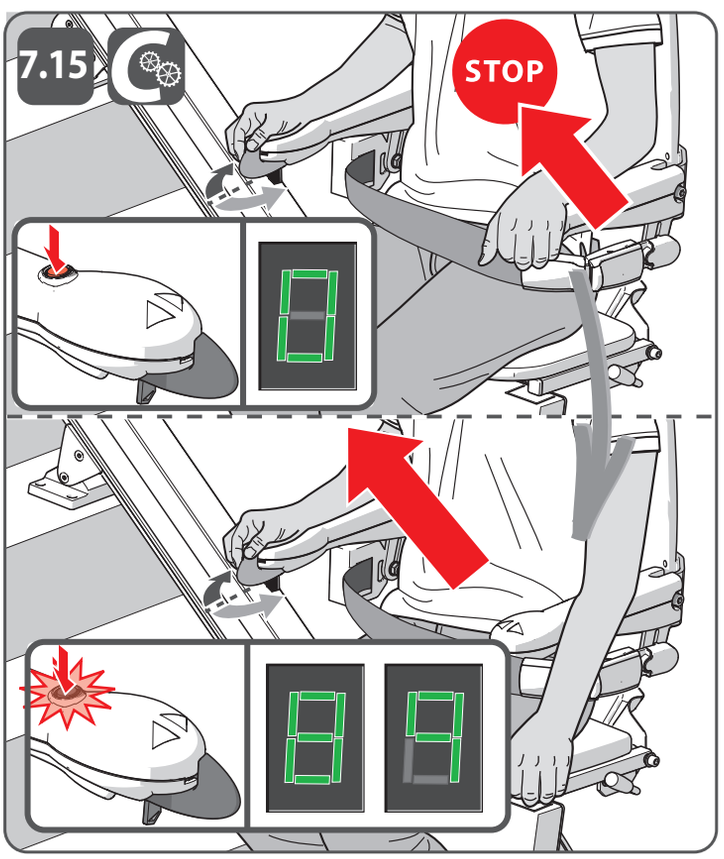
7.15
STOP

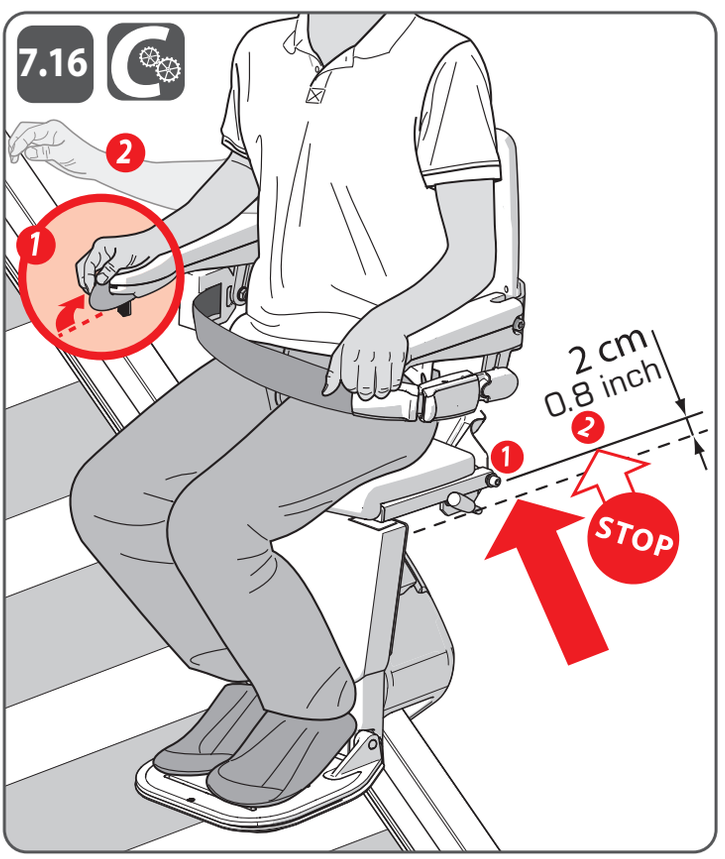
7.16
2
1
2 cm
0.8 inch
2
1
STOP


7 TEST SAFETY/FUNCTIONS

SAFETY LIMITS - STOPS

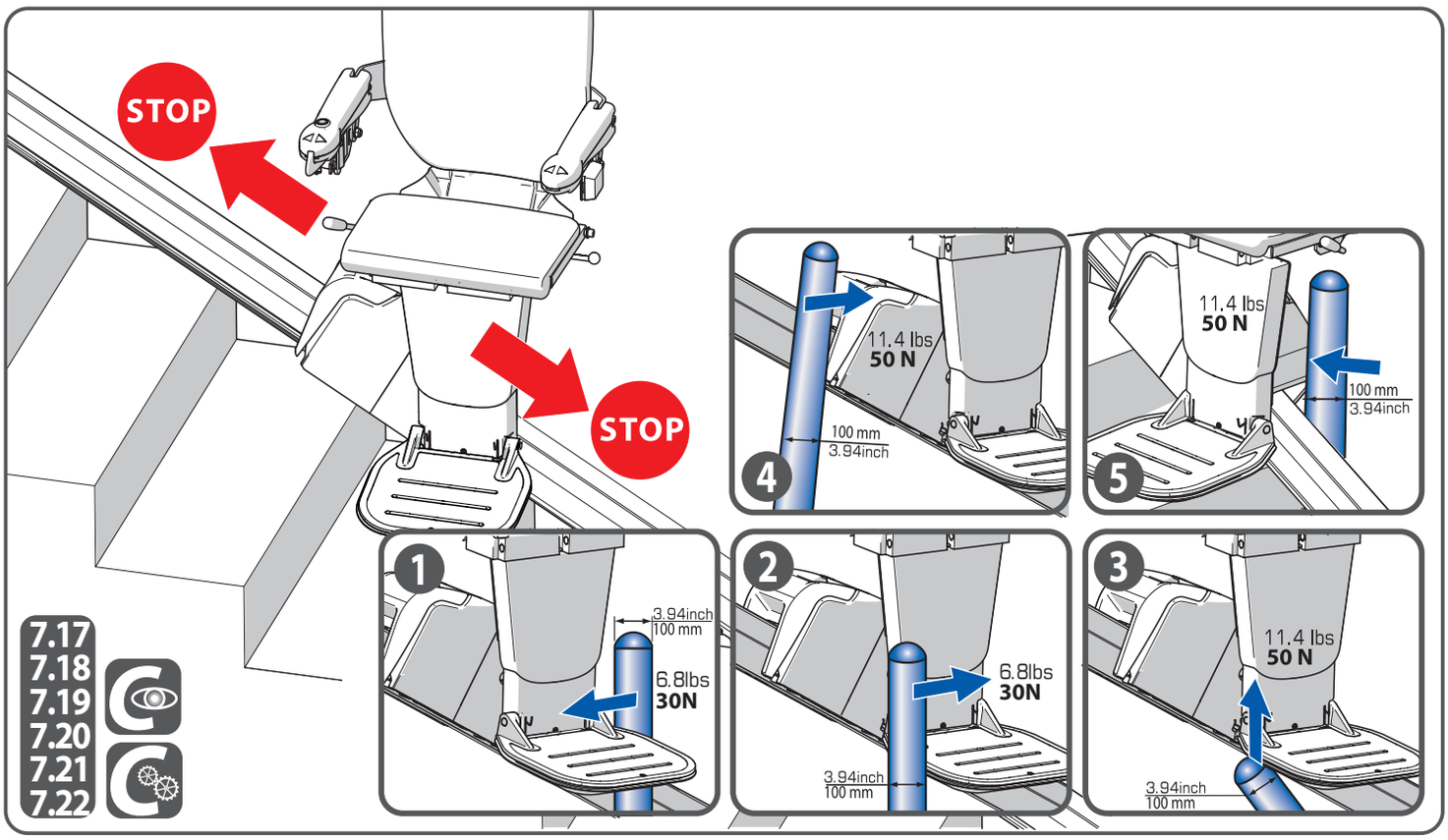
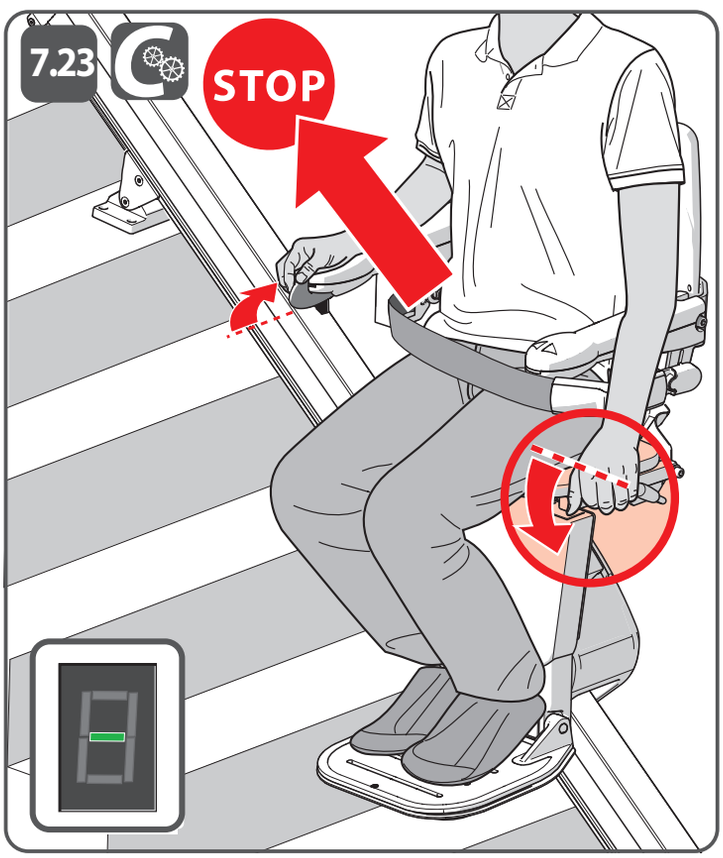
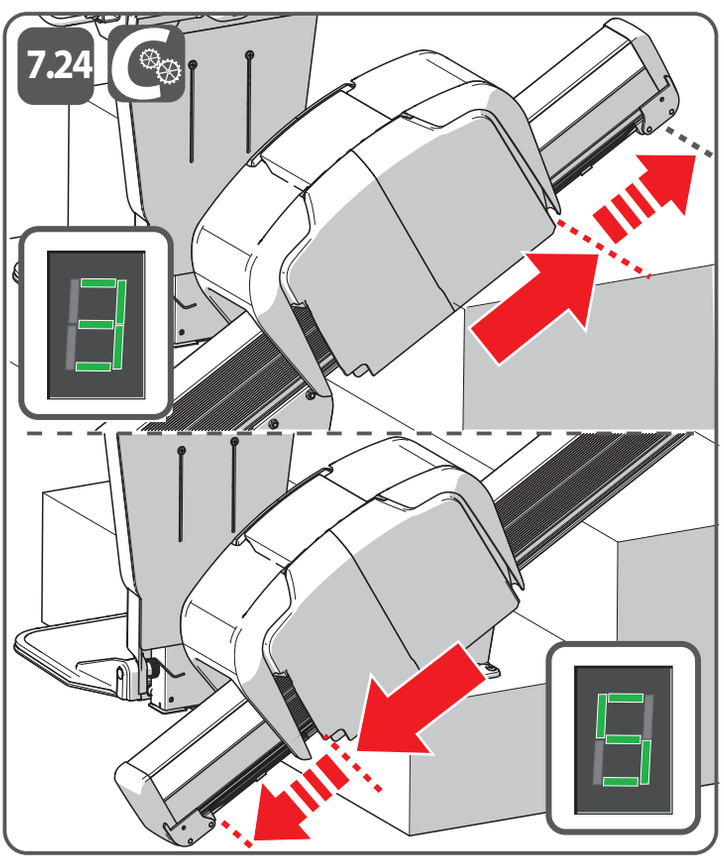
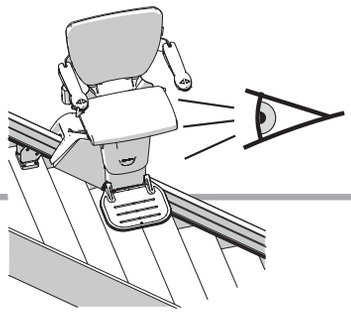
8 SERVICE LIFT OPTIONS


POWERED FOOTREST

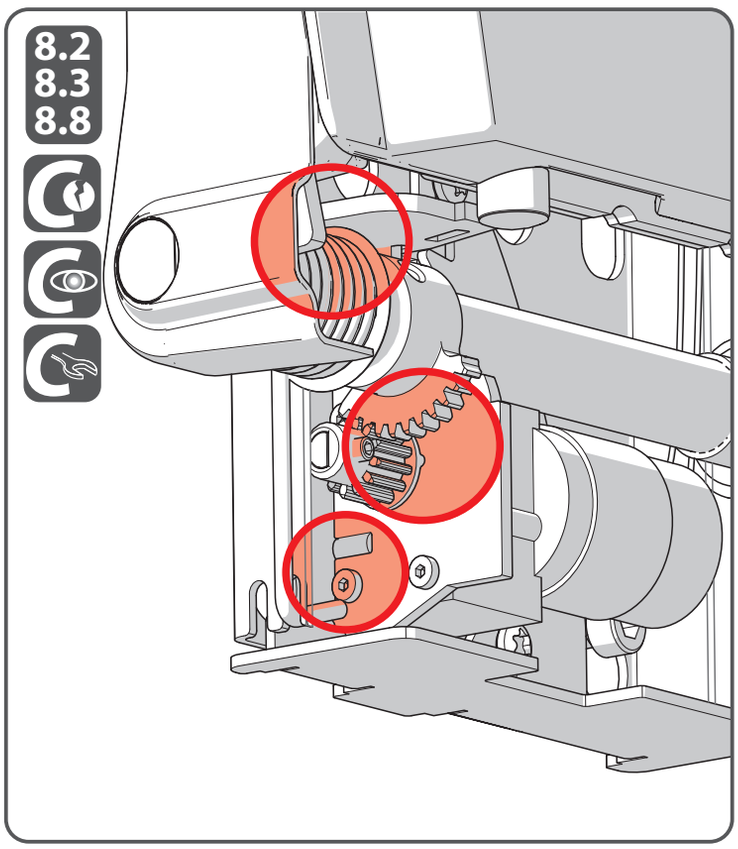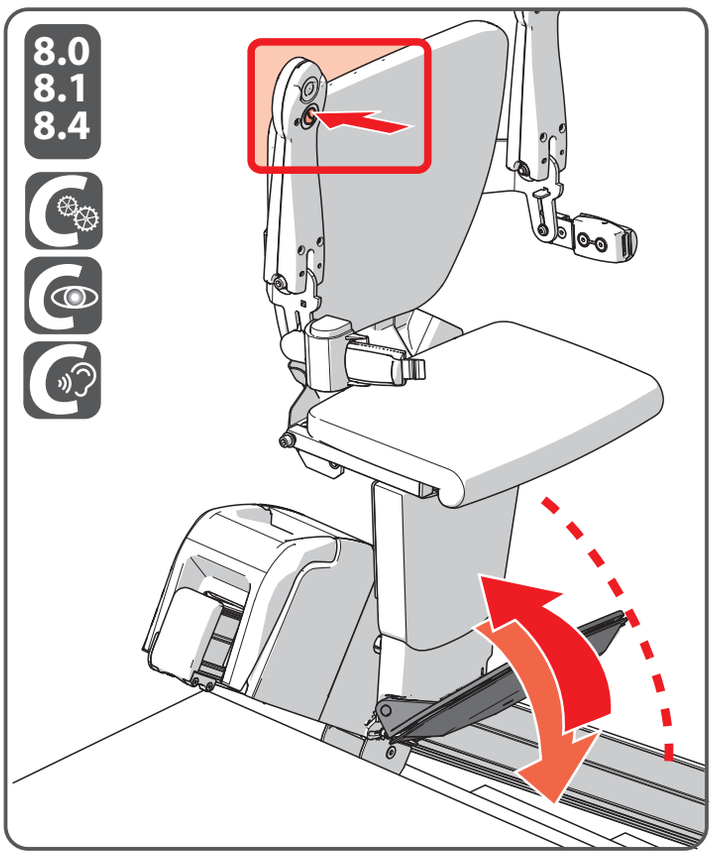
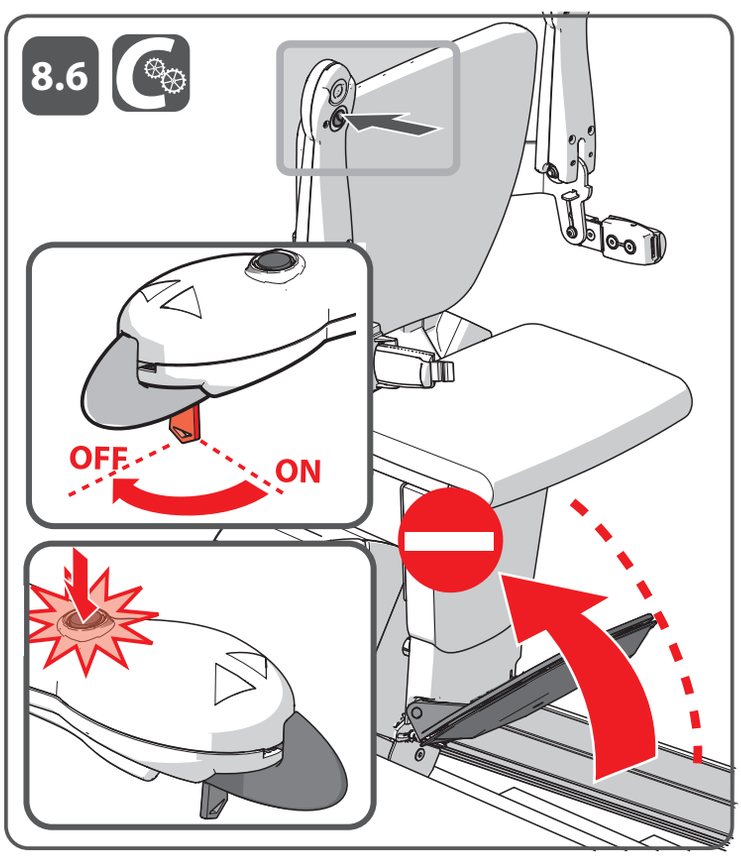
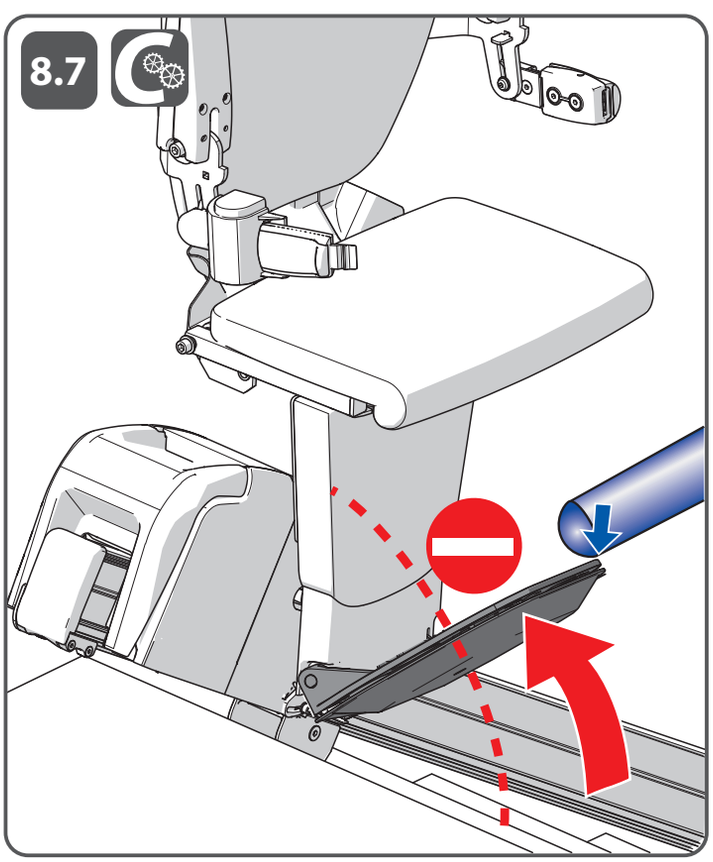
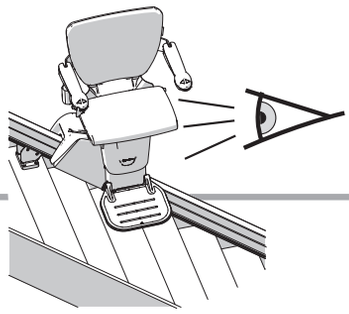
9 SERVICE LIFT OPTIONS


AUTOSWIVEL

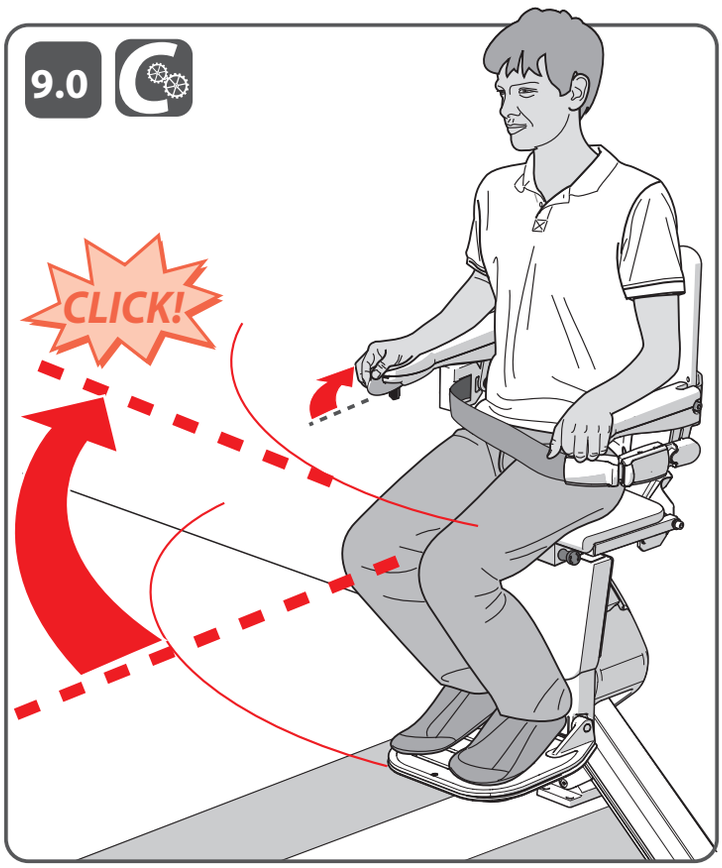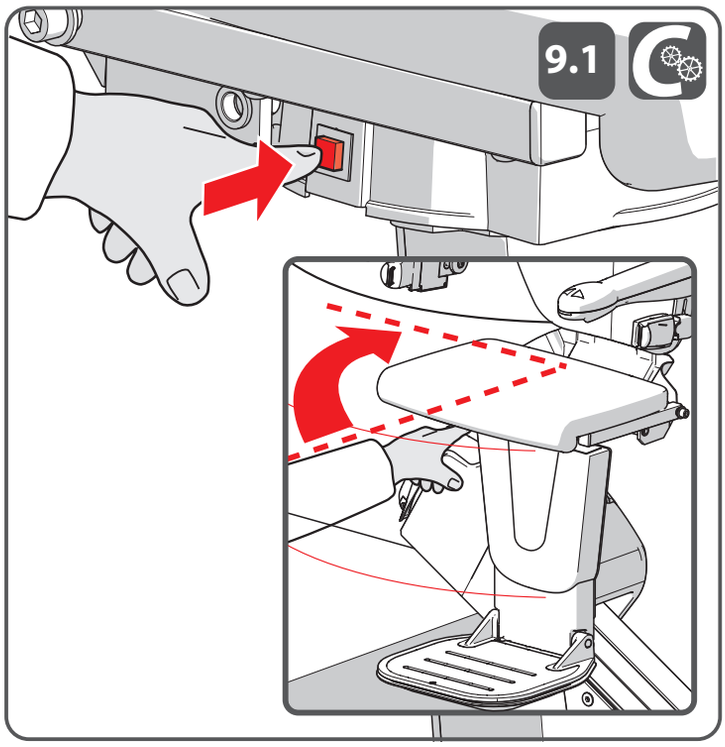
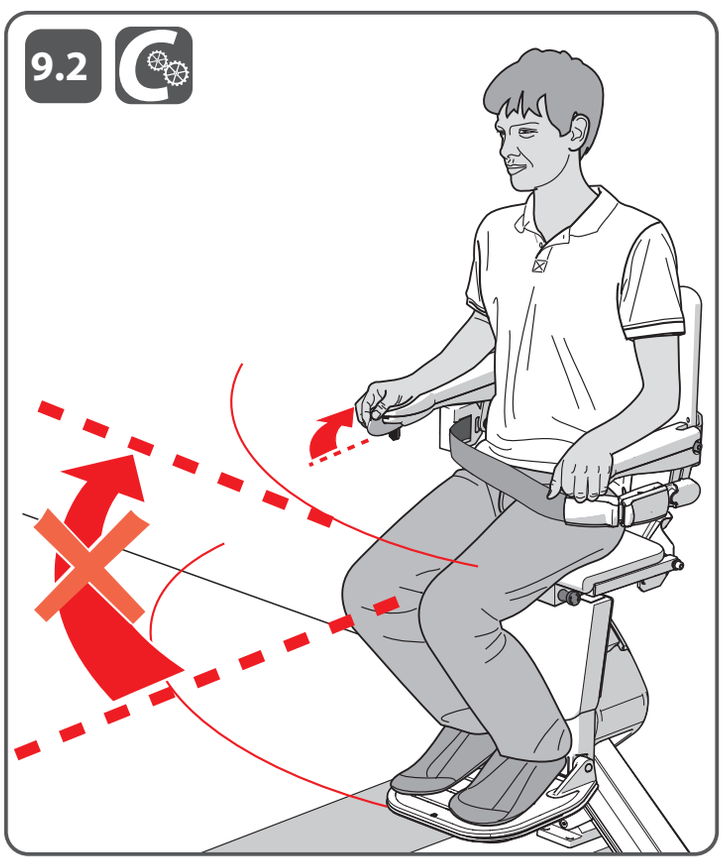
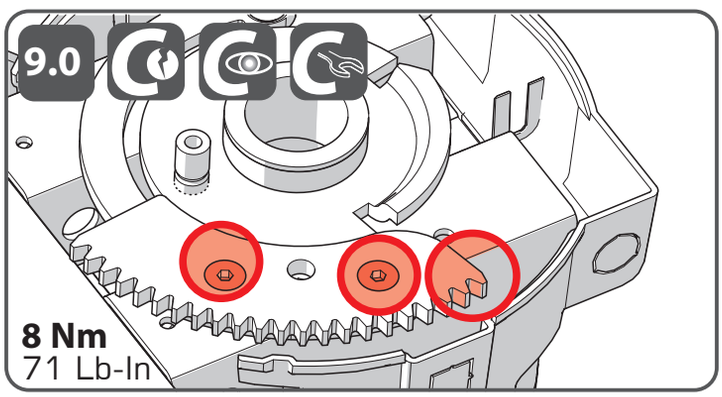
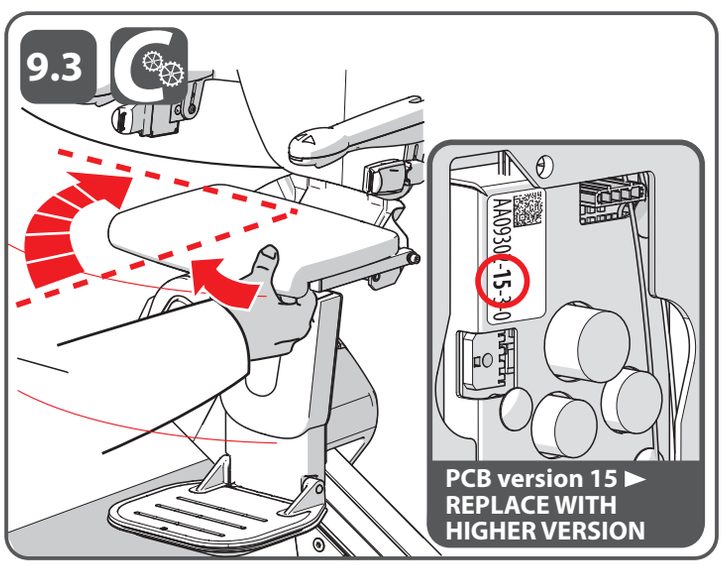
10 SERVICE LIFT OPTIONS


EASY HINGE

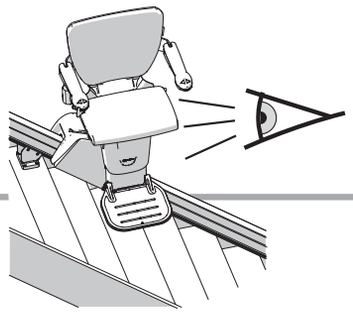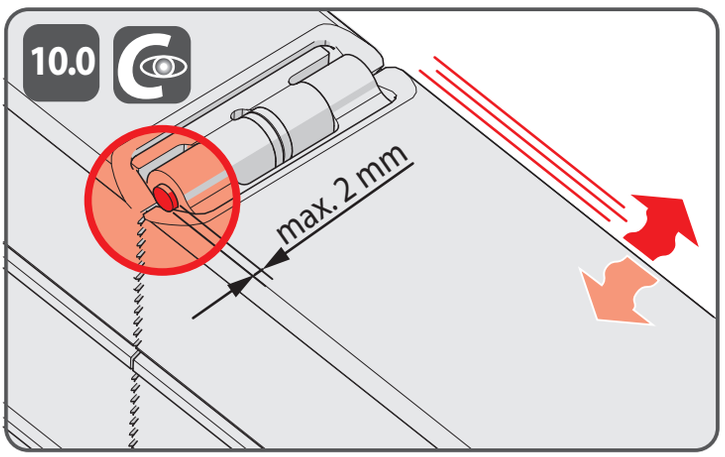
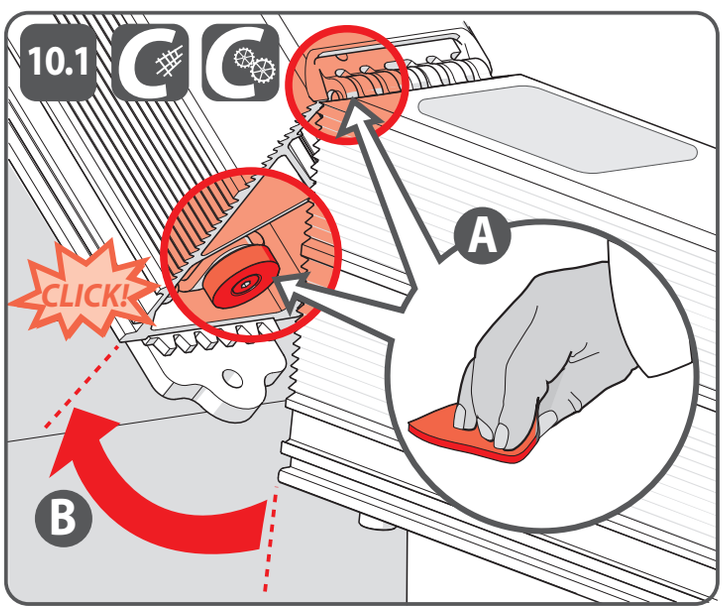
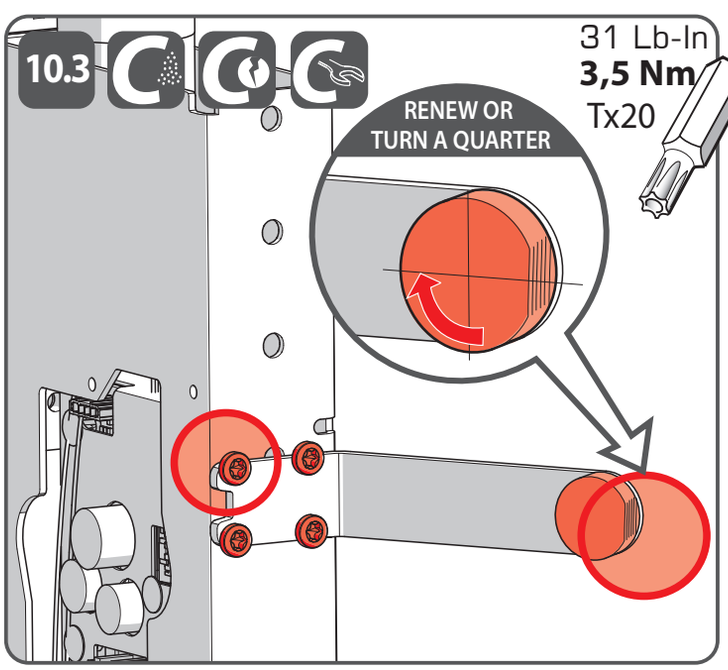
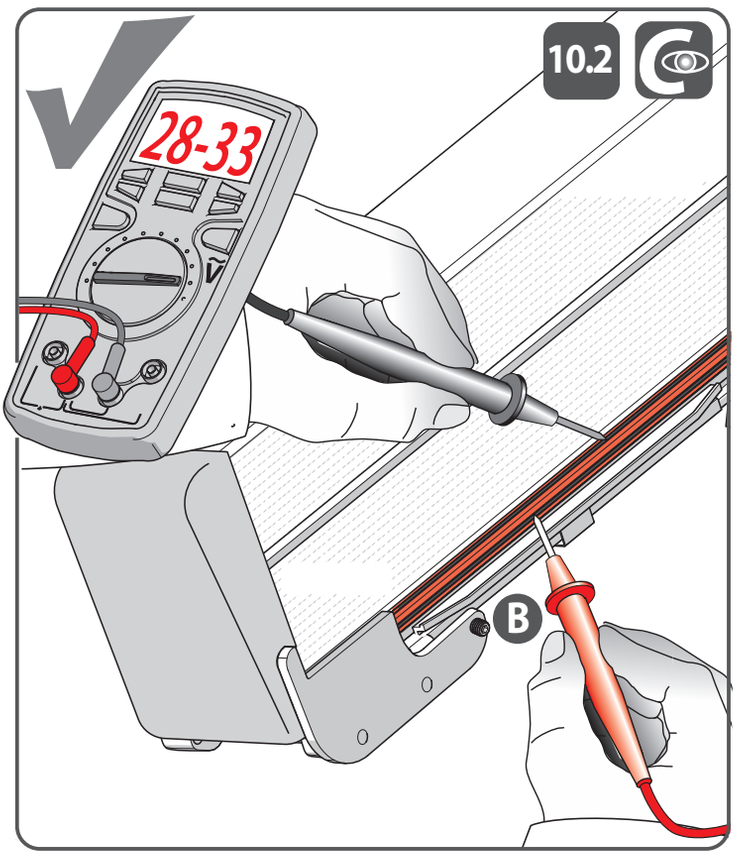
11 SERVICE LIFT OPTIONS


SLIDE TRACK

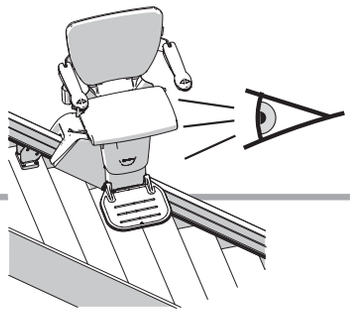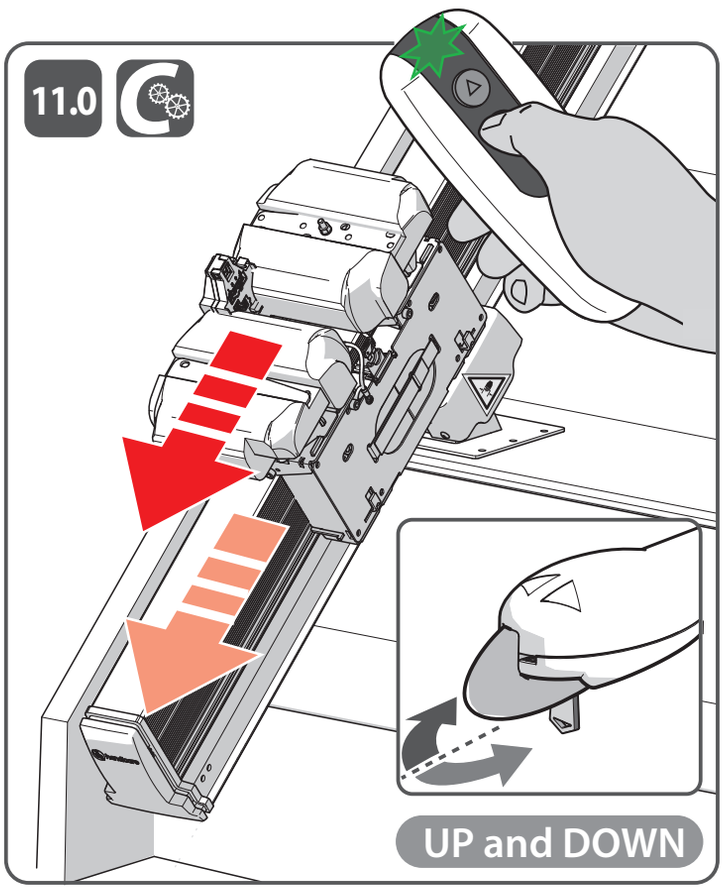
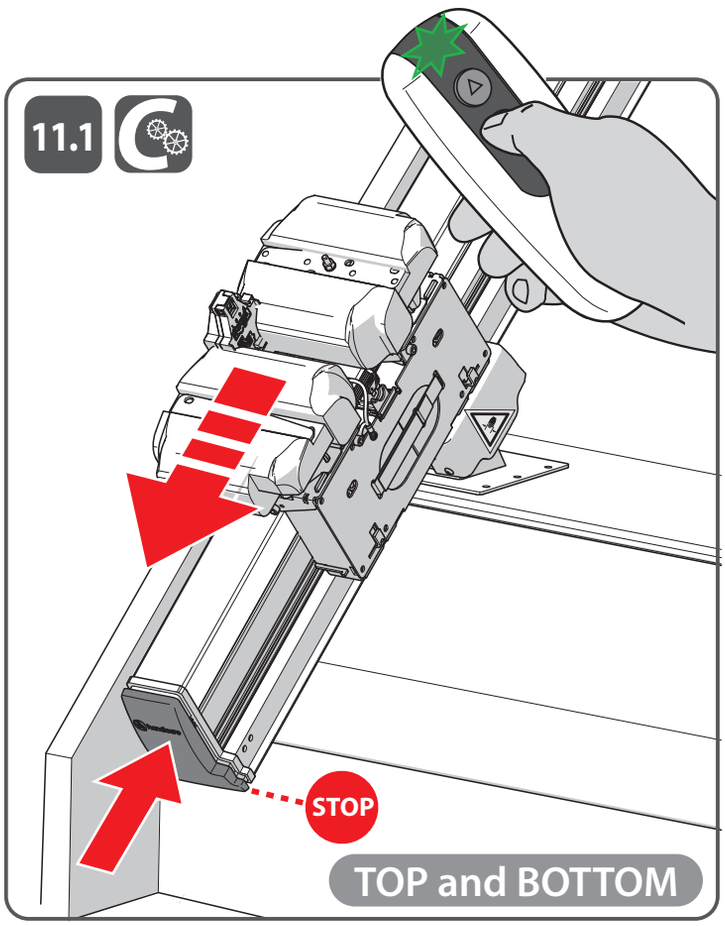
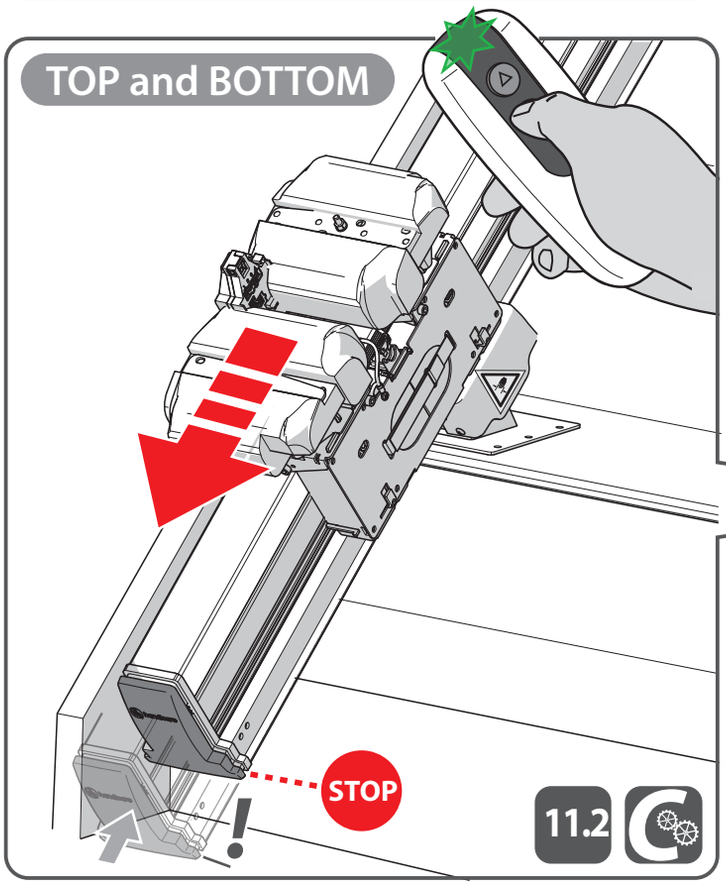
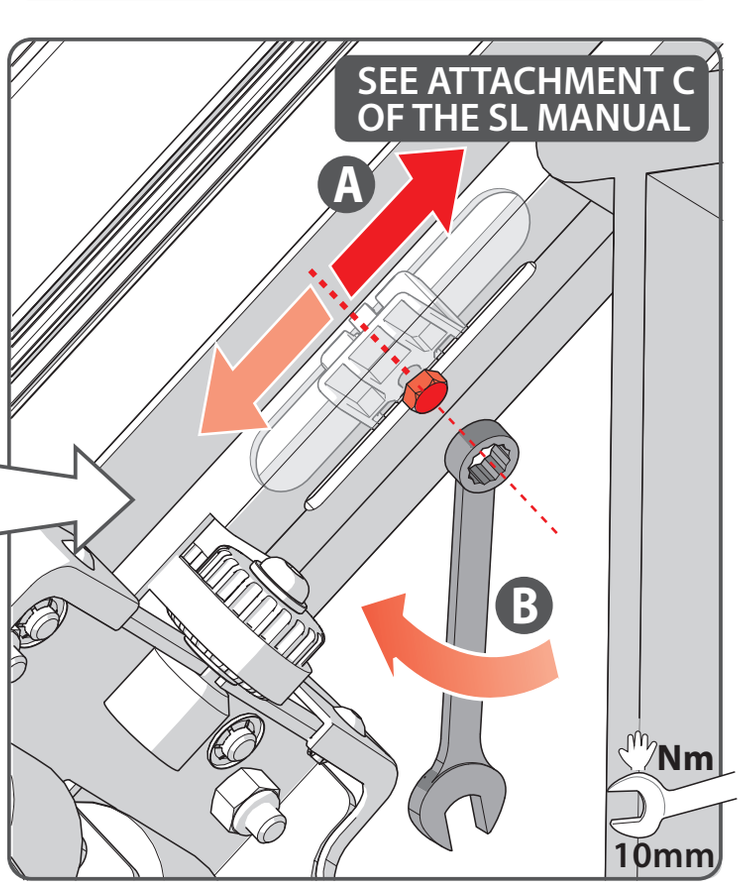
11 SERVICE LIFT OPTIONS


SLIDE TRACK

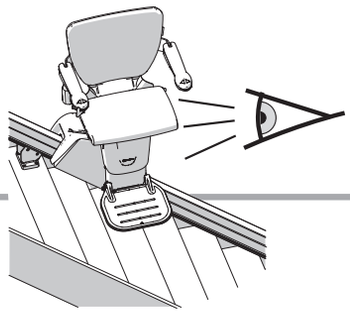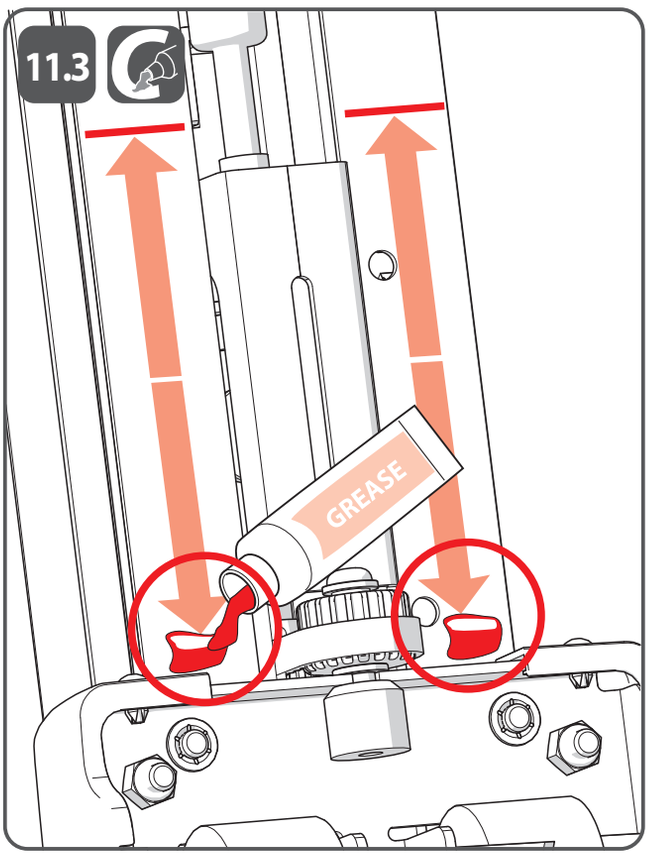
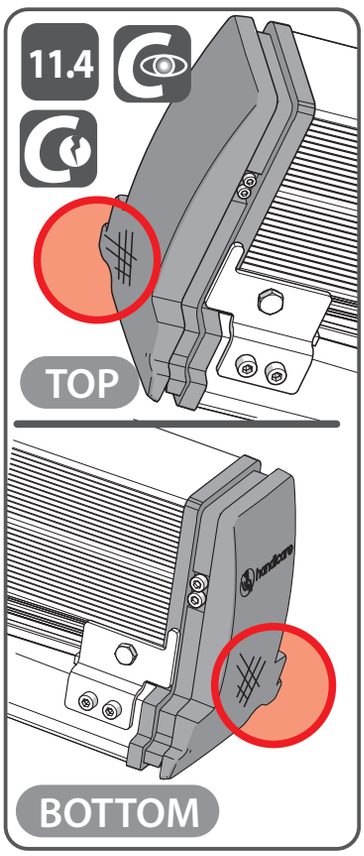
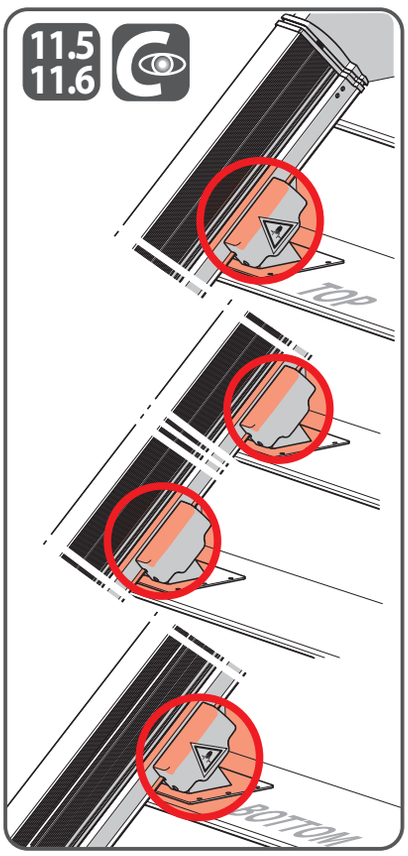
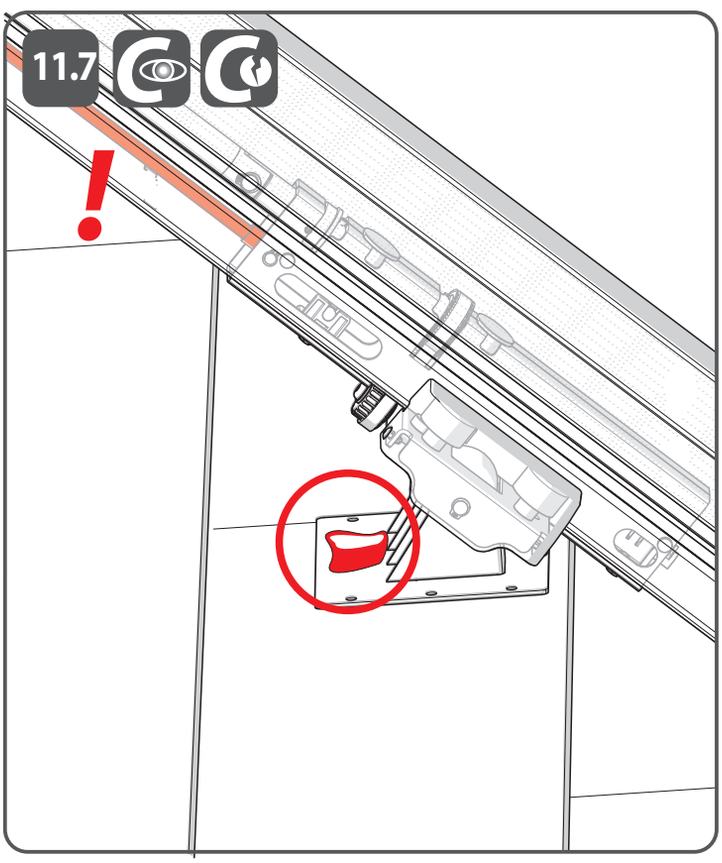
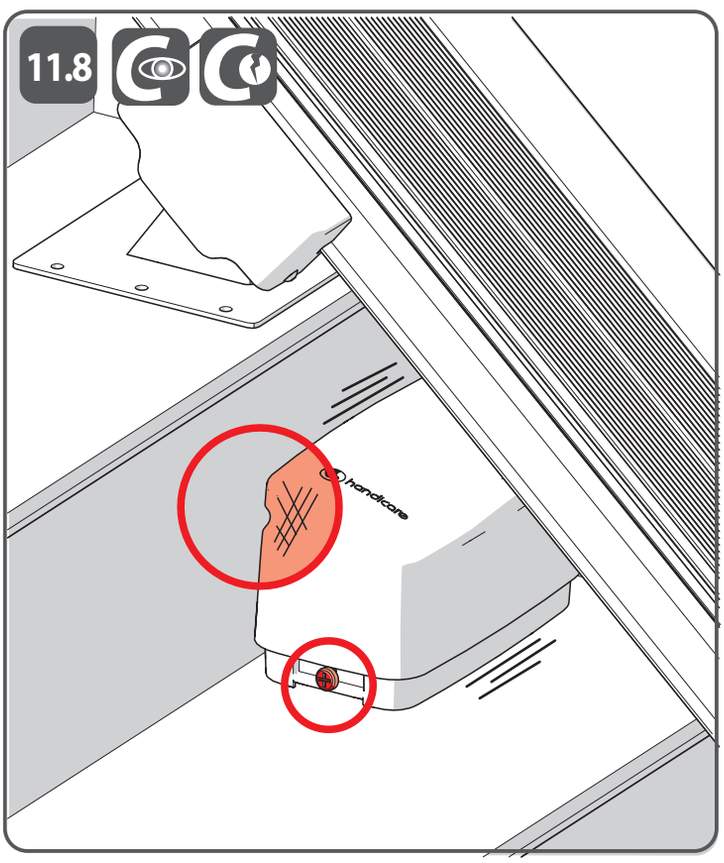
12 FINISH

FINAL CHECK

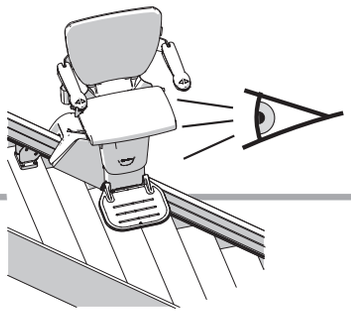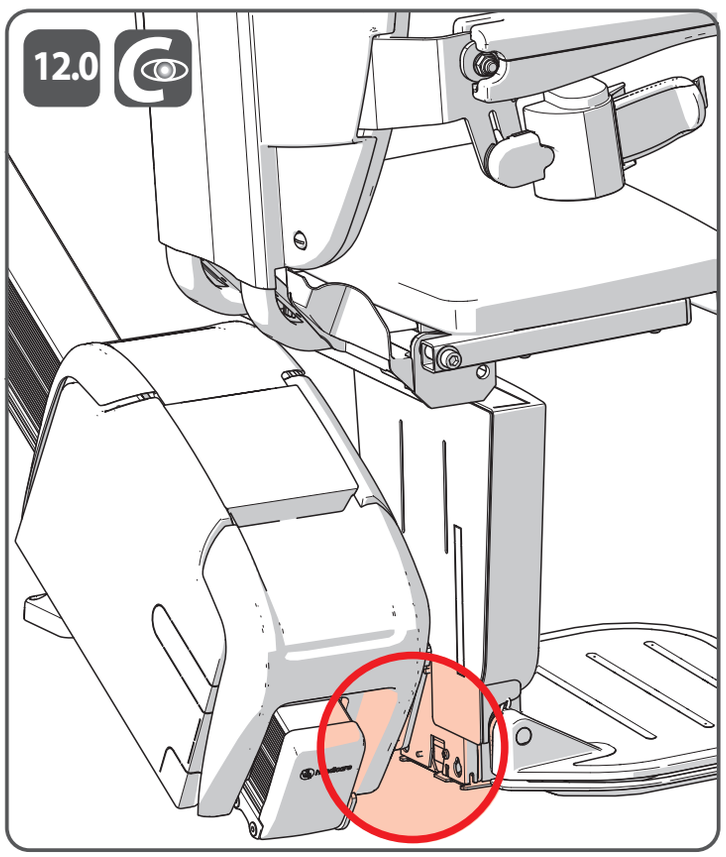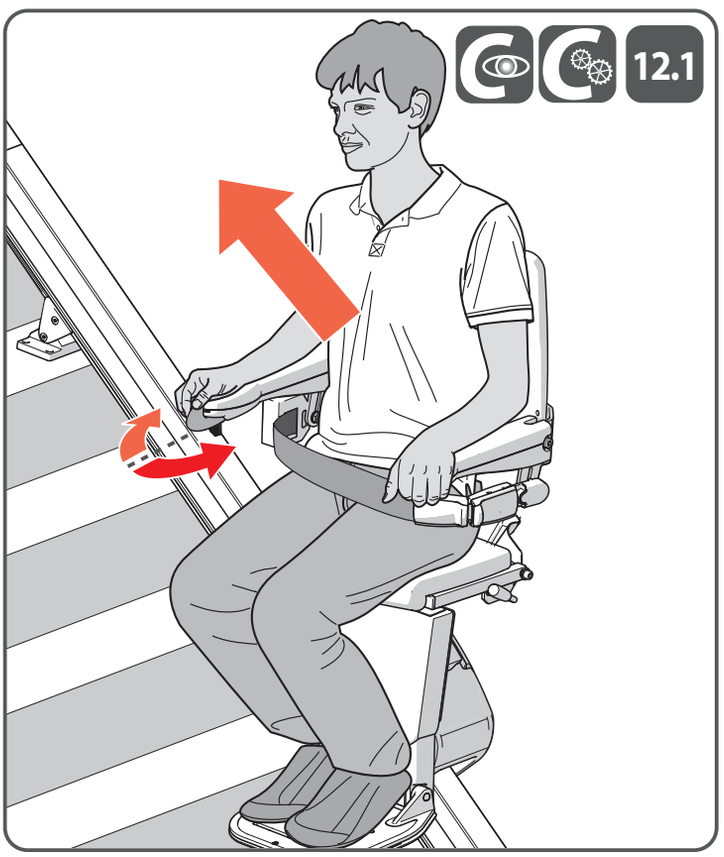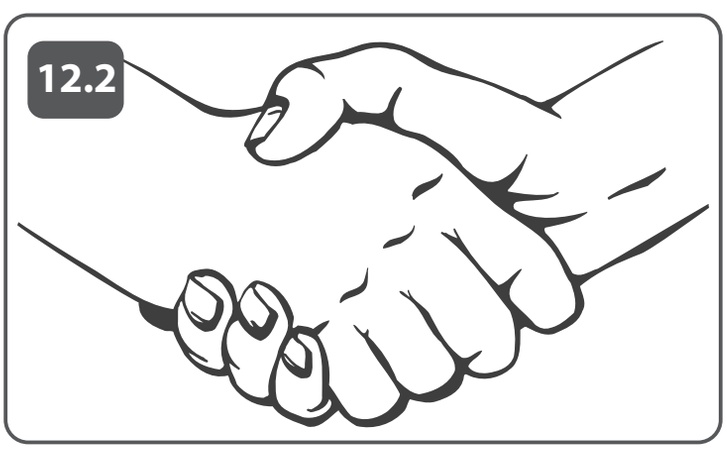
CERTIFICATE
CERTIFICAAT
CERTIFICADO
ZERTIKAT
CERTIFICAT
CERTIFICATO

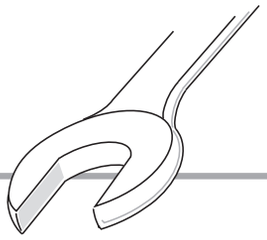CERTIFICATE FOR PERIODIC EXAMINATION HANDICARE 1100

Name:		Additional details:
Adress:		
Zipcode:		
City:		
Telephonenumber:		
E-mail:		
Rail Type:	1100	Projectnumber:
Type chair:	Manual seat / Powered seat	Unitnumber:
Options:	Non / Powered footrest / Powered swivel	Printnumber:
Rail:	Standard / Slide / Hinge	Year:

0. START		✗	✓
0.0	Meet customer and ask for his **remarks** about the lift		
0.1	Do a **test ride** under load and check if the lift is functioning okay		
0.2	Remove **the key** from the switch		


1. STANDARD RAIL		✗	✓
1.0	Check the rail for **damages**		
1.1	Clean the rail for **dirt**		
1.2	Check the **grooves** for damages/irregularities		
1.3	Check all **fixings** of the rail		
1.4	Check if all rail parts are **aligned**		
1.5	Check the fixings of the **supports** on the stairs		
1.6	Check the fixings of the **supports** on the rail		
1.7	Check if the fixation of the supports are still **tight**		
1.8	Check the presence of both **final limit** cams and mechanical endstops		
1.9	Check the charger **output** to rail		
1.10	Clean the **rail caps** for dirt		

2. CHARGER		✗	✓
2.0	Check the function of the **charger**		
2.1	Check the **wiring** for damages		
2.2	Check the connection to the **charging strip**		
2.3	Check that the charger does not make **noise**		
2.4	Check that the plug is in the socket and the **safety label** is placed.		

3. REMOTES		✗	✓
3.0	Replace the **batteries** in the remotes (Penlight AA)		
3.1	Check remotes for **damages**		
3.2	Check if the light in the remotes **illuminates**.		
3.3	Check if the remotes are fitted correctly to the **wall**		
3.4	Check the **function** of all remotes		

For examination of this section the chair and the covers of the unit should be removed and the unit should to be taken of the rail

4. BATTERIES		✗	✓
4.0	Check batteries for damages and corrosion on the **connections**		
4.1	Check **wiring** battery-connections for damages		
4.2	Check **fuseholder** for corrosion		
4.3	Check if the **display** shows a "**2**" (and the lift will beep) when not charging, and shows a "**-**" when charging		
4.4	Check the **capacity** of the batteries with a battery tester (battery 1=..... , battery 2=.....)		
4.5	Check the charge **voltage** from the PCB to the batteries		

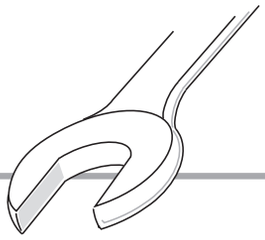
5. UNIT / FRAME		✗	✓
5.0	Check the **unit** for required replacement (see attachment D)		
5.1	Check the **unit frame** for damages		
5.2	Check batteries **support** for damages		
5.3	Check the fixation of the **strap/bracket** on the batteries		
5.4	Check torque on **bolts** (backside of unit) (14Nm \| 124 Lb-In)		
5.5	Check the b**rushes on the power pack** for damages		
5.6	Clean the space around **rollers** and bearings and remove all dirt		
5.7	Check rollers and **bearings** for wear		
5.8	Check if **rollers** and bearings run smoothly		
5.9	Check the charging **springs** for shape, wear and functionality (spring force)		
5.10	Check the shape of activation **rods**		
5.11	Check if activation **rods** move smoothly and activate the switches of the print circuit board		
5.12	Check the fixation of the **adaptor plates**		
5.13	Check the **print circuit board** for damages		
5.14	Check the function of the **final limit** switches in the unit		
5.15	Clean the brushes on the power pack		

For examination ot this section reinstall all removed parts of the lift

6. WIRING		✗	✓
6.0	Check the **wiring** for damages		
6.1	Check if the wiring is **fitted** correctly		
6.2	Check the **earth wire** connections on the frame, unit and rail		

7.CHAIR		✗	✓
7.0	Check the presence of the **CE sticker**		
7.1	Check the chair for **damages**		
7.2	Check for too much **noise**		
7.3	Check the **armrest** and the armrest hinge for damages		
7.4	Check the function of the **joystick** and for damages		
7.5	Check the function of the **keyswitch** and the presence of the **spare key**		
7.6	Check the fixation of the safety belt bracket and function of the **safety belt**		
7.7	Check the **mechanical lock** of the seat		
7.8	Check the **footrest** for damages		
7.9	Check the fixation of the **footrest**		
7.10	Check presence and fixation of the bolts and **securing bolt** under the seat		
7.11	Check the **fixation** of the chair		
7.12	Check all welding of the chair frame **pivot**		
7.13	Check the **perpendicularity** of the chair with a spirit level		
7.14	Check functionality of the **armrest switches**		
7.15	Check functionality of the **stop button**		
7.16	Check **accuracy** of the joystick		
7.17	Check the function of the **chair support** detection switch		
7.18	Check the chairsupport detection is less than **50 N** in the up direction		
7.19	Check the function of all **cover** detectionswitches		
7.20	Check the cover detection is less than **50 N** in both directions		
7.21	Check the function of all **footrest** detection switches		
7.22	Check the footrest detection is less than **30 N** in 3 directions		
7.23	Check the function of the **chairposition** switch under the seat		
7.24	Check the function of the **final limit** switches		
7.25	Check for ***charging problems***		

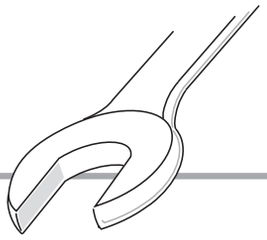	***8. OPTIONS >** POWERED FOOTREST*	✗	✓
8.0	Check that the **button** of the powered footrest is functioning		
8.1	Check that the footrest **(un)folds** smoothly		
8.2	Check the **gear and rack** segment of the folding mechanism for wear		
8.3	Check the **spring** for damage		
8.4	Check the sound of the **motor** (replace when the noise is too much)		
8.5	Check the folding **mechanism** for damage		
8.6	Check that the footrest is not functioning if the **emergency button** is pressed or the **keyswitch** is off.		
8.7	Check that the **safety** is functioning (when the footrest is (un)folding, it must stop when blocked)		
8.8	Check that the **bolts** of the motorframes are still fixed		

	***9. OPTIONS >** AUTOSWIVEL*	✗	✓
9.0	Check the function of the **autoswivel**		
9.1	Check the function of the **release button**		
9.2	Check that the autoswivel is not functioning if the **emergency button** is pressed or the **keyswitch** is off		
9.3	Check the ***resistance*** of the autoswivel		


	***10. OPTIONS >** EASY HINGE*	✗	✓
10.0	Check the **joint**		
10.1	Clean all parts of the hinge and check the function of the **magnet**		
10.2	Check the charger **output** to the hinged rail part		
10.3	Check the **collision bracket** (on wear and thightness and straightness)		
10.4	Check that the **safety** is functioning (the lift must stop when it hits the hinged rail)		

	***11. OPTIONS >** SLIDE TRACK*	✗	✓
11.0	Check the function of the **Slide Track**		
11.1	Check that the **safety edges**(top and bottom) are functioning (the track must stop when it hits something)		
11.2	Check that the track **stops** at the right position (top and bottom). Adjust the magnet positions if necessary.		
11.3	**Lubricate** the underside of the Slide Track rail near the bottom foot (described in the installation manual)		
11.4	Check that the **safety edges** (top and bottom) are intact and securely fitted to the ends of the rail.		
11.5	Check that all **roller foot covers** are present and are securely fitted to the roller feet.		
11.6	Check that **warning labels** are intact and securely fitted to the first and last roller foot covers.		
11.7	Check if the gas strut is leaking (significant amount of oil in the area of the first (fixed) mounting foot).		
11.8	Check that the **battery box** is intact and securely fitted to the mounting bracket		

	12. FINISH	✗	✓
12.0	Check that all removed parts are **reinstalled** proper		
12.1	Make again a **test ride** and check if all functions are still okay		
12.2	Thank the customer for their time and **mention** them that they can call the dealer if there are any problems		

The stairlift is safe / unsafe	Remarks
What repairs, renewals or alterations are required to enable the stairlift to continue to be used safely?	
State specified time limit for these repairs:	Date of examination:
Signature Handicare lifts employee:	Signature customer:

DEALER IS RESPONSIBLE FOR STORAGE (AT LEAST 5 YEARS) OF THIS CERTIFICATE [FILLED IN AND SIGNED WHEN TESTED]

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


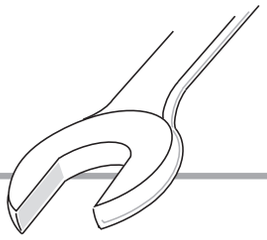
A UNINSTALL LIFT


DOWNWARDS

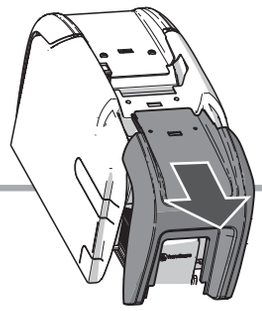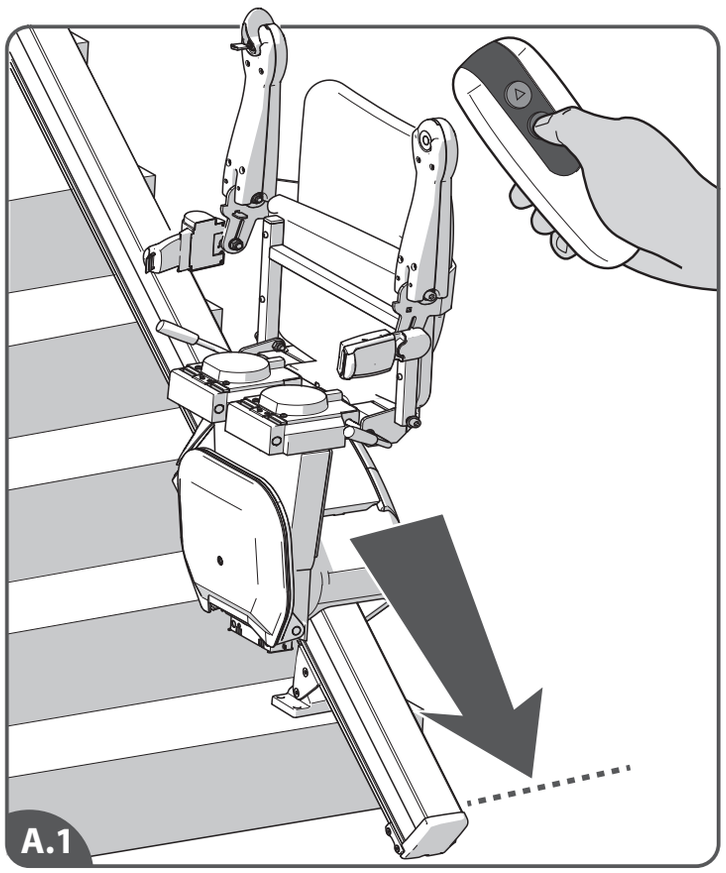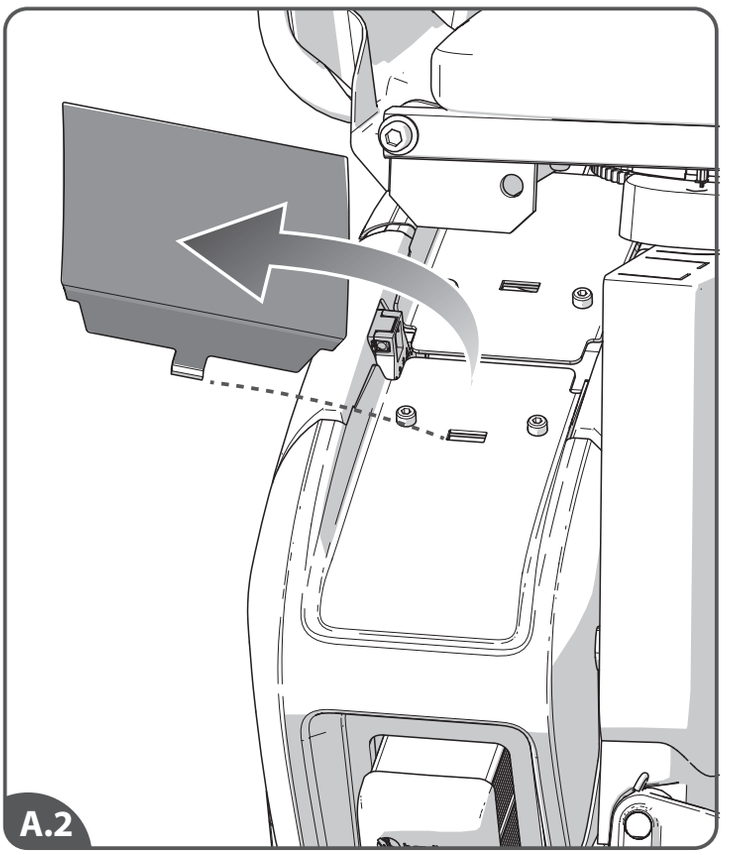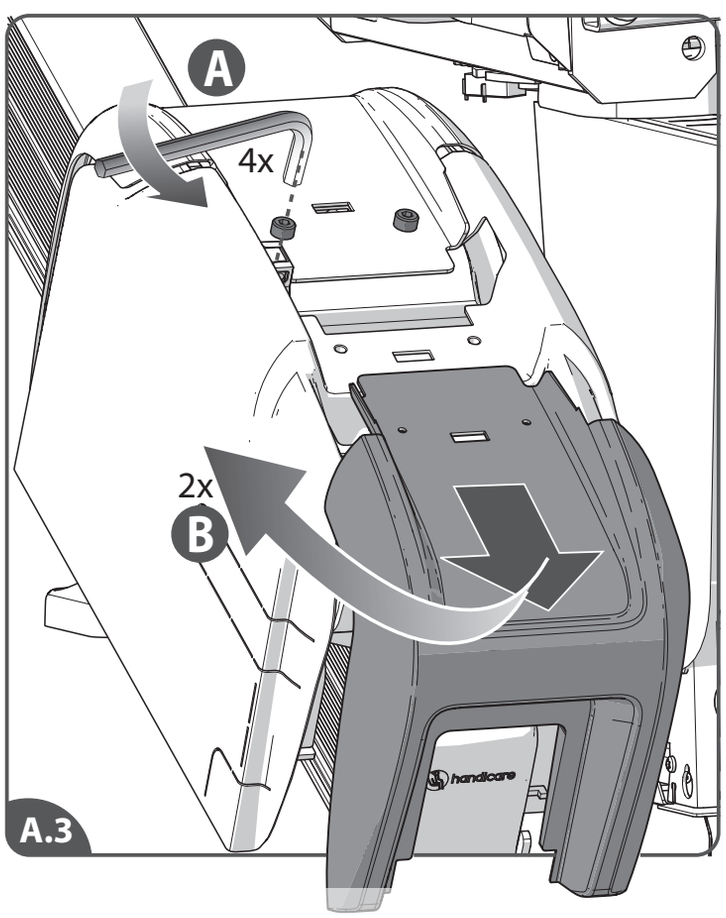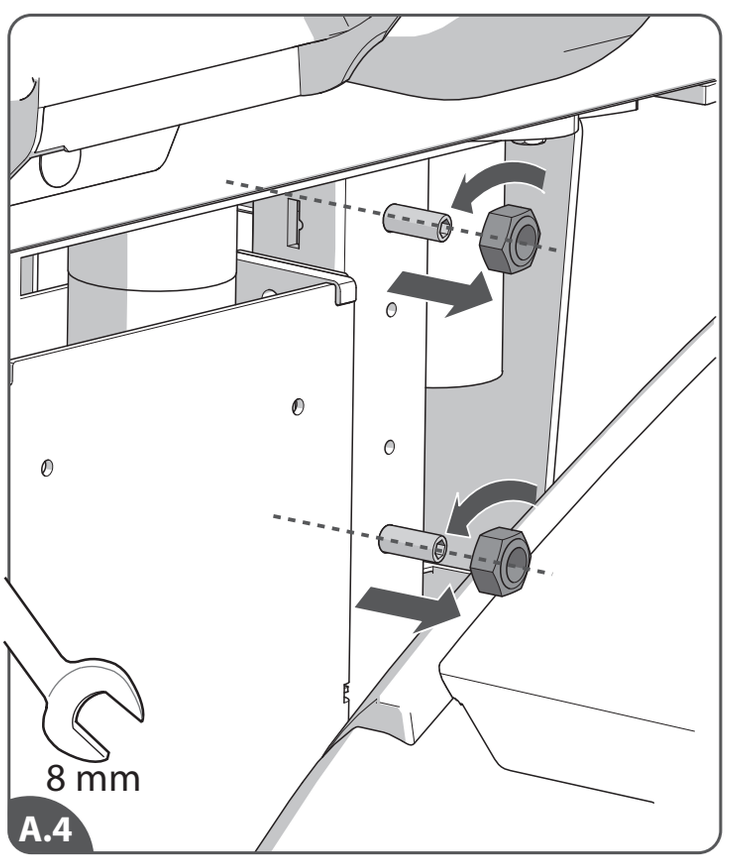[A.1]Ride the lift down.
[A.2]Remove the diagnostics display cover.
[A.3]Unscrew the safety guards{A}. Slide the safety guard over the rail and remove from the power pack{B}.
[A.4]Unscrew the two nuts.

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



IT


A UNINSTALL LIFT COVERS

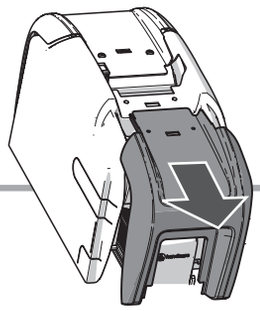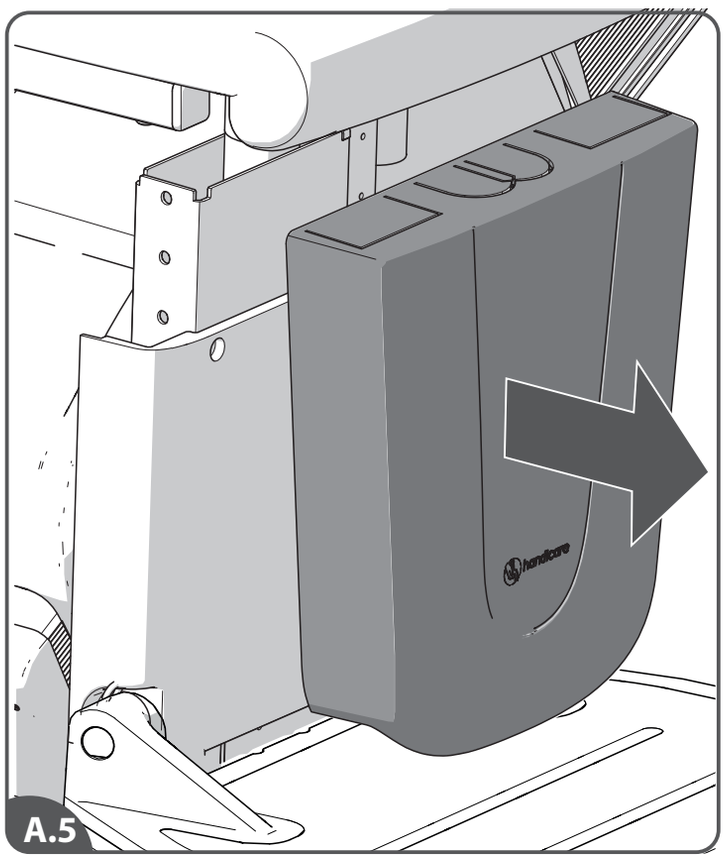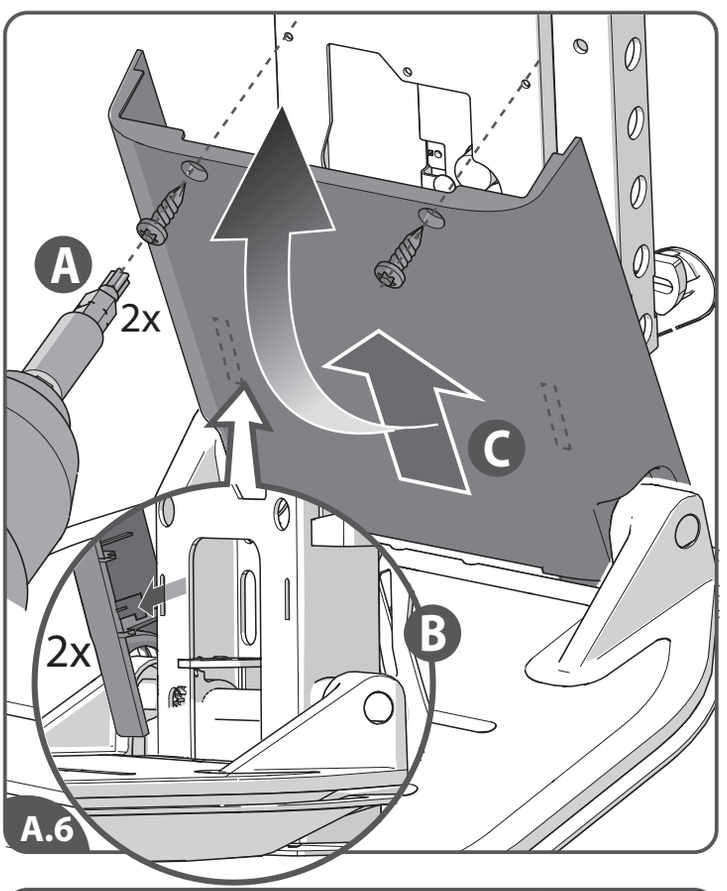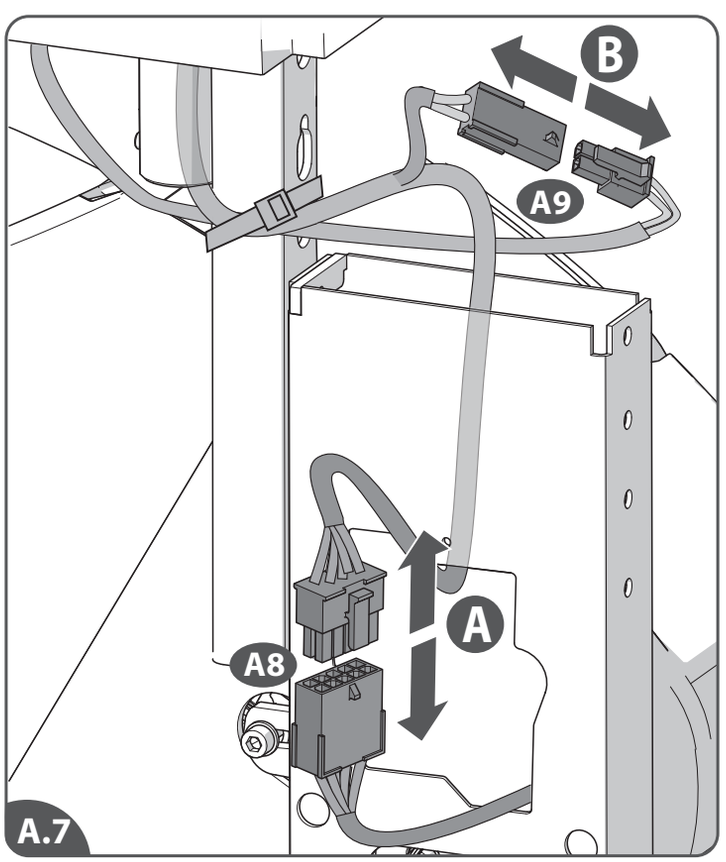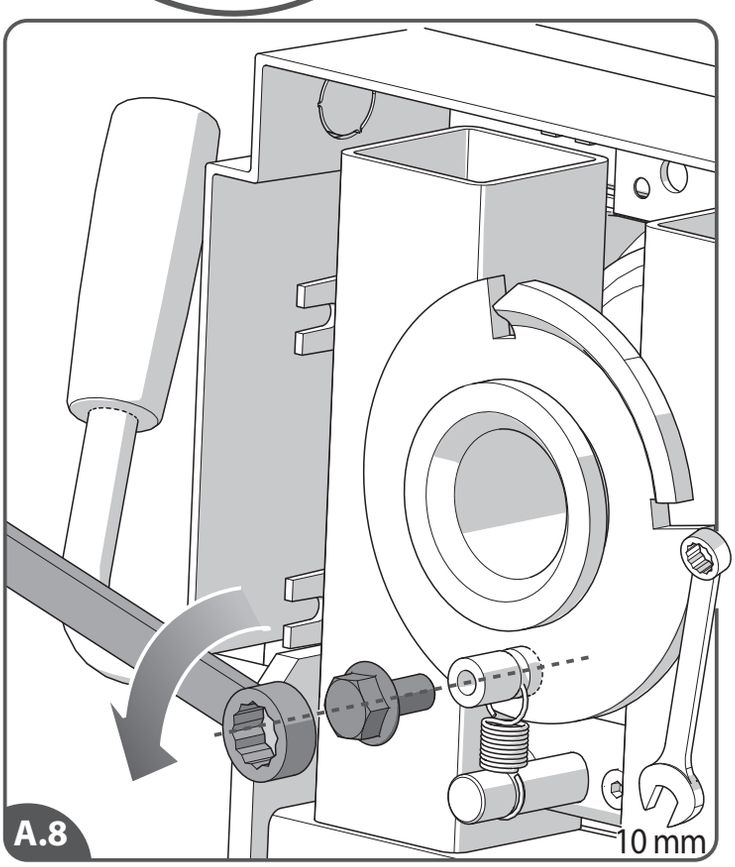[A.5]Remove the upper tilt chassis cover.
[A.6]Unscrew the bottom tilt chassis cover {A}. Remove the bottom tilt chassis cover {C}. Slide the two ribs on the cover out of the recesses on the frame {B}.
[A.7]Disconnect all other connections of the seat looms.
[A.8]Remove the securing bolt from the bottom of the seat.

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



IT


A UNINSTALL LIFT SEAT

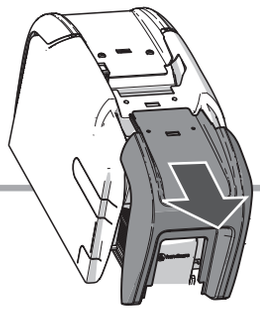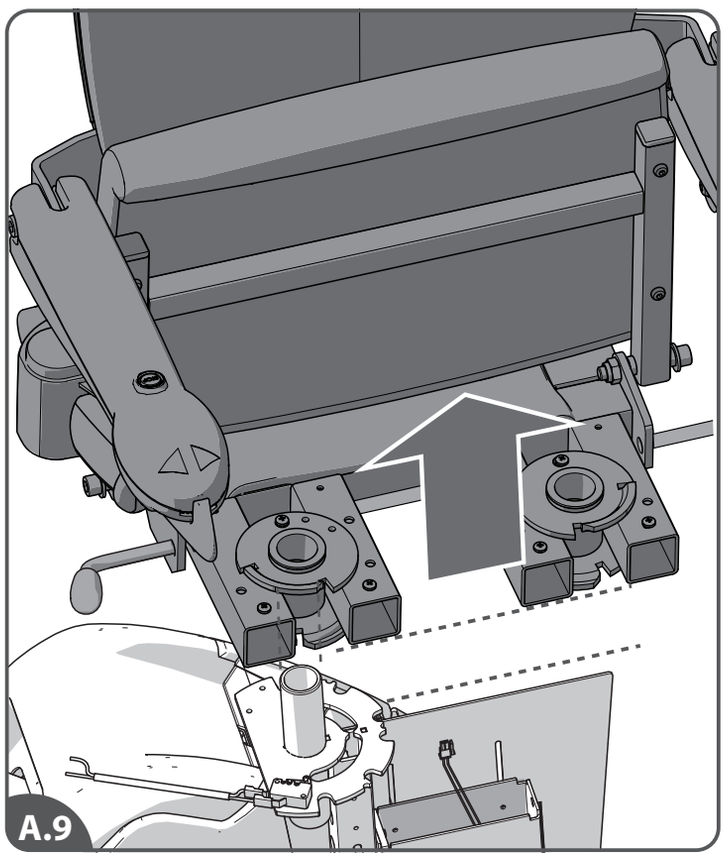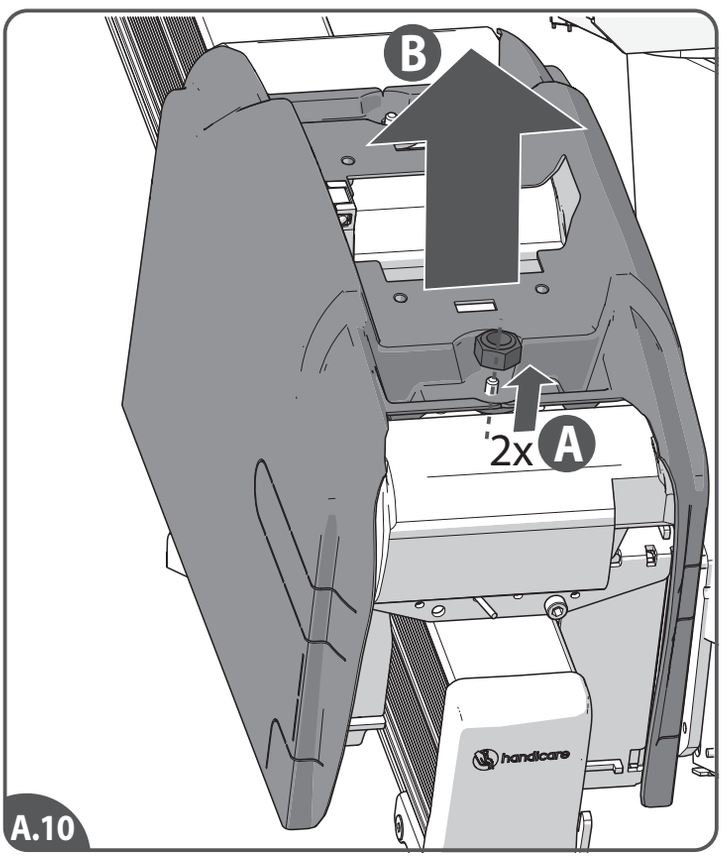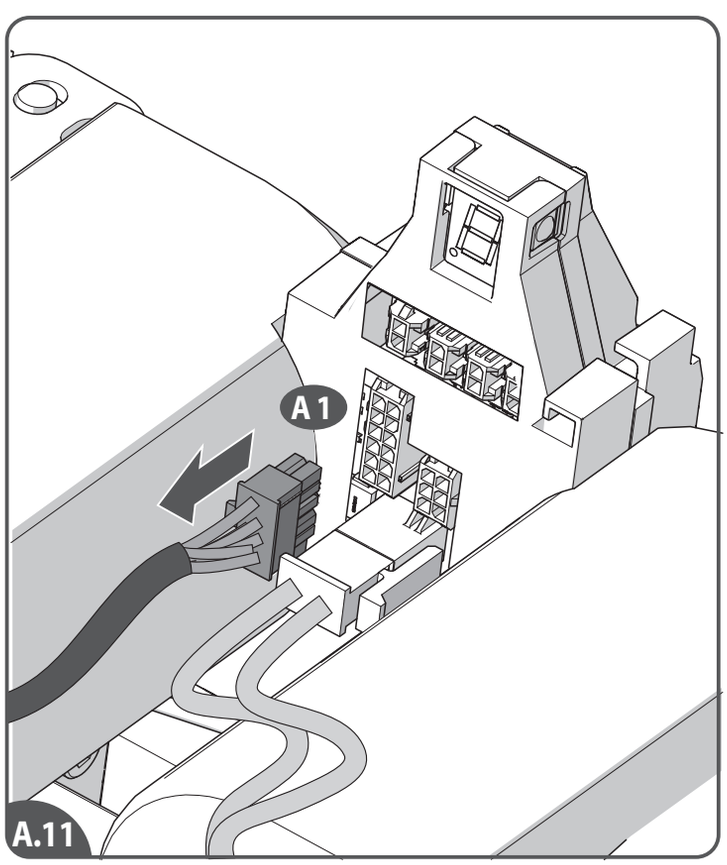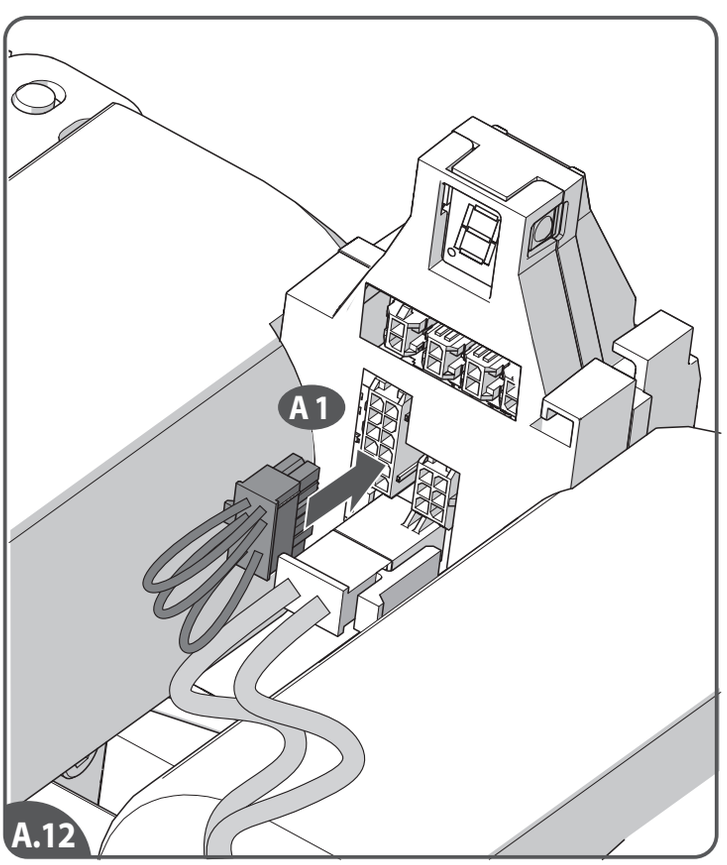[A.9]Remove the seat and lay it down in a safe place.
[A.10]Remove the two nyloc nuts{A}.Remove power pack cover from the unit{B}.
[A.11]Disconnect the main loom from the control box.
[A.12]Connect the shortcut connector to the PCB.

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
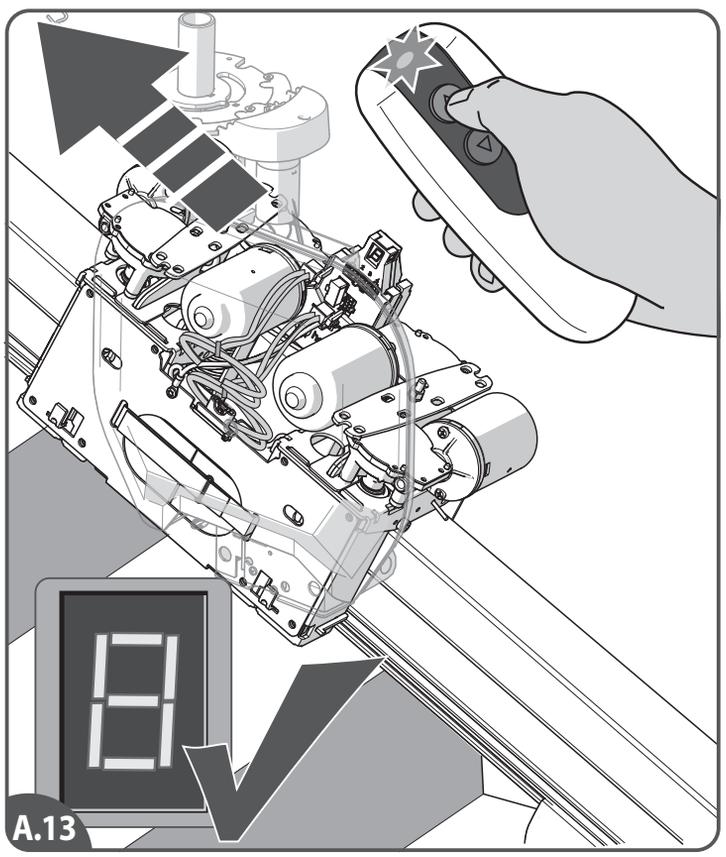
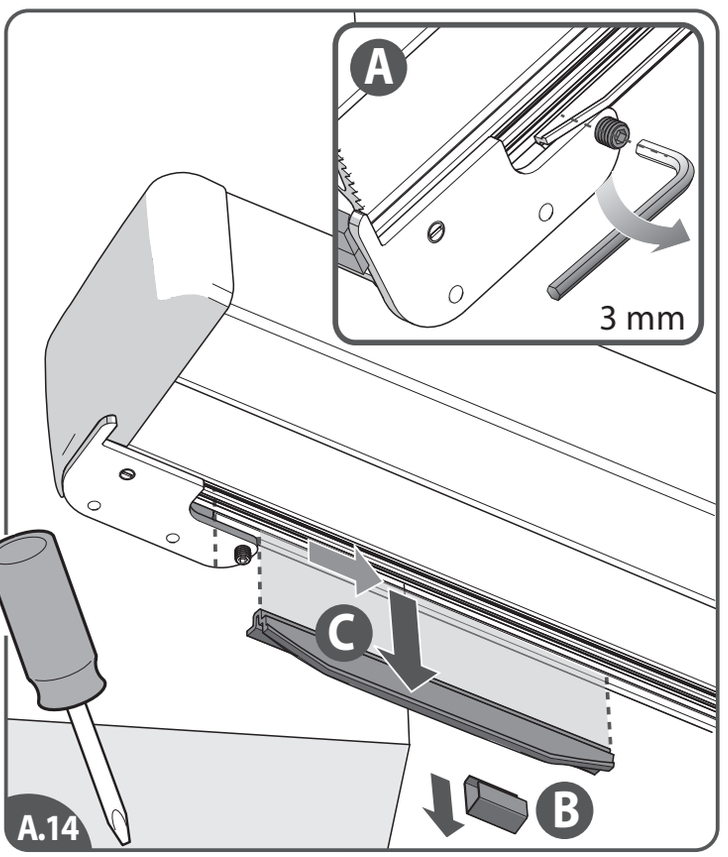
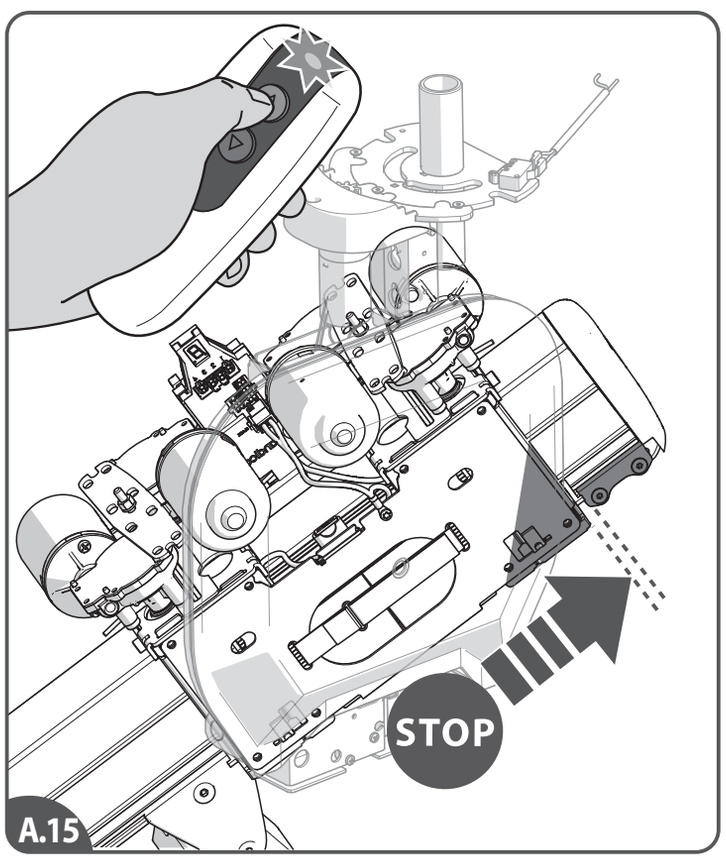
STOP
A.15

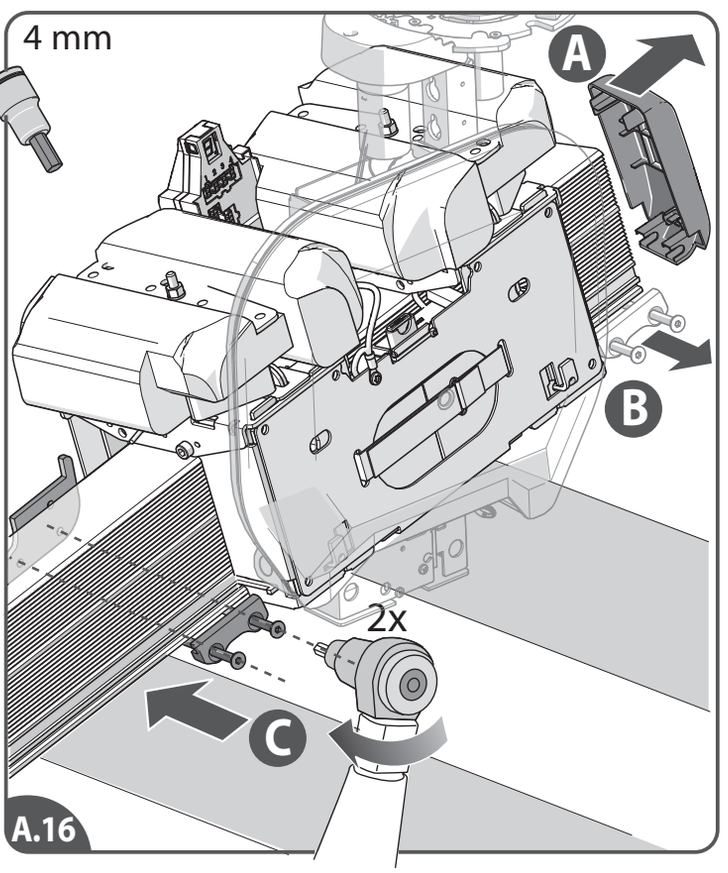4 mm
A
B
2x
C
A.16

[A.13]Press the remote control to drive the power pack up the rail.
[A.14]Remove the limit actuator and the circlip on top of the stairs.
[A.15]Ride the lift to the top of the stairs just before touching the mechanical endstop.
[A.16]Remove the end cap{A}. Remove the mechanical endstop assembly{B}. Install the mechanical endstop assembly just below the power pack (to stop the power pack when the tension is released) and fasten the bolts {C}.

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



IT

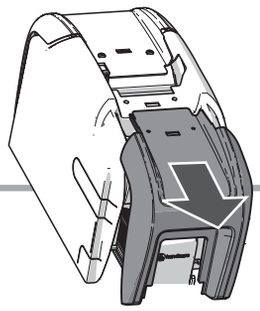UNINSTALL LIFT UNIT

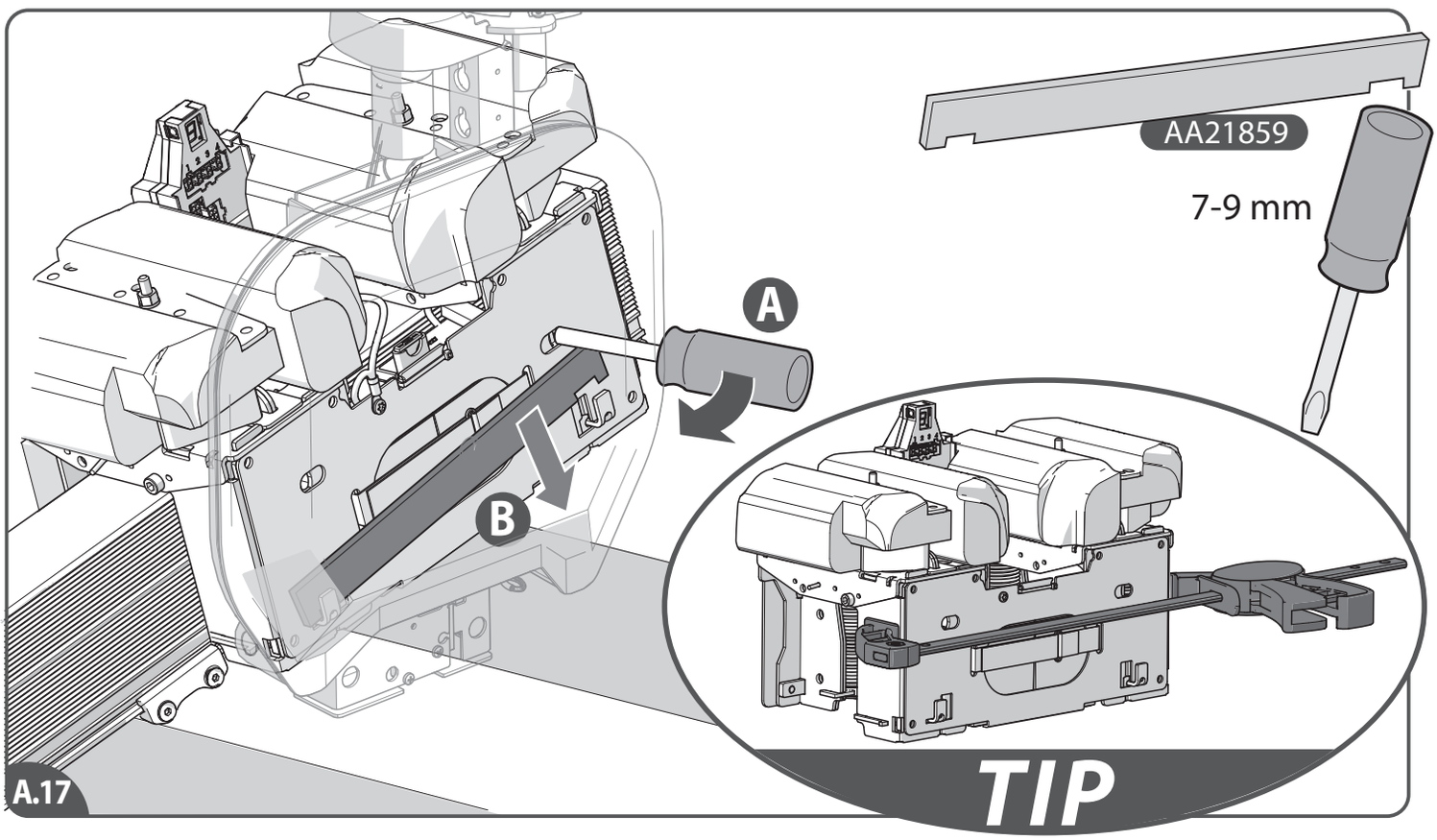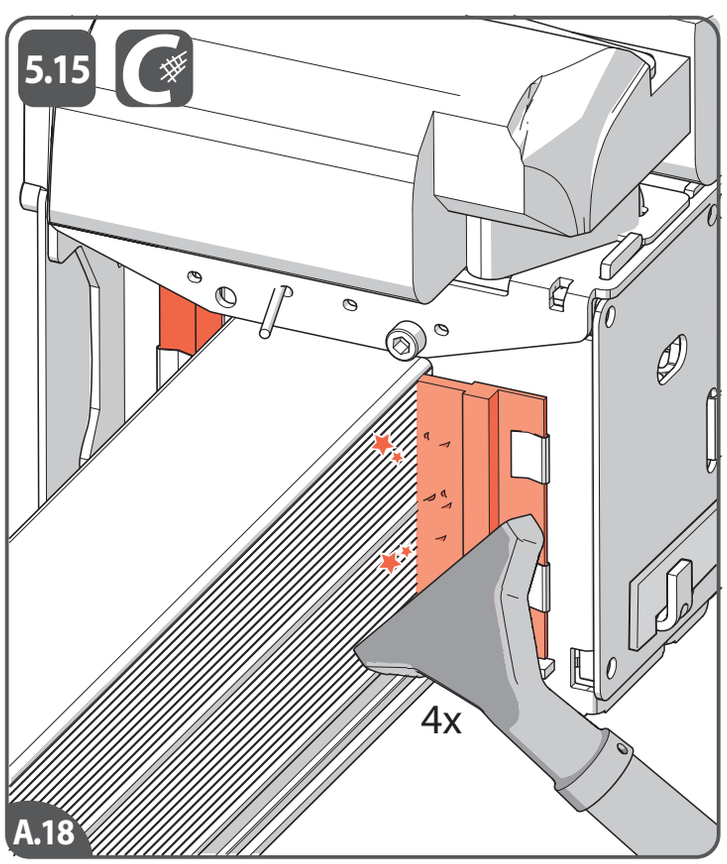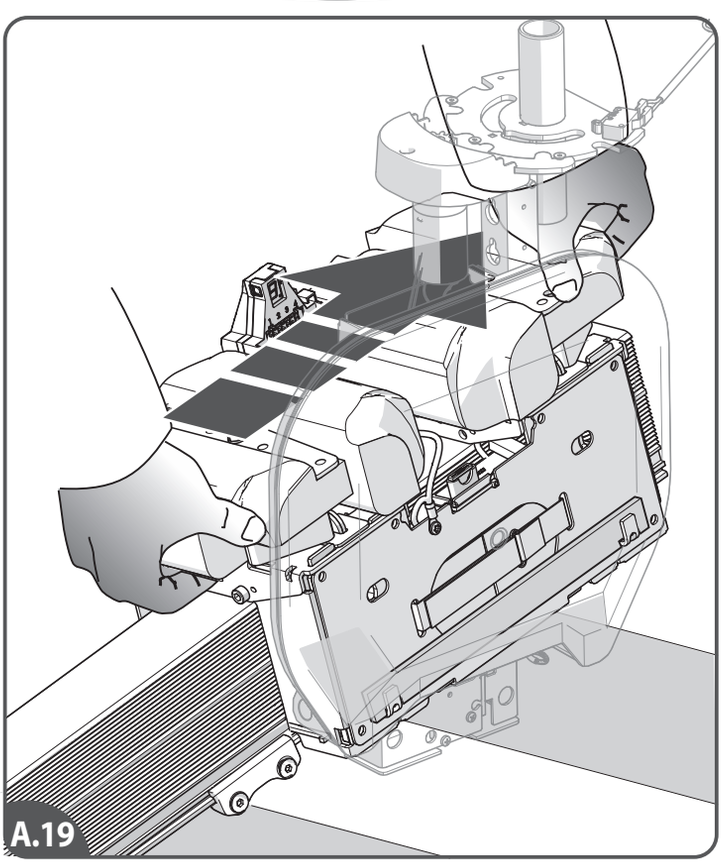[A.17]Use a flat screw driver and push the upper wheel set inwards{A}. Take away the tension by placing the bar{B}.
[A.18]Vacuumclean the brusches on the power pack.
[A.19]Remove the power pack and place it on a mat.

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

B REMOUNT TILT CHASSIS/STEM

FOOTREST_CHASSIS


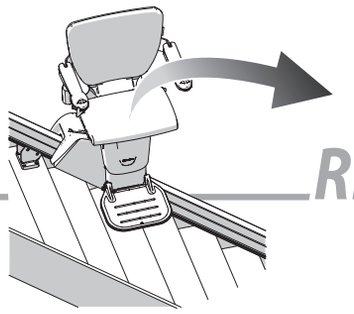
NOT ALLOWED AT STANDARD SERVICE WORK.


NECESSARY KIT IN CASE THE TILT CHASSIS OR SEAT STEM WILL BE REMOVED !

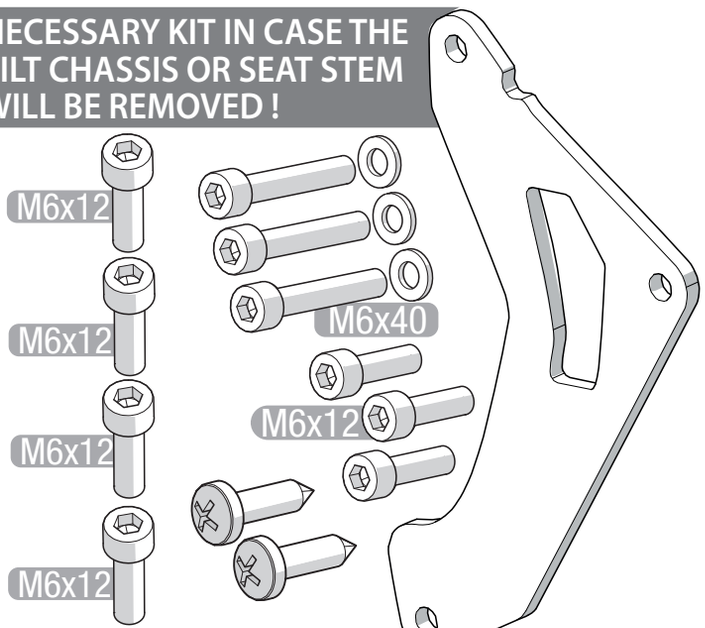REMOUNT TILT CHASSIS KIT:

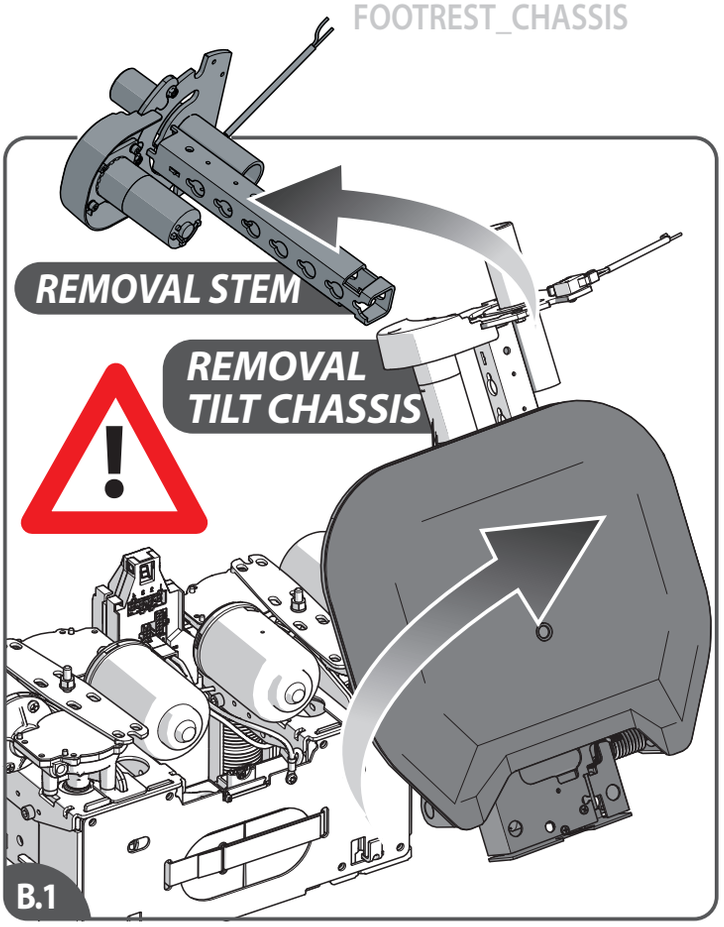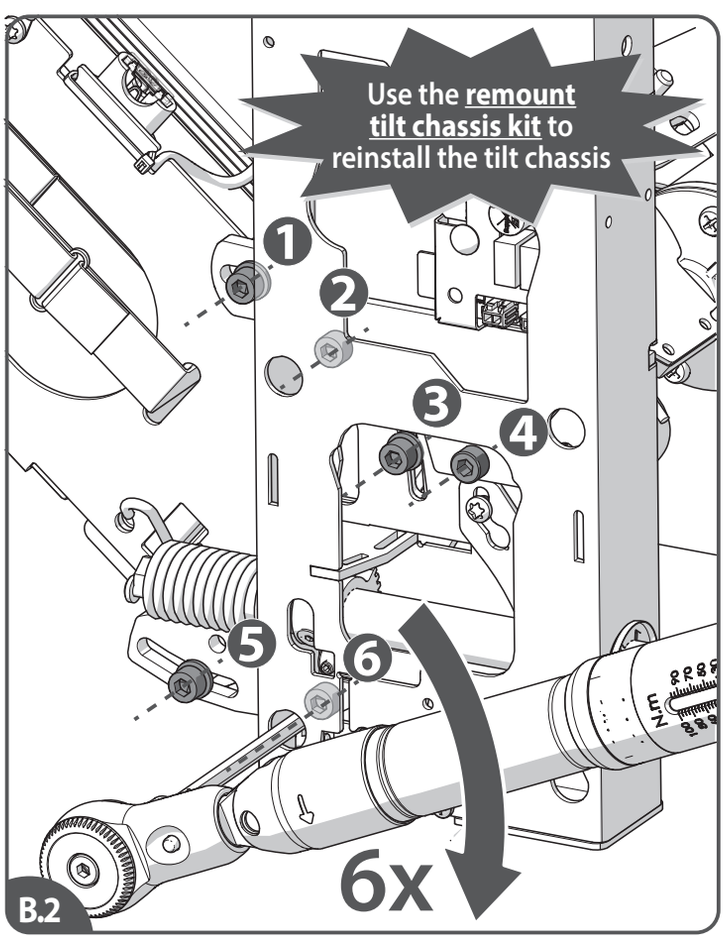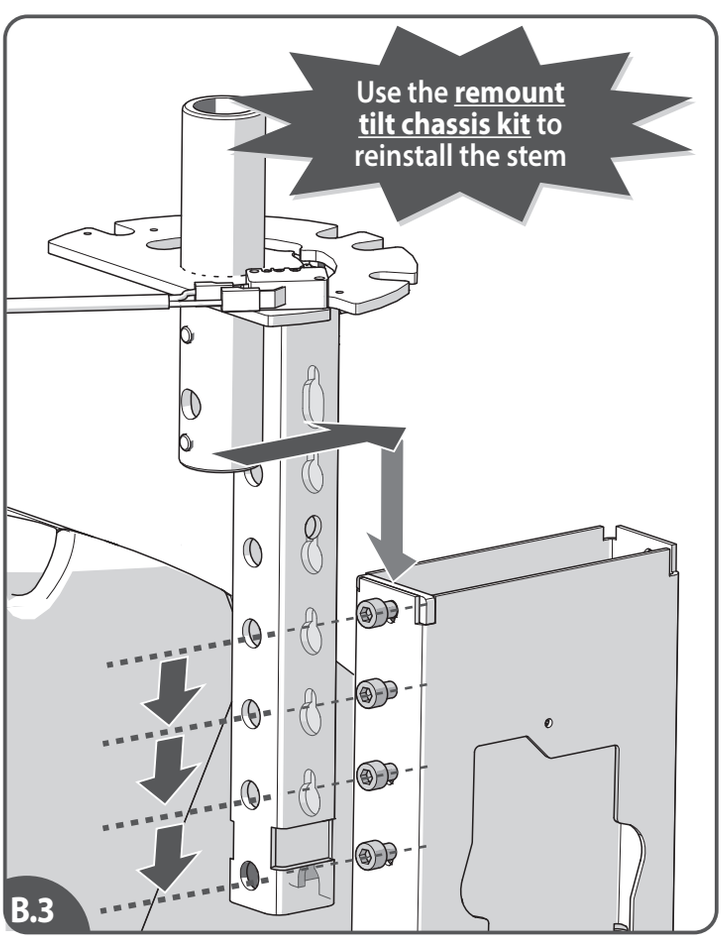NOT ALLOWED AT STANDARD SERVICE WORK !
[B.1]In case of removal of the tilt chassis or the seat stem renew the bolts and the adaptorplate.
[B.2]Use the Remount tilt chassis kit to reinstall the chassis.
[B.3]Use the Remount tilt chassis kit to reinstall the stem.

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



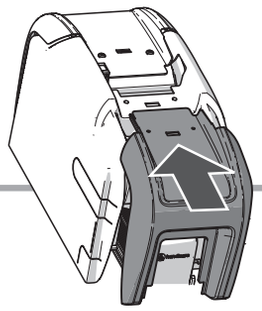
C REINSTALL LIFT

POWER PACK

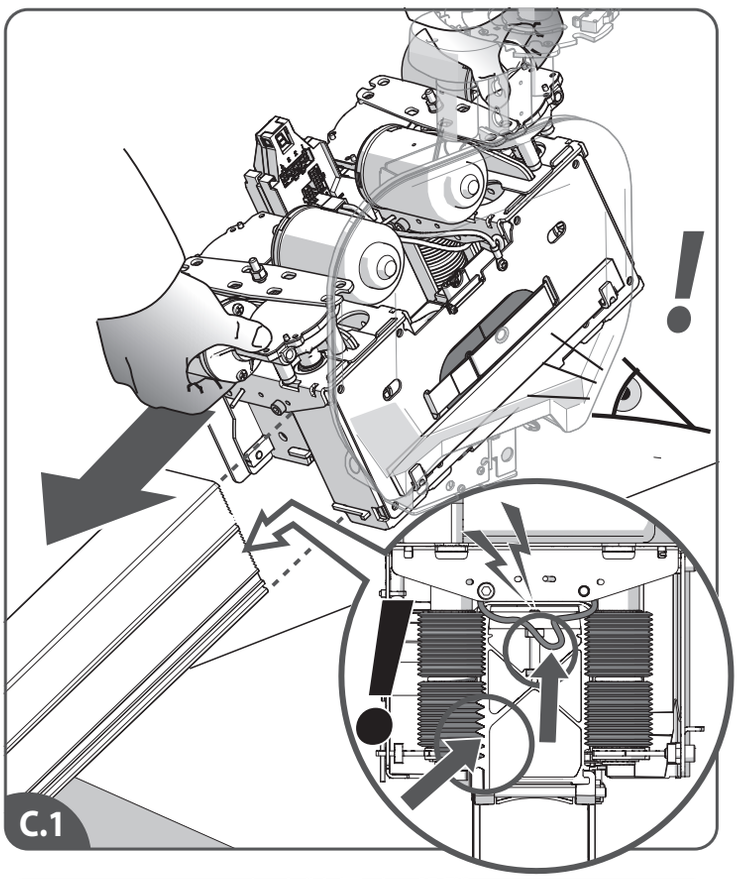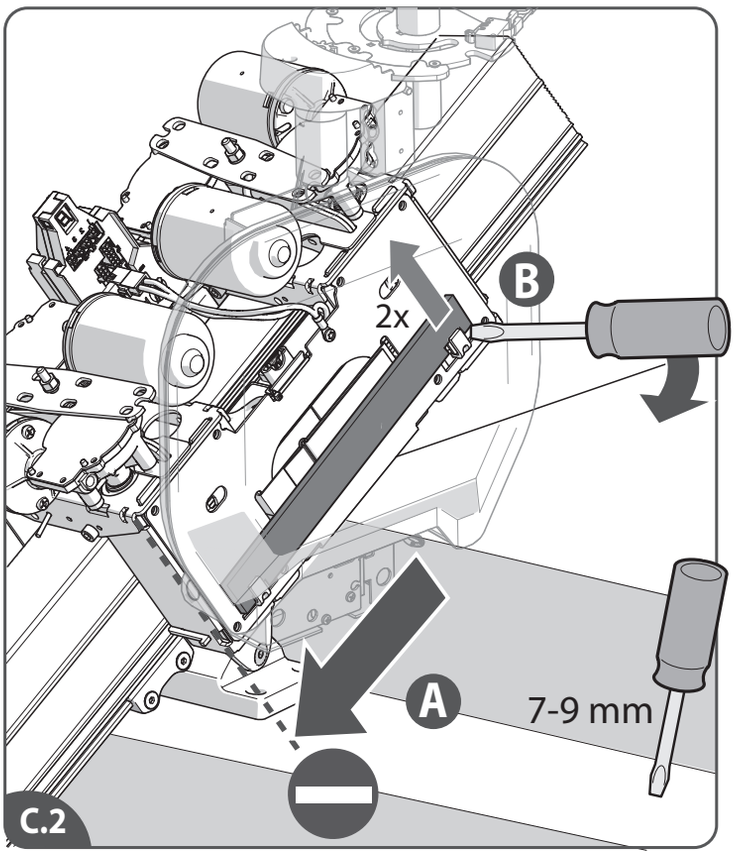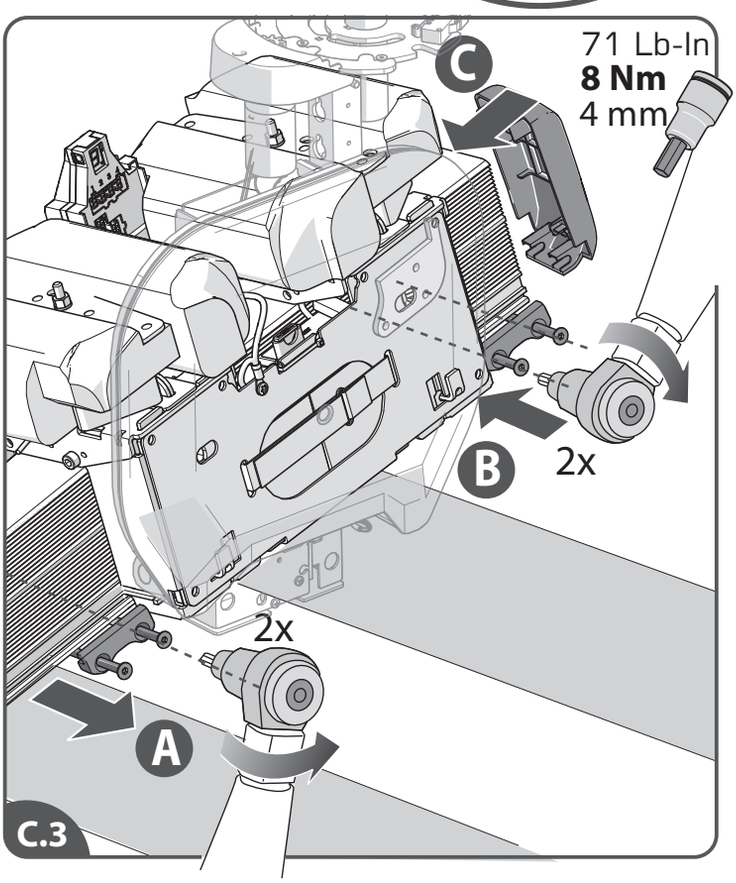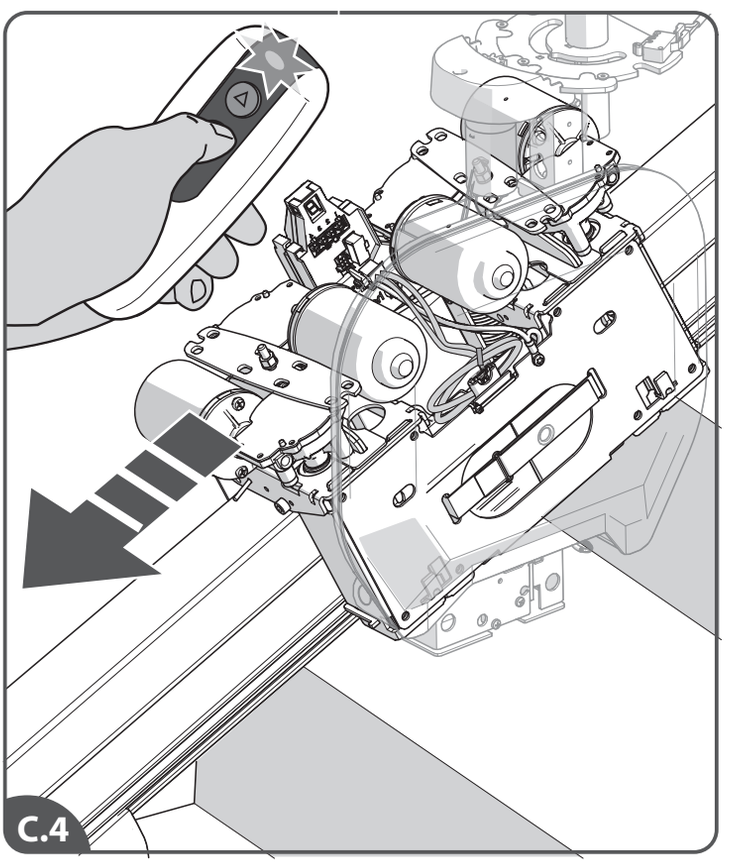[C.1]Reinstall the power pack on the rail. Ensure that the motorised rollers engage properly with the rail and the looms do not pinch between the rail and the unit.
[C.2]Lower the power pack down the rail until the power pack hits the endstop{A}. Remove the bar, so the power pack clamps the rail {B}.
[C.3]Reinstall the endstop to the rail (align precisely with the rail end) {A}. Tighten the bolts securely {B} {8 Nm/71 lb-in}. Fit the end cap {C}.
[C.4]Ride the power pack down.

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




IT

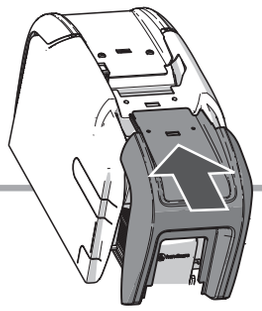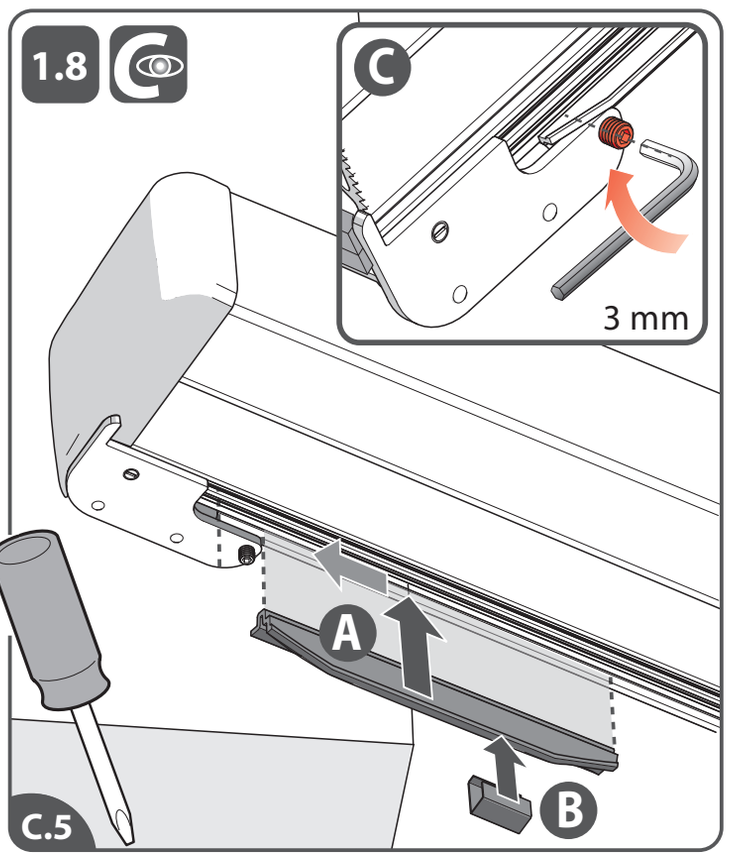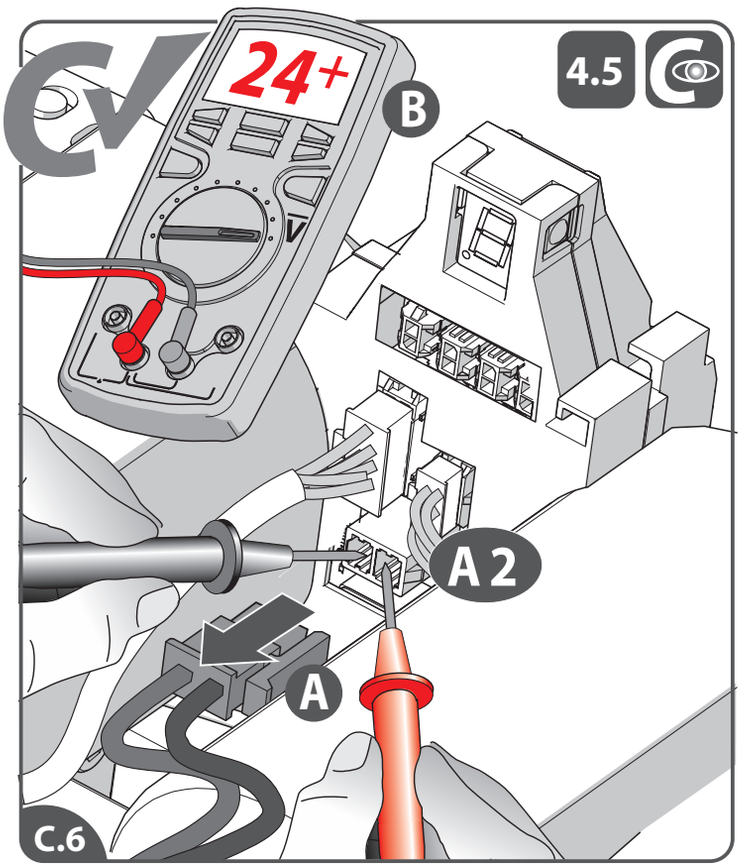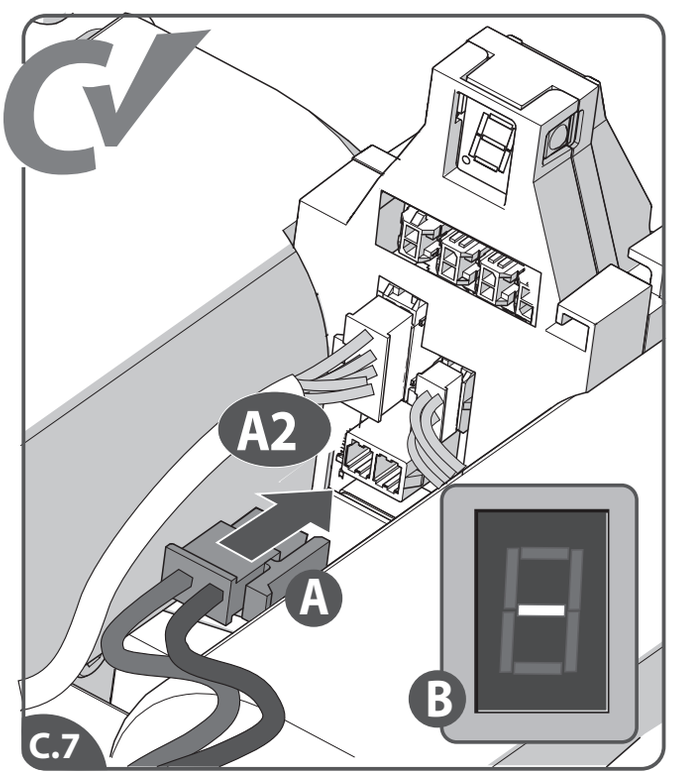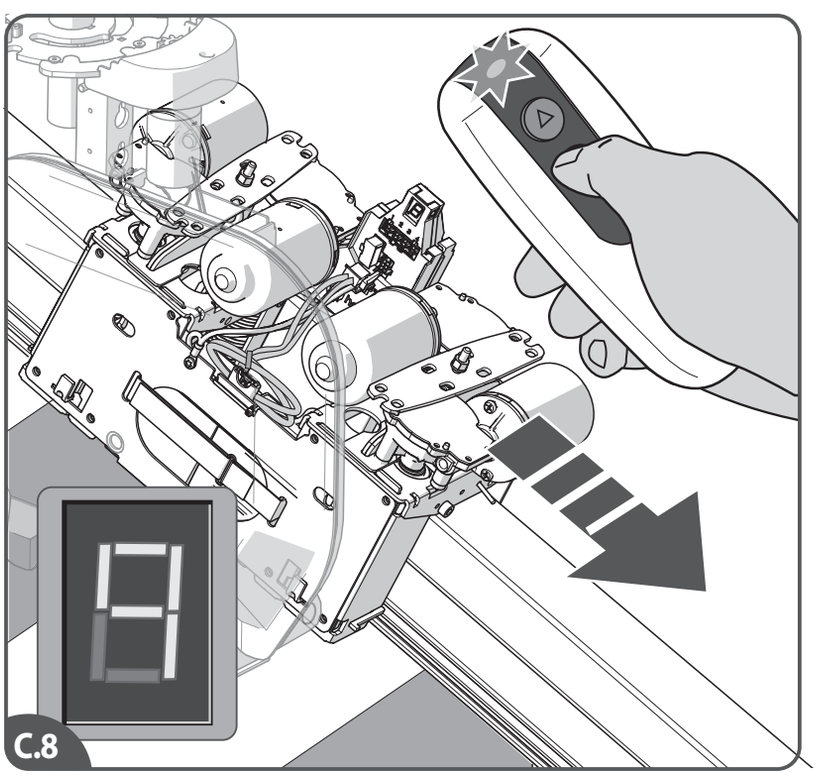SEE THE INSTALLATION MANUAL FOR THE FURTHER REINSTALLATION OF THE LIFT

[C.5]Slide the limit actuator into the rail, up against the end stop (slanting sides at the top!) {A}. Fit the circlip over the limit actuator{B}. Secure the adjusting screw in the end stop[C].
[C.6]Remove the battery loom connector from the control box(red/black wires){A} and use a Volt meter to measure voltage (should be 24+ V){B}.
[C.7]Reconnect the battery loom{A} and check that the lift is charging{B}.
[C.8]Press the remote control to drive the power pack all the way down.
SEE THE INSTALLATION MANUAL FOR THE FURTHER REINSTALLATION OF THE LIFT

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


IT

REQUIRED REPLACEMENT OF THE UNIT

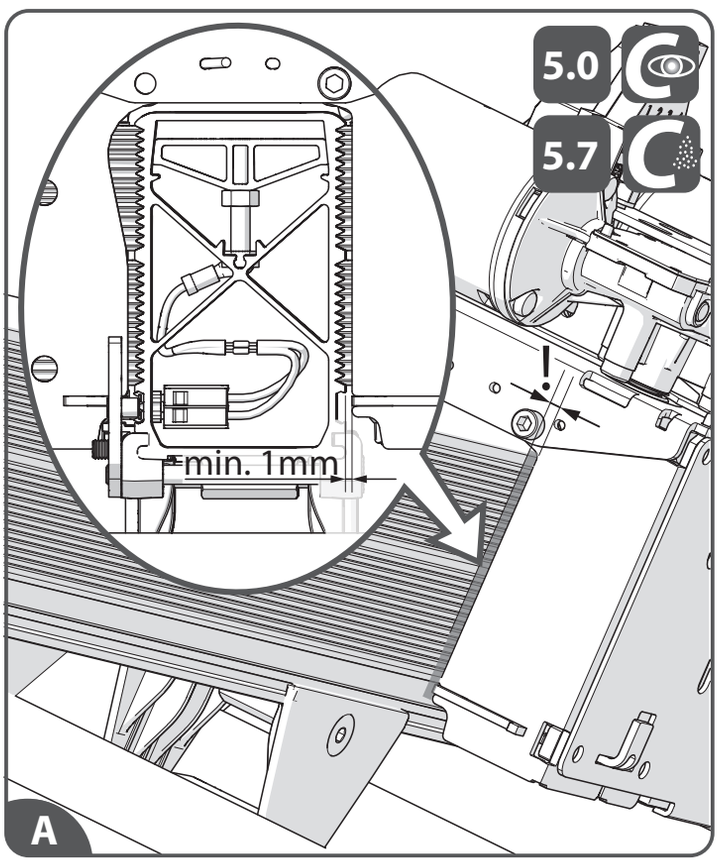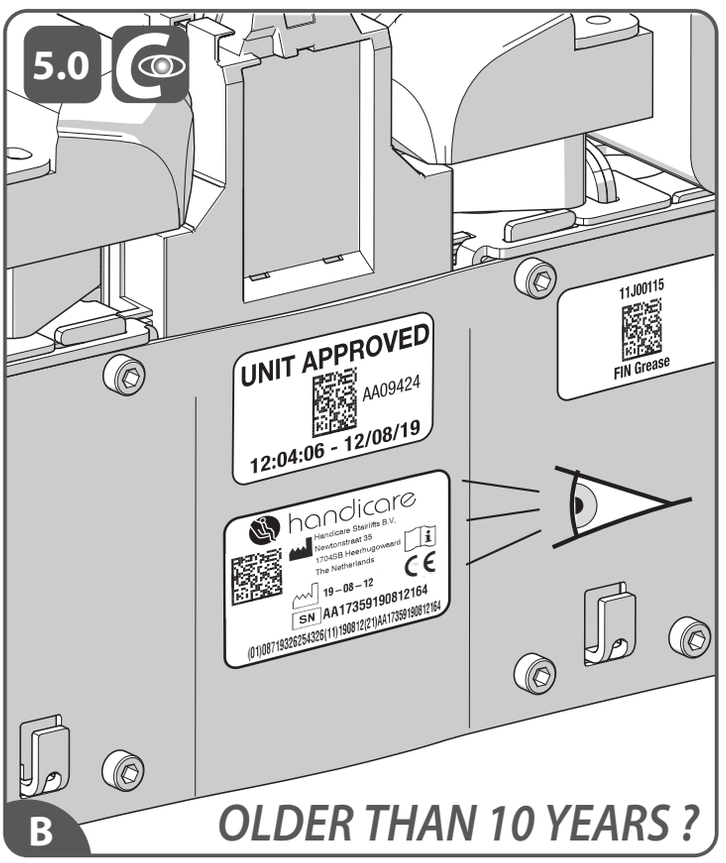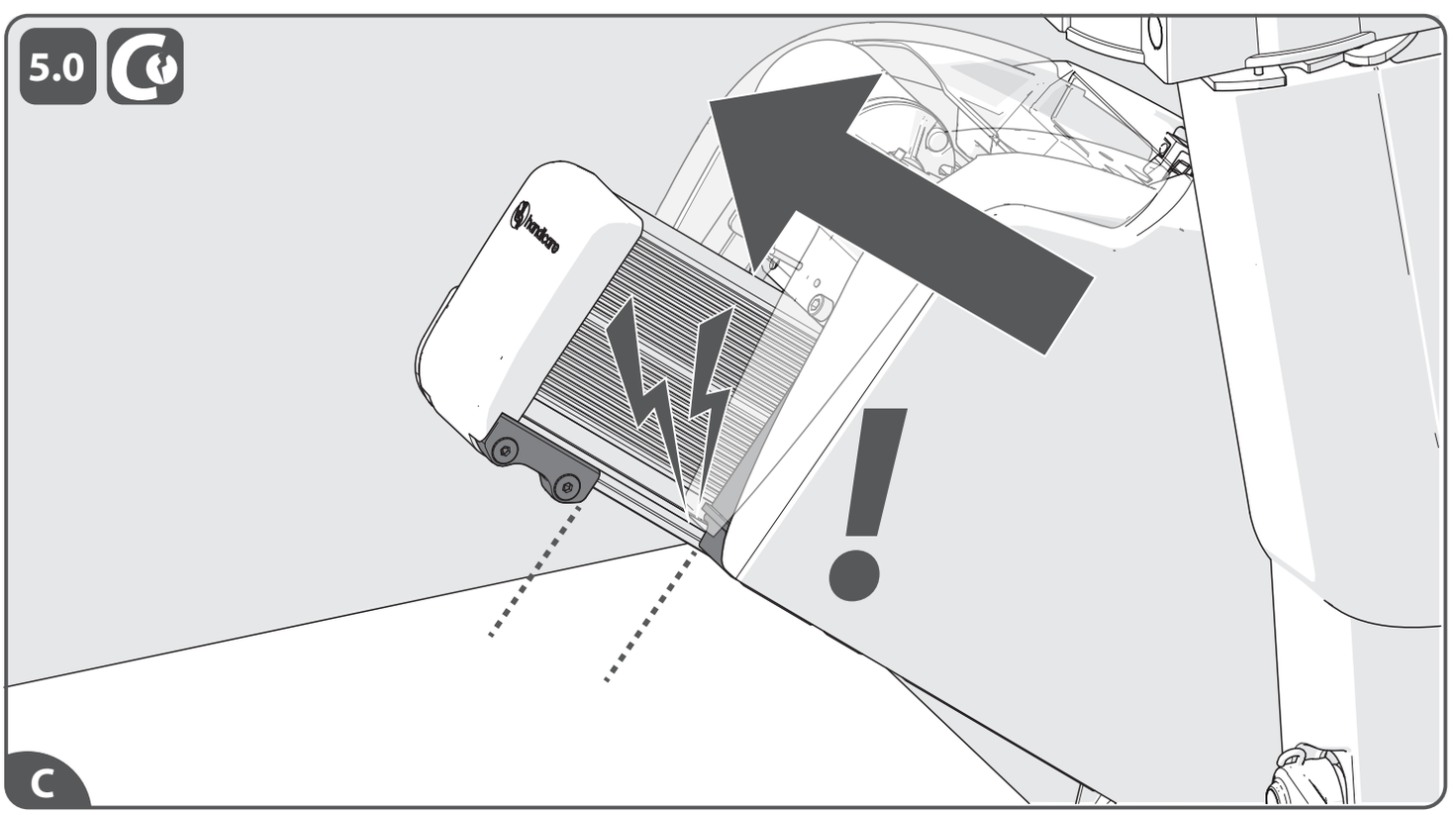Replace complete 1100 unit in case of:
[A]If due to wear of the rollers the horizontal distance between rails and steel front plate of the power pack at position of the lowest (roller which bears most load) roller becomes less than < 1mm.
[B]If power packis more > 10 years old.
[C]If the lift has run on the mechanical endstop.

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
LABEL ON CHARGE STRIP

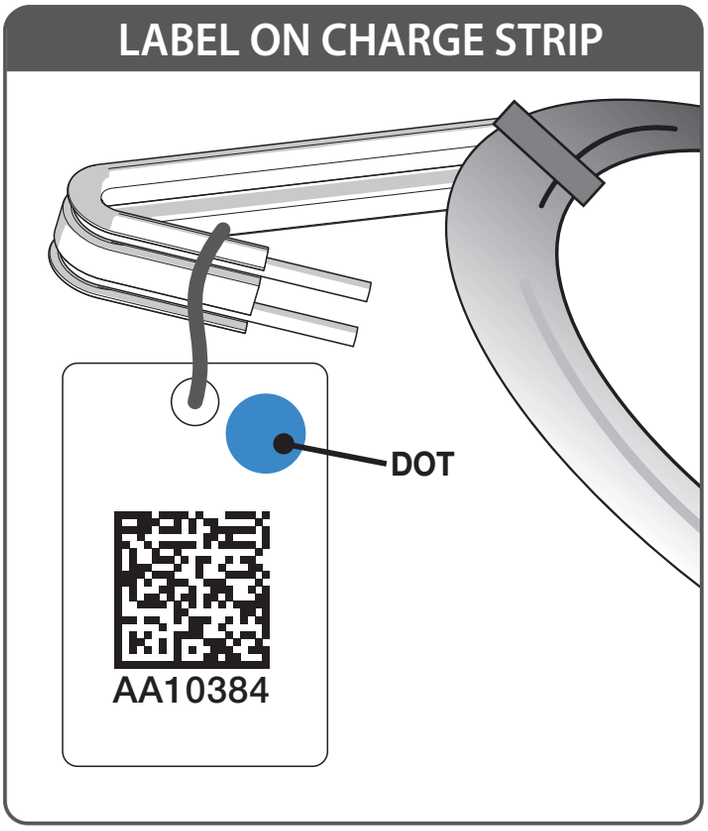STICKER ON FITTING KIT

Label on charge strip	Sticker on fitting kit	Matches with rail
No dot	Rail fitting kit	Batch 1xx
No dot	Rail fitting kit	Batch 1-5
Green dot	Rail fitting kit - CCP	Batch 1-7
Yellow dot	Rail fitting kit - CCP-Y	Batch 1-10
Blue dot	Rail fitting kit - CCP-B	Batch 1-11

EN

The dimensions of this slot may vary between different rail batches. You can find the batch number of the rail on the top of the box. In a single delivery, the charge strips match the delivered batch of rails. If the charge strips and rails are seperated you can check if a charge strip is suitable for a particular rail thanks to a sticker on its label. There is also further guidance on the fitting kit's sticker.Please see below an overview of the different charge strips and their labelling.
By the end of July 2020 we will change the slots of the rail in such a way that matching of fitting kits and rails is no longer necessary. We ensure that sliding of the charging strip into these rail slots will be far more easy. These rails will be marked as "batch 100".

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
CHARGING PROBLEMS

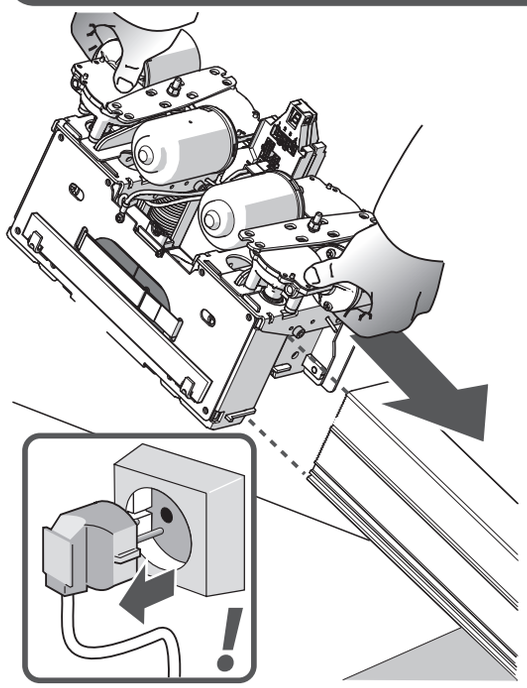Check that the charger is unplugged before reinstalling the power pack.

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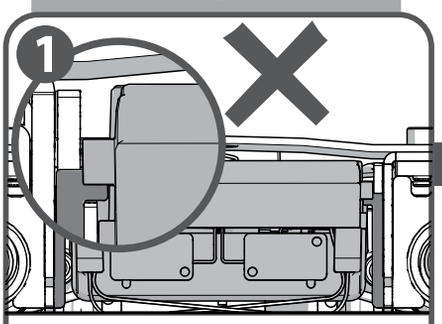
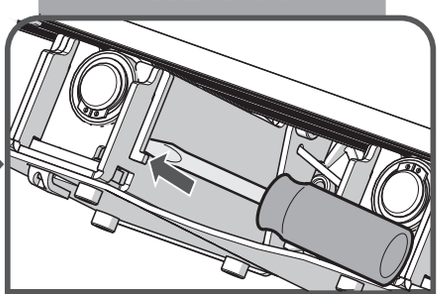
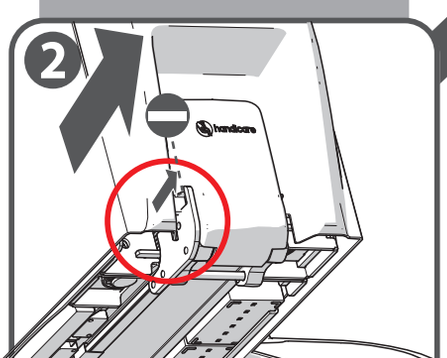
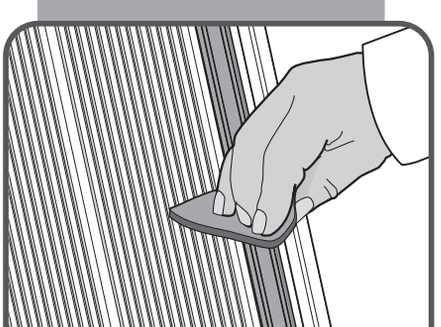
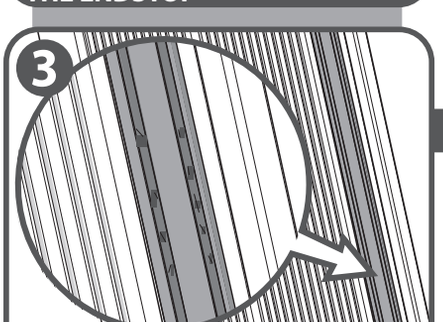
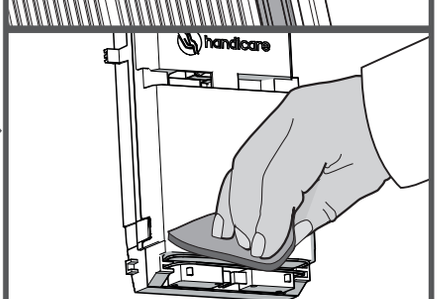
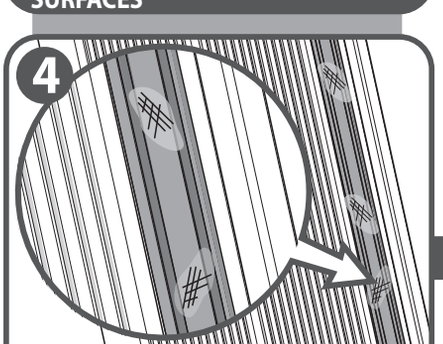
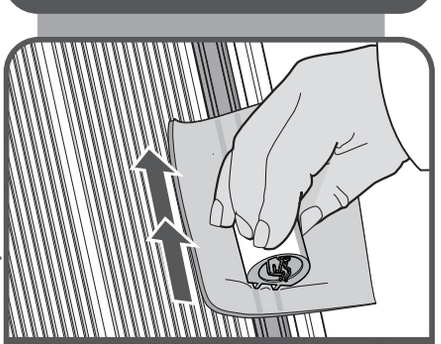
CHECK !	REPAIR
1  THE PCB IS NOT FITTING WELL/ STICKING IN HIS BRACKETS	 STRAIGHTEN THE BRACKETS OF THE PCB BY BENDING THEM A BIT
2  OVERRUNNING: LIFT STOPS BECAUSE UNIT CAP IS HITTING THE ENDSTOP	
3  CORROSION ON THE CHARGING SURFACES	 CLEAN THE CHARGING SURFACES WITH A SCOTCH-BRITE HANDPAD
4  POLLUTION ON THE CHARGING SURFACES	 CLEAN THE STRIP WITH A TISSUE WITH A DEGREASING LIQUID (USE KIT AA23773)

CHARGING PROBLEMS
[1]The PCB is not fitting well/sticking in his brackets > straighten the brackets of the pcb by bending them a bit.
[2]Overrunning: lift stops because unit cap is hitting the endstop > straighten the brackets of the PCB by bending them a bit.
[3]Corrosion on the charging surfaces > clean the charging surfaces with a scotch-brite handpad.
[4]Pollution on the charging surfaces > clean the strip with a tissuewith a degreasing liquid (use kit AA23773).

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
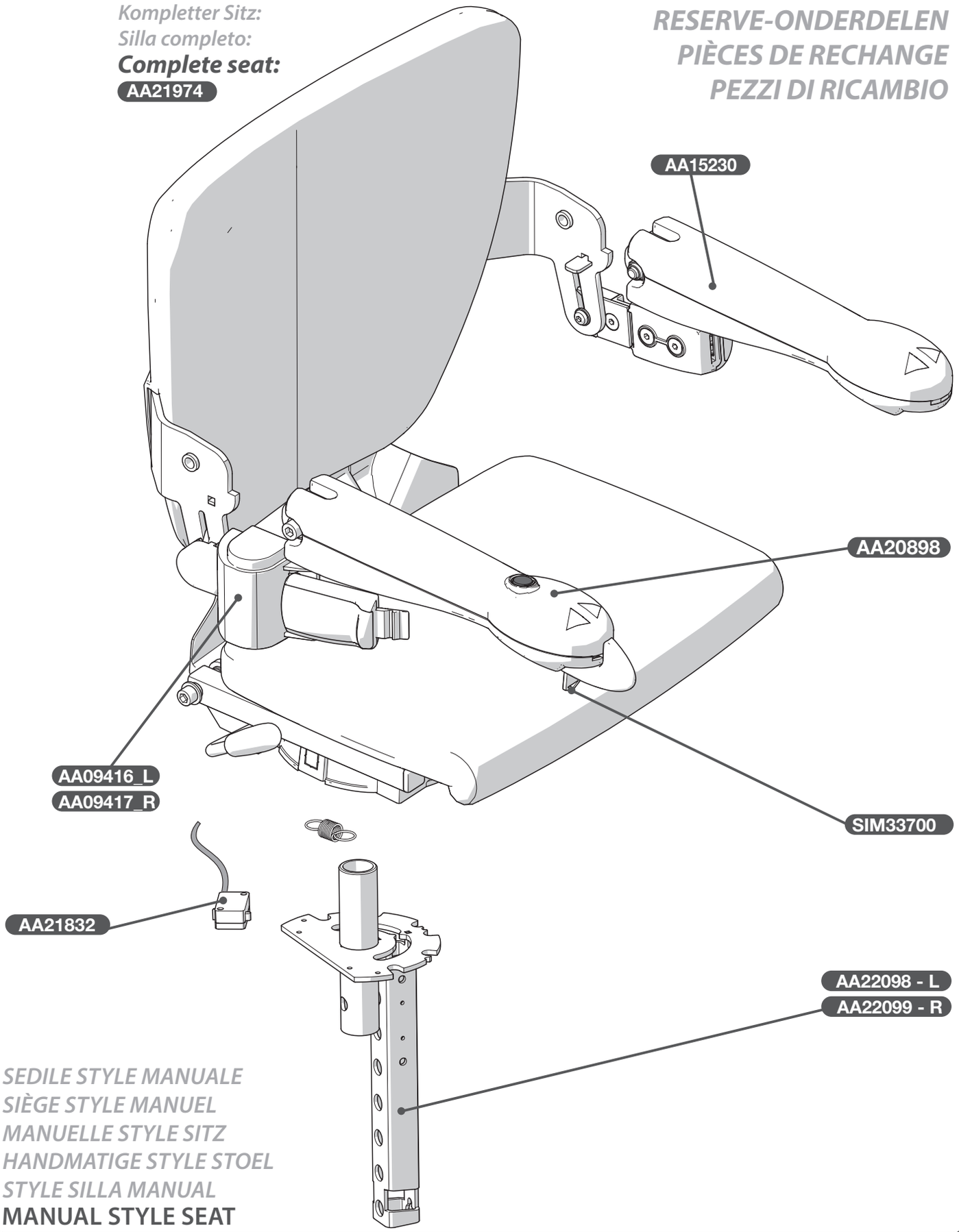
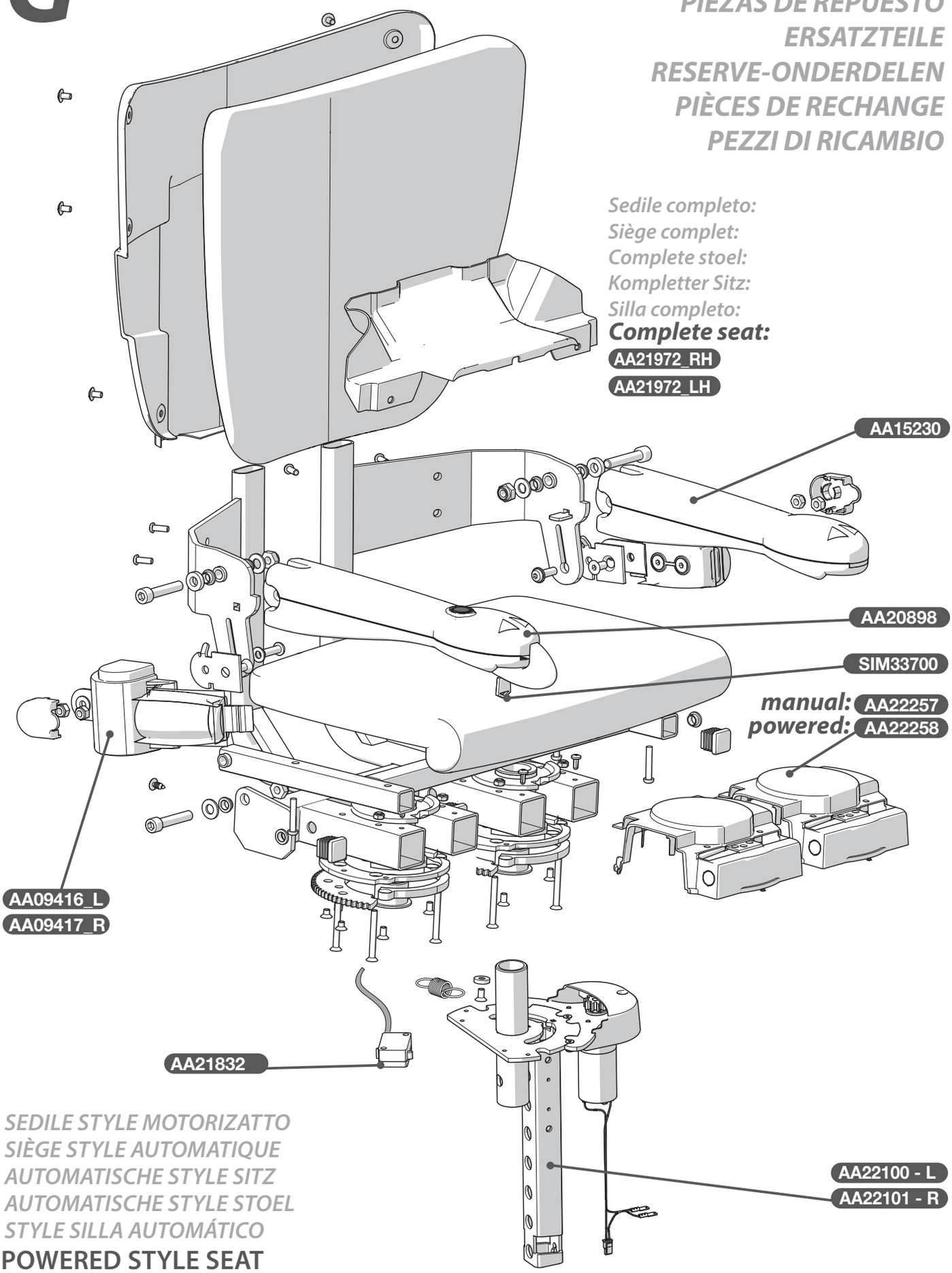
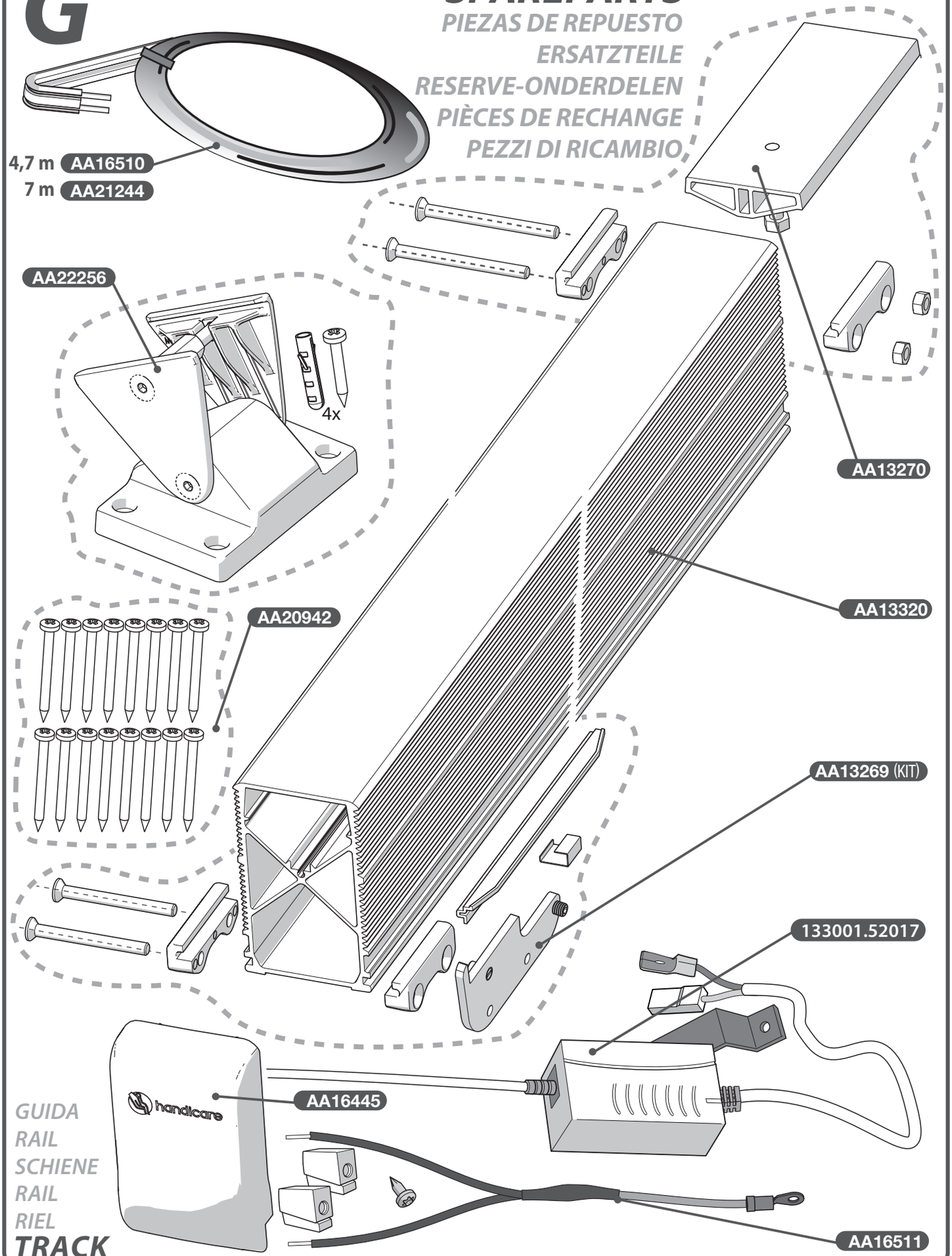
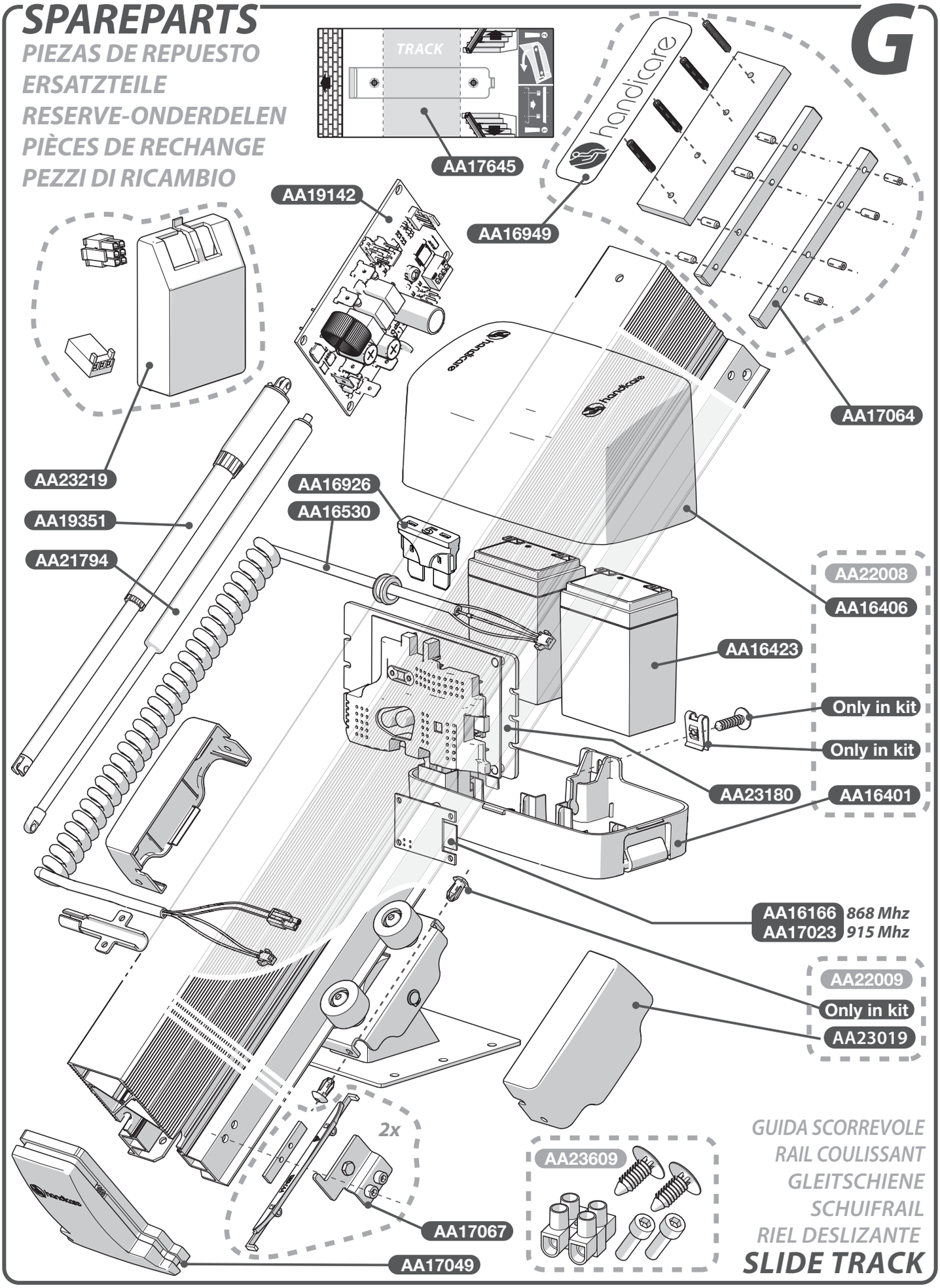
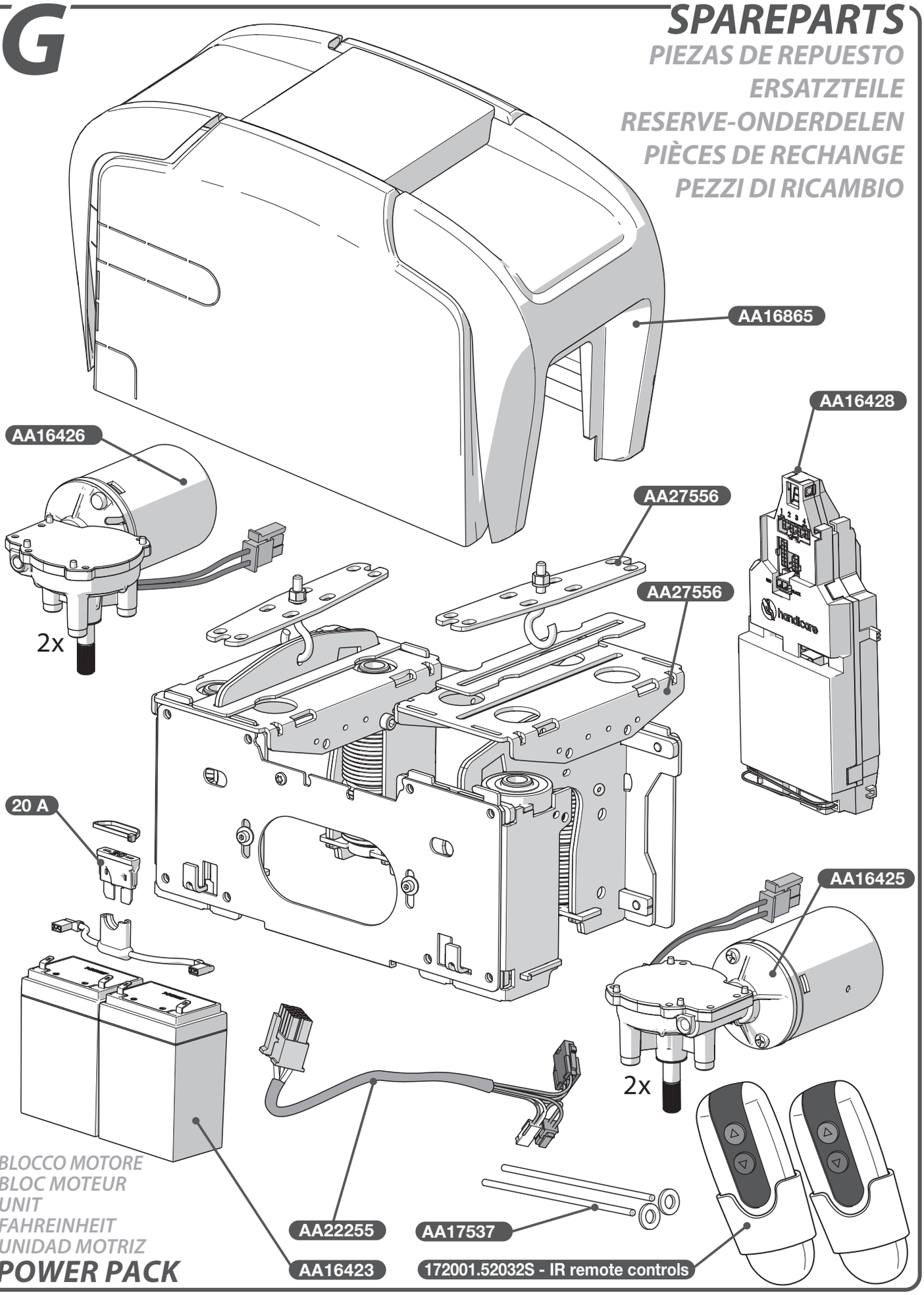
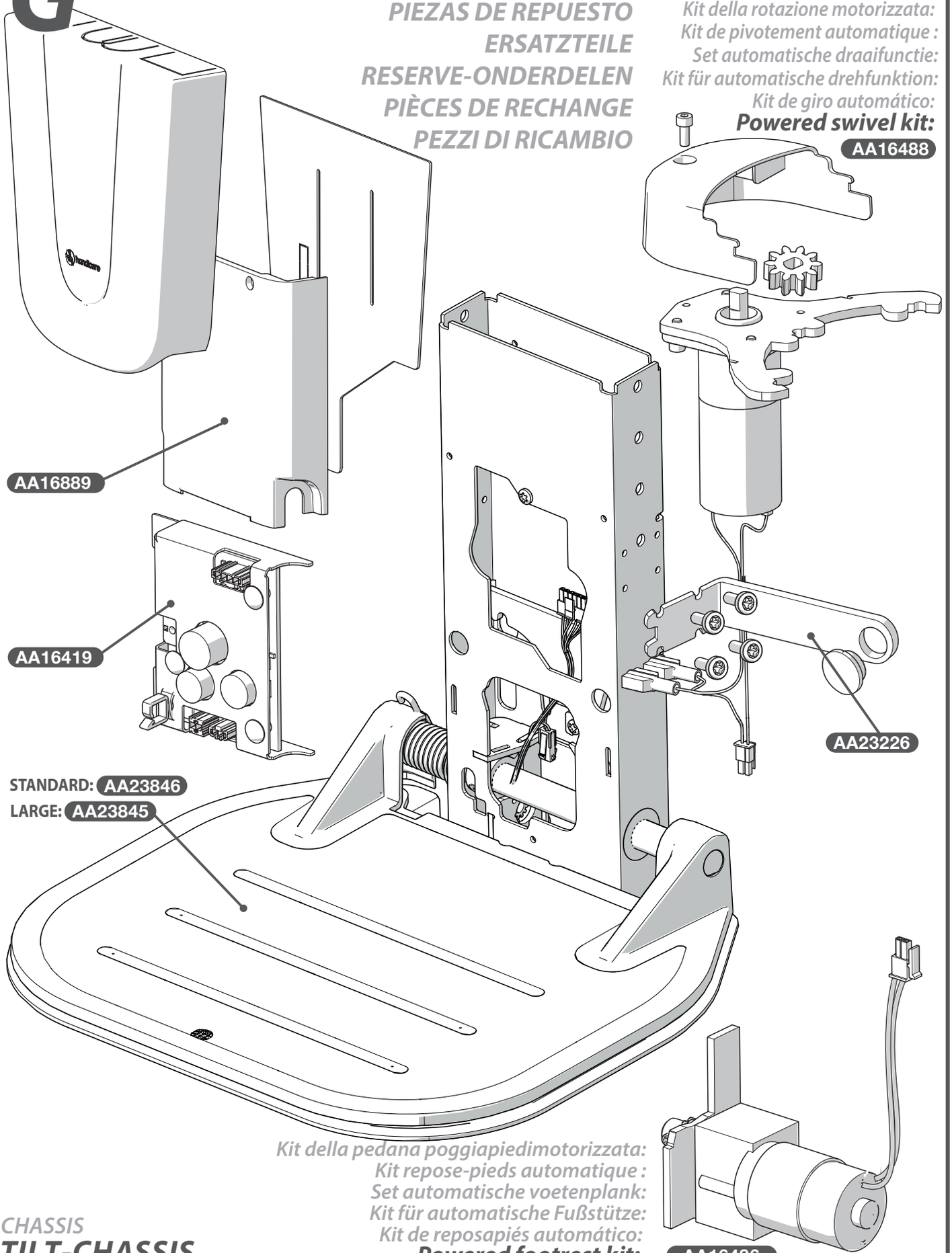
CHASSIS
TILT-CHASSIS

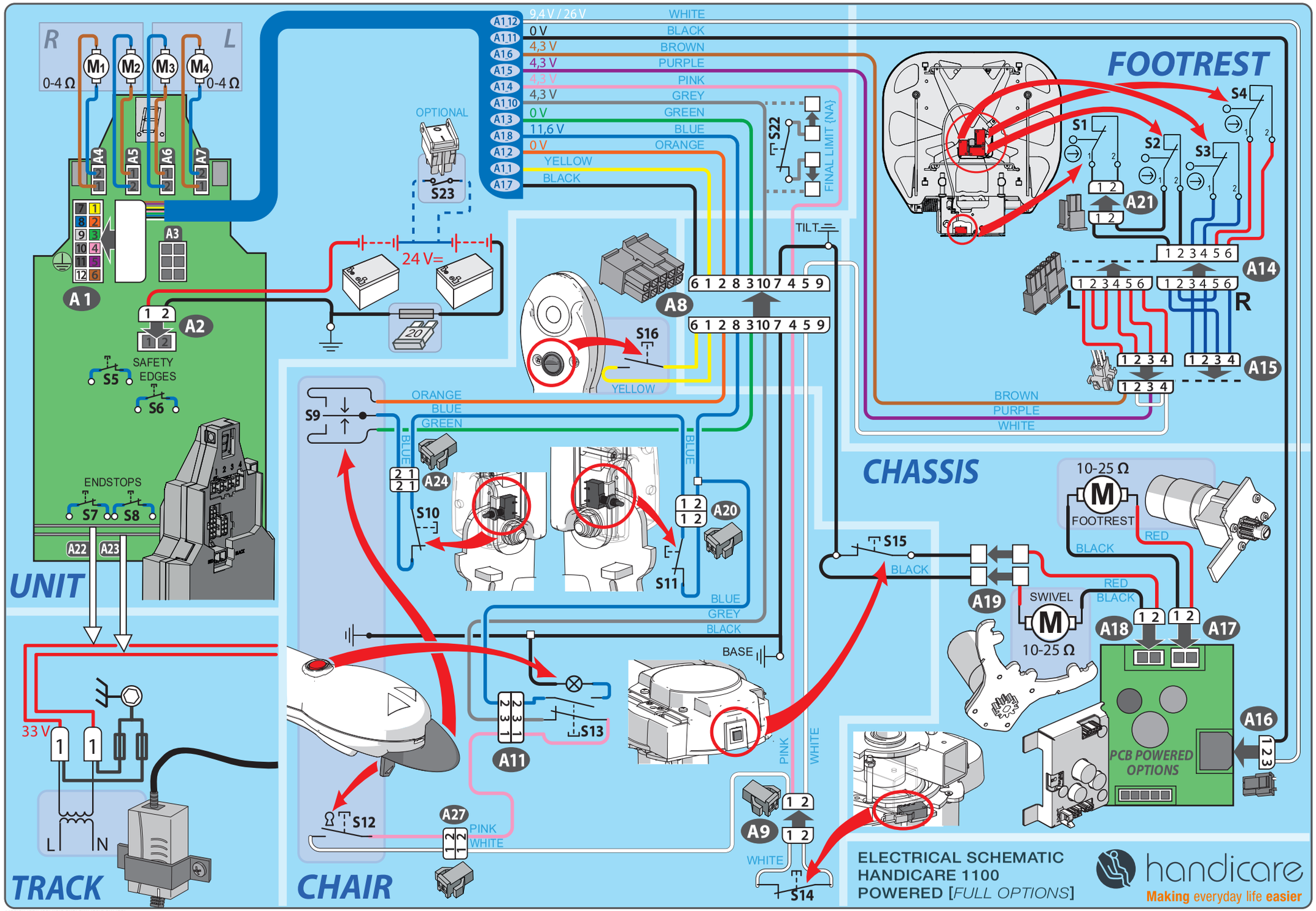For the latest version of this document check the Handicare partnersite.

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Handicare Canada
81 Romina Drive,
Concord, ON L4K 4Z9
Canada
T : 905.850.9003
F : 416.260.0256

www.handicare.com

SERVICE MANUAL
SERVICEHANDLEIDING
MANUAL DE SERVICE
KUNDENDIENSTANLEITUNG
MANUALE PER SERVIZIO
MANUAL DE SERVICIO
MANUAL DE MANUTENCAO
VEDLIKEHOLDSHANDBOKEN
UNDERHALLSMANUAL
INSTRUKCJA

AA17833_rev 21_090622