www.drivetech.info


RTL10353 INSTALLATION

www.drivetech.info


1. Verify that the tip will fit the cane.



2. Remove the Allen head wedge from the locking wedge.

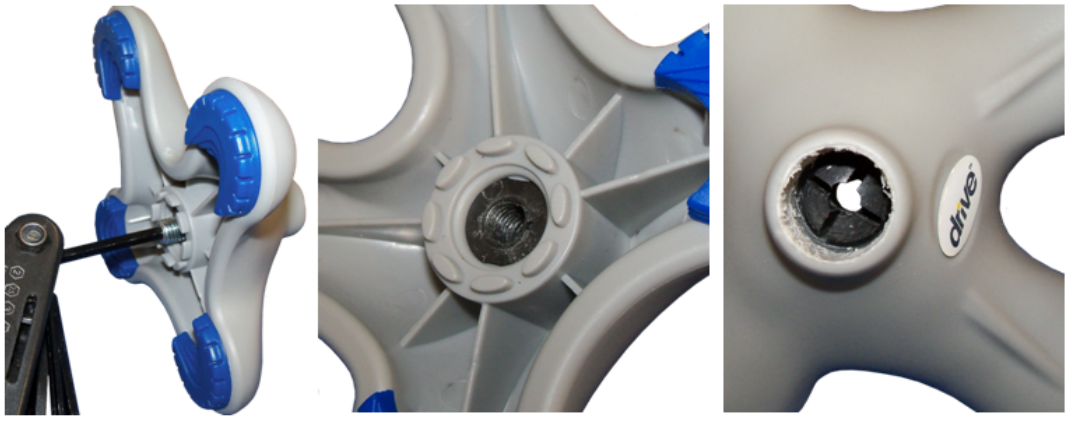

3. Remove the existing tip from the cane. Mark off in pencil 1 1/8” from the bottom of the cane.

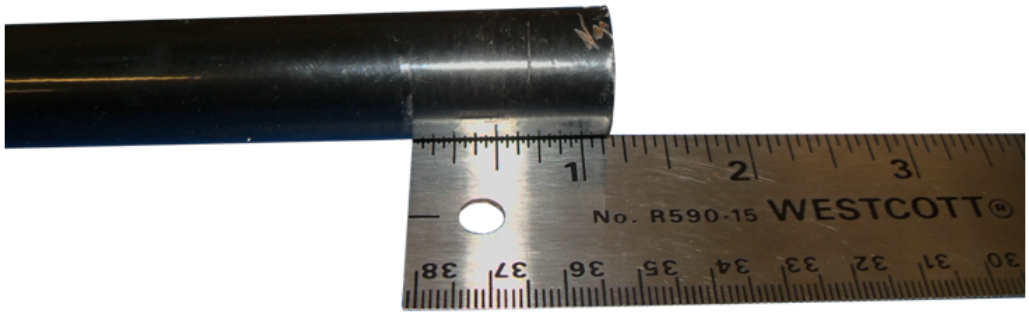

4. Place the Quad Tip on the ground and insert the can into the tip.



5. Push the cane into the tip by leaning into the cane. The cane has to go in at least 1”. Verify this by checking the pencil make that was put on the end of the cane.





6. Insert the Allen head wedge into the locking wedge and tighten.



