



Buzzaround GB120



The Buzzaround CarryOn is small in stature and big on fun! This super portable scooter features a new sporty blacked-out style and upgraded ergonomic seat for ultimate comfort. It features a 300 lb. weight capacity, an infinitely adjustable tiller, with a convenient USB port, drink holder, cane holder, front and rear LED lighting, and flip-up armrests. Enjoy long range coverage up to 22 miles on the airline friendly 12.5 AH and spare 6.5 AH lithium-ion batteries, which can be removed and taken as a carry on while traveling by air, sea, or train. A kickstand with a rear basket allows for vertical storage and serves as spare battery holder. The heaviest piece is 47.5 lbs. (without seat and batteries) and the height when folded is only 12". It disassembles easily into two pieces and is designed to store in the smallest compartments. It comes standard with red color shrouds installed and an extra set in blue. Optional color shroud kits let you express a style all your own! Pack your dreams and go anywhere with the CarryOn Scooter!



Standard Colors:



Velvet Red
(Installed)



Coastal Blue
(Included)

Optional Colors:



Crisp White



Canary Yellow



Stone Grey



Orchid Pink



Army Green
(Matte)



Desert Tan

Sporty Fold-Flat Scooter! Airline Friendly! Optional Travel Bag!

Standard Features:

- 300 lb. Weight Capacity
- 4-Point Carry Handles for Easy Transport
- 12" Height When Folded (without seat)
- 4" Ground Clearance
- 12.5 AH Lithium-ion Battery with Travel Case & 6.5 AH Spare
- Flip-Up Armrests
- Infinite Adjustable Tiller
- USB Port, Drink Holder & Cane Holder
- Disassembles Easily

**FREE Spare 6.5 AH
Lithium-ion Battery for Up
to 22 Mile Combined Range
with 12.5 AH Battery!**

2-Year Battery Warranty!




Model	GB120
Class I Medical Device	Yes
Weight Capacity	300 lbs
Drive Wheels	Rear
Maximum Speed	3.6 mph
Operating Range ¹ 12.5 AH	14.4 miles
Operating Range ¹ 6.5 AH Battery	7.5 miles
Total Combined Range ¹	21.9 miles
Ground Clearance Under Deck	4"
Ground Clearance to Transaxle/Motor	1.5"
Turning Radius	48.75"
Battery Type (Lithium-Ion)	12.5 AH & 6.5 AH
Colors	Red, Blue, White, Yellow, Grey, Pink, Green & Tan
Freewheel Mode	Yes
Electronic Speed Control	Yes
Electro-Mechanical Brakes	Yes
Charger	Off Board
Controller	50 Amp
Dimensions	
Length	38"
Width	19.5"
Height (ground to top of back of seat)	34.5"
Maximum Height (folded without seat)	12"
Ground to Top of Seat	23.5"
Ground to Top of Deck	6"
Deck to Top of Seat	17.5"
Front Axle to Rear Axle	28.7"
Weight of Unit: (assembled)	65.5 lbs
Main Body	47.5 lbs
Seat w/Arms	12 lbs
Batteries	5 lbs each
Tires:	Solid/Low Profile
Front	6"
Rear	7"
Rear Anti-Tip Wheels	2"
Standard Seat	Stadium Style
Back Height - no headrest	13.25"
Width x Depth	15" x 13.75"
Color	Black Vinyl
Warranty:	
Frame	Lifetime
Drive Train	2 Years
Electronics ²	2 Years
Batteries	2 Years

¹ Battery range at 200 lbs., and will vary due to rider weight, drive surface, terrain, and battery type.

² Electronics warranty excludes batteries.

LIGHTWEIGHT PORTABLE SCOOTER

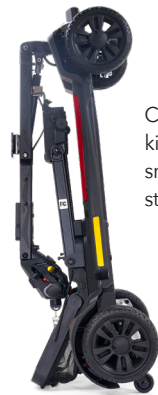

Convenient kickstand for small space storage



Folds easily for a compact profile of only 12" high (without seat)



Spare battery storage in rear basket



Front & rear LED lighting



Front carry handle makes travel easy



LED battery indicator on the console



Optional Travel Bag

Literature is current at time of printing. Golden Technologies reserves the right to make changes to the product or literature at any time.



Phone: 800-265-9615
www.goldentechnologies.ca

Find us on
Facebook
www.facebook.com/goldentechcanada

LIT-SHT-GB120-CAN-0326