

HEADQUARTERS

Apex Medical Corp.
No. 9, Min Sheng St., Tu-Cheng,
New Taipei City, 23679, Taiwan
T +886 2 2268 5568
F +886 2 2268 6525
marketing@apexmedicalcorp.com

PAN-ASIA

TAIWAN

Apex Medical Corp.
No. 9, Min Sheng St., Tu-Cheng,
New Taipei City, 23679, Taiwan
T +886 2 2268 5600 /0800 054 133
F +886 2 2268 6526
customer@apexmedicalcorp.com

STURDY INDUSTRIAL

168, Sec.1, Zhongxing Road, Wugu District,
New Taipei City 24872, Taiwan
T +886 2 2981 5245
F +886 2 2984 8208
foreign@sturdy.com.tw
www.sturdy.com.tw

CHINA

Apex Medical (Kunshan) Corp.
No. 1368, ZiZhu Road, Kunshan,
Jiang Su Province, China
T +86 512 8617 1660 Ext 1150
F +86 512 8617 8663
customer.cn@apexmedicalcorp.com

THAILAND

Apex Medical (Thailand)
No.111/152 Soi Pho Kaeo Yak 19,
Nawamin, Beung Kum,
Bangkok 10240,
Thailand.
T +662 018 2140 Ext 4210
sales.thai@apexmedicalcorp.com

MIDDLE EAST • AFRICA

29, Misr Liltameer Building,
Stree#10, AREA#8, Sheraton Building
Heliopolis Cairo, Egypt.
T +202 2 2671303
F +202 2 2678649
raafat.wahby@apexmedicalcorp.com

PAN-EUROPE

SPAIN

Apex Medical S.L.
Elcano 9, 6ª planta 48008
Bilbao. Vizcaya. Spain
T +34 94 4706408
F +34 94 4706409
info@apexmedical.es

UK

Apex Medical Limited
Unit 33, Great Western Business Park
McKenzie Way, Worcester ,WR4 9GN
United Kingdom
T +44 (0) 190 5774 695
sales@apexmedicalcorp.co.uk

FRANCE

Apex Medical France
4 Boulevard de la Chanterie
49124 Saint-Barthelemy d'Anjou
France
T +33 2 72 79 27 13
F +33 2 85 73 07 67
info.france@apexmedicalcorp.com


GERMANY

SLK Vertriebsgesellschaft mbH
Am Herdicksbach 18, 45731 Waltrop
Germany
T +49 (0) 231-925360-0
F +49 (0) 231-925360-29
info@slk-gmbh.de
www.slk-gmbh.de

AMERICAS

USA

Apex Medical USA
927 Mariner Street,
Brea, CA 92821, USA
T +1 714 671 3818
F +1 714 494 8173
sales.usa@apexmedicalcorp.com

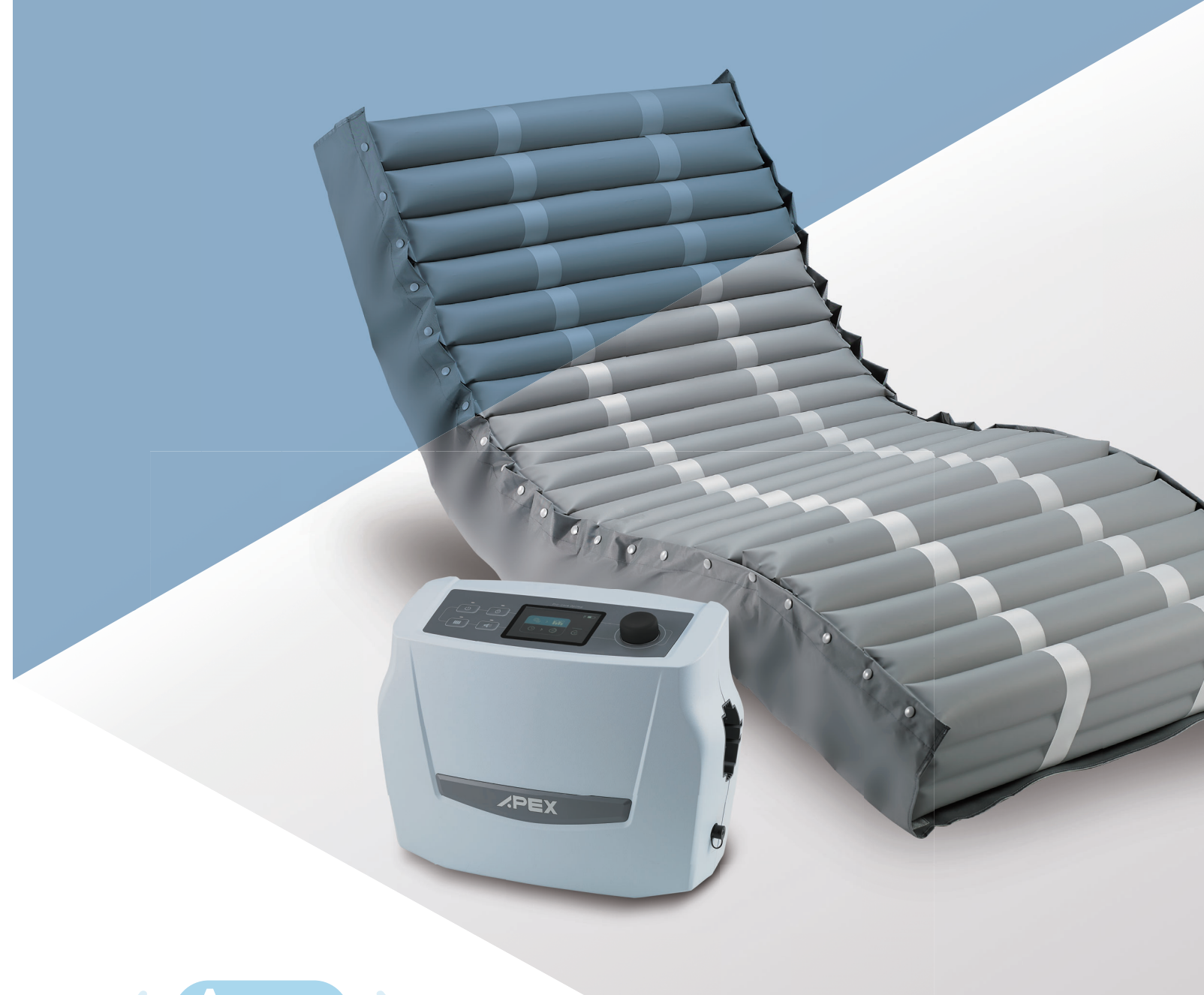

**AUTO
SENSING**

Pro-care Optima

Dynamic Pressure Redistribution Air Mattress

- Suitable for patients in high to very high risk of pressure injuries
- Auto-sensing technology with intuitive one touch set up
- Multi-zone air pressure for greater pressure redistribution and optimal comfort
- Designed for improving work efficiency while performing the best pressure care to patients



Distributed by


Pro-care Optima

Dynamic Pressure Redistribution Air Mattress

"Simplicity Redefined"


Multi-zone Air Pressure

It demonstrates significant pressure redistribution to fit different patient contours with stability and comfort by providing **different air pressure levels** in static head pillow, torso and lower leg sections

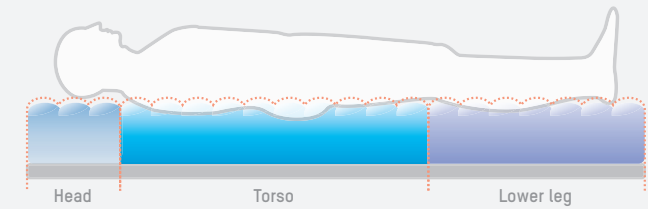


Auto Seat Inflation
Automatically initiates when the head section has been raised by ≥ 30 degrees to prevent bottoming out


80%

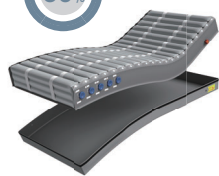

80% of Time Saving
During decontamination, air cells with quick release detachable base can be easily removed from mattress

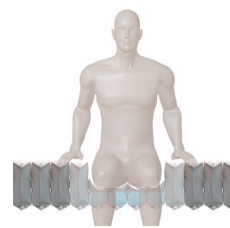

6 Cell-in-Cell Air Cells
Providing excellent support in sacral area for patient when in sitting position or sitting on the side before exiting from bed

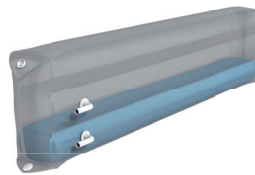

Micro Low Air Loss
Provide good ventilation and help reduce the accumulation of heat and moisture



Visual and Audible Alarm

- Low Pressure
- Power Failure
- Tech Support



CPR Knob

Mattress can be quickly deflated in less than 15 seconds for performing CPR to the patient



Cable Management

Reducing tripping hazards for a safer healthcare environment



4-way Stretch Cover

- Low friction and shear
- Vapour permeable
- High-flame retardant
- Water resistant



Water Resistance



4-Way Stretch




Vapour Permeable



Flame Retardant



Heel Relief Function
Single air-cell deflation can float the heels with **zero pressure**

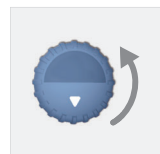


Activating heel protection with the one-handed control knob



Transport Mode/ Built-in Battery (Optional)

Non-stop dynamic pressure relief for more than **6 hours** during power failure or transport (24 hours for static mode without battery)

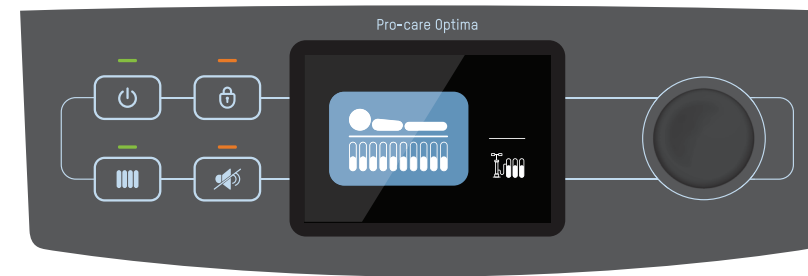

One Touch Set Up

- System automatically runs on default programme (Alternating therapy in 10-minute cycle)
- Just one touch, easily set up therapy or parameters

Built for Clinical Environment

The robust pump unit with smooth surface and less cavities is easier for daily cleaning and disinfection process in busy clinical environment



Intuitive Auto-Sensing Technology

- **Auto-recognition:** automatically initiates the corresponding pressure setting to different mattress
- **Auto-adjustment:** ensures that the pressure is constantly adjusted in accordance with the patient's weight and position

Specifications

Model	Pro-care Optima	8" Replacement
Pump	Dimensions	13.4" x 6.5" x 10.2" / 34.1 x 16.5 x 26.0 cm
	Weight	4.5 kg / 9.9 lb, 5 kg / 11.0 lb (With battery)
	Case Material	Flame retardant ABS
	Supply Voltage	220-240V / 50 Hz ; 110-120V / 60 Hz
Mattress	Operating Cycle	10, 15, 20, 25 min.
	Dimensions	78.7" x 36"/34"/32" x 8" / 200 x 80/85/90 x 20 cm
	Cell Height	8" x 21 cells
	Weight	9.2 Kg / 20.2 lb
	Cover Material	4-way stretch PU / Polyester cover (Optional: Nylon / PU quilted cover)
	Cell Material	TPU
	CPR Knob	1 Piece
Maximum Patient Weight	250 kg / 550 lb	
Flame Retardant Standards	EN 597-1 & EN 597-2	

Specifications are subject to change for improvement without notice