

Model Number	3W: S66 4W: S67
Weight Capacity	400 lbs.
Maximum Speed <sup>1</sup>	Up to 5 mph
Ground Clearance <sup>3</sup>	3W: 3.25" at side deck 4W: 3.5" at side deck
Turning Radius <sup>3</sup>	3W: 40" 4W: 52"
Overall Length <sup>3</sup>	3W: 46" 4W: 47"
Overall Width <sup>3</sup>	20.75"
Seat-to-Ground Height	21"-23.5"
Seat-to-Deck Height	15"-18"
Front Tires	3W: 3" x 8" 4W: 3" x 9"
Rear Tires	3" x 9"
Suspension	Front and rear (CTS)
Range Per Charge <sup>1,2</sup>	3W: Up to 17.4 miles (200 lbs.) Up to 13.2 miles (400 lbs.) 4W: Up to 17.8 miles (200 lbs.) Up to 12.4 miles (400 lbs.)
Total Weight Without Batteries <sup>5</sup>	3W: 120 lbs. 4W: 132 lbs.
Heaviest Piece When Disassembled <sup>3</sup>	3W: 48 lbs. (rear section) 4W: 57 lbs. (front section)
Standard Seating	Type: Compact, foldable Dimensions: 20" width, 18" depth
Drive System	Rear-wheel drive sealed transaxle with 24-volt DC motor
Dual Braking System	Regenerative and electromechanical
Battery Requirements <sup>4,6</sup>	Size: (2) U1 batteries; 25 lbs. each
Battery Charger	3.5-amp, off-board
Warranty	Frame: lifetime limited Drivetrain: 2-year limited Electronics: 2-year limited Battery: 6-months

<sup>1</sup> Varies with user weight, terrain type, battery amp-hour (AH), battery charge, battery condition, and tire condition. This specification can be subject to a variance of (+ or -) 10%.

<sup>2</sup> Tested in accordance with ANSI/RESNA, MC Vol2, Section 4 & ISO 7176-4 standards. Results derived from theoretical calculations based on battery specifications and drive system performance. Test conducted at maximum weight capacity.

<sup>3</sup> Due to manufacturing tolerances and continual product improvement, this specification can be subject to a variance of (+ or -) 3%.

<sup>4</sup> Battery weight may vary based on battery manufacturer.

<sup>5</sup> Includes standard seat weight.

<sup>6</sup> AGM or gel cell type required.

**PRIDE**<sup>®</sup>  
MOBILITY  
LIVE YOUR BEST<sup>®</sup>

Beamsville, Ontario • 888-570-1113

pridemobility.ca |    

NOTE: All specifications subject to change without notice. The information contained herein is correct at the time of publication; we reserve the right to alter specifications without prior notice. Speed and range vary with user weight, terrain type, battery charge, battery condition and tire pressure. Due to manufacturing tolerances and continual product improvement, length, width, turning radius and ground clearance can be subject to variance of + or - 3%. Special order seats or non-standard dimensions will not be fold-flat and depth-adjustable.

0001667-070620

**REVO<sup>®</sup> 2.0**


**RUGGED**  
WITH  
**FEATHER-TOUCH**  
**DISASSEMBLY**



**LIVE YOUR BEST<sup>®</sup>**  
pridemobility.ca |    

# TOUGH AND BUILT TO LAST

20" x 18" Pinchless hinge, black vinyl seat with 8-point swivel for increased comfort

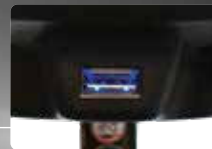

USB port in tiller

Standard front basket

Standard under seat storage

Easy-to-grip tiller adjustment knob

Exclusive stylish, lightweight, non-scutting, black, low-profile tires  
(Optional pneumatic tire on the 4W: S67)



Pride's CTS Suspension



Feather-touch disassembly



LED headlight in tiller

Experience the **Revo® 2.0**, the most impressive mid-size scooter that combines an array of features into one great package. Tough and built to last, the Revo 2.0 offers rugged dependability you expect from luxury mid-size scooters. **Comfort-Trac Suspension (CTS)** ensures a smooth and comfortable ride over varied terrain. Experience convenient, feather-touch disassembly for portability that is comparable to our **Travel Mobility** line. Enjoy standard under seat storage and speeds up to 5 mph.

- Exclusive one-hand feather touch disassembly
- Easily charge the batteries on or off the unit
- ECO mode/sport mode selector switch offers an optional ECO mode at 4 mph with improved range per charge or the sport mode at 5 mph
- Delta tiller with ergonomic wraparound handles lets you operate the scooter with one hand and rest your wrist
- Easy-to-read LED battery meter
- Ultra heavy-duty drivetrain with 70-amp controller
- User-friendly console
- Easy-access tiller-mounted XLR USB charging port for mobile devices

Available Colours



Grey Street



True Blue

## Accessories

- Cane/crutch holder (single, double)
- Cell phone holder
- Cup holder
- Forearm crutch holder
- Oxygen tank holder
- Quad cane holder
- Rear basket
- Rearview mirror
- Saddlebag
- Safety flag
- Walker holder
- Weather cover
- Wishbone crutch holder
- Lap belts: 50", 60" and 70"