

# Folding Wheelchair

Ultra Lightweight





**FUTURE MOBILITY  
FOLDING WHEELCHAIR**



Stellato II is an innovative, low maintenance, high performance wheelchair built for comfort and affordability. It's bursting with standard features that make it the most desired wheelchair in its category. At a weight of 31 lbs., the Stellato II Wheelchair is a strong, lightweight option for multi-purpose use.



## Features

- Multi position axle plate
- Composite flip-up footplates
- Wheel locks (push or pull to lock)
- Elliptical heavy duty cross brace tubing
- Quick-flip anti-tips
- Seat upholstery slip-on back upholstery
- Angle-adjustable push handles
- Heel loops
- Height adjustable-angle caster housings with - precision bearings
- Lightweight folding frame with swing away footrests
- Full-length, fixed-height, flip-back vinyl wrap around arms

Canadian Made  
folding wheel-  
chair with quality  
and durability  
that defines  
Future Mobility

## ADJUSTABLE BACK CANES

Set height of adjustable back canes between -5 and 21 degrees.





## UNIQUE CROSS-BRACE

Single aluminium cross-brace for optimum rigidity and ease of folding. Elliptical heavy-duty cross brace tubing creates a structurally sound frame.



## HEIGHT ADJUSTABLE AXLE PLATE

The axle plate ensures you are getting the maximum adjustment range you need without changing parts or adding unnecessary weight.

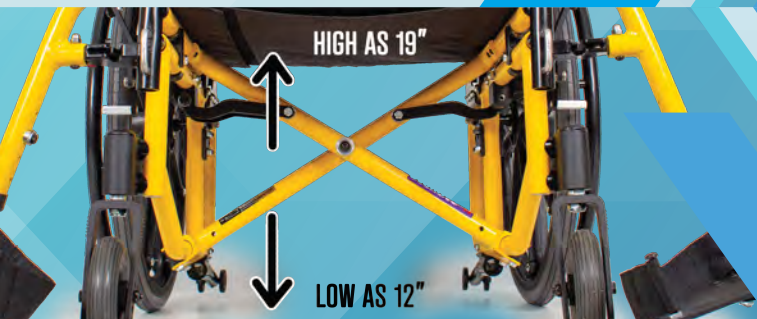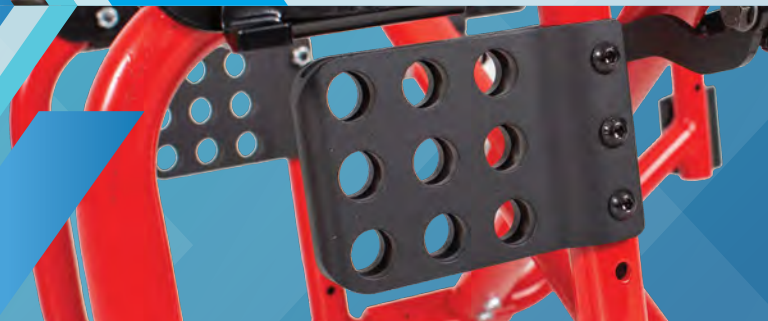

## HEIGHT ADJUSTABLE JOURNAL

Allows for a seat-to-floor range of 12"-19" on the frame. It features dual bolt adjustment, which allows for easy adjustability. Journals can be positioned in a leading or trailing position, allowing for a stable footprint or a tighter turning radius.

Future Mobility is a North American manufacturer of wheelchairs, backs, cushions and accessories with our parts manufactured in house. This gives us the ability to customize anything needed for specific customer requests.

**CUSTOMIZE  
ANY PART!**



**Seat Width**  
14" - 20"



**Seat Depth**  
15" - 20"



**Adjustable Back Angle**  
-7°, 0°, 7°



**Caster Size**  
4", 5", 6", 8"  
Urethane



**Front Rigging**  
Swingaway, Dual Swingaway,  
or Elevating Legrest



**Back Style**  
19" - 22" Low  
22" - 25" High  
Slip on Back Upholstery



**Weight Capacity**  
275 lbs. (125 kg)



**Seat to Floor**  
13" - 20"



**Rear Wheel**  
20", 22" or 24"  
Composite Urethane



**Product Weight**  
31 lbs. (14 kg)



**Overall Width**  
+ 9"

**ADP Code**  
WM0001188

**HCPCS Code**  
K0004

## Warranty:

Future Mobility Healthcare Inc. warrants the Stellato II wheelchair to be free from defects in materials and workmanship for two (2) years from the date of purchase. The frame has a lifetime warranty and the upholstery for 90 days, upon normal usage by the original purchaser. If within this warranty period the product shall be proven to be defective, such product shall be repaired or replaced, at Future Mobility Healthcare Inc. discretion.

## USA | CANADA

1 Buffalo River Place,  
Buffalo  
NY 14210  
716-783-9130  
www.futuremobility.com

6750 Professional Court  
Mississauga  
ON L4V1X6  
1-888-737-4011  
www.futuremobility.ca